



Innovation and quality in the field
of foot care and podiatry

Catalogue 2017

20 - 23

Complete Offers

24 - 47

Treatment chairs

28 - 35

Podo Xdream series

36 - 47

Pioneer series

48 - 89

Furnitures

50

Cabinets

70

Work chairs

78

Magnifying lamps

86

Foot baths

90 - 133

Treatment systems

97

Studio Line

100

Suction and wet technology

104

Compressed air technology

111

Flexi Line

112

Suction technology

118

Wet technology

129

Accessories

134 - 151

Mobile foot care

136

Trolley system SÜDA mobility

140

SÜDA cases

148

Case lights

149

Legrests

150 - 151

Accessories

152 - 201

Instruments

154

Cutter and grinder

178

Nippers and scissors

174, 198

Accessories

202 - 217

Hygiene

204

Ultrasound equipment

206

Sterilization technology

208

Hot Air-Sterilizatores

210

Autoklaven

214

Accessories

216

Hygiene Sets



SÜDA

Your partner and pioneer for podiatry and foot care

A light-footed walk is emblematic of our physical wellbeing, the foot as a whole is always in the tradition of our company. We will be a lifelong competent dialogue and consulting partner when it comes to support you with the most innovative and multiply awarded technology in your work: from the beginning.

SÜDA dedicated nearly 100 years of development of exclusive technologies for maintaining the health of the foot and provides a comprehensive high quality product range and a concerted practice facilities. Maintaining the health of the foot is the center of our research. We are worldwide in conjunction with innovative tradition for durable modernity: Made in Germany.

It's easier with SÜDA.





From the way to the destination

Freedom, well-being and independence, which means mobility into old age. Who is doing something good for his foot, strengthen his immune system and his health. Chiropody and podiatry are indispensable part of the health care and serve each individual and the society as a whole. We support you and your customers in a careful handling of one's own body to increase well-being. Let us give our feet special attention!

SÜDA accompanies you.

Vision

SÜDA dedicated to holistic health, the environment and society. We encourage young talent through education, training and competitions and by the perception of social and corporate responsibility.

Design, development and production of SÜDA devices were made directly by our company and are perfectly, not only to the specific needs of our customers, but also integrate industry standards and norms.

SÜDA produces energy-saving and clean and our particularly climate-friendly shipping bears the corporate responsibility bill.



German production & quality service

From raw materials to finished treatment chair, from the board to the finished device: our products are manufactured by our employees at the company's own production site in Leipzig. We assure the quality in all areas of our company through the annual certification to the internationally recognized quality management systems ISO 9001 and ISO 13485 for medical devices. Delivery of goods at home and abroad rely on a daily basis, our modern logistics center as well as a variety of bulk goods.

Our all-round carefree package includes on demand: supply, installation, and professional training in the functionality of the goods at your premises by employees in six offices nationwide. We take a two year warranty on our technical branded products and offer reliable security. Our customer service department is available for comprehensive service and fast, competent management in maintenance, examination or repair. And for a failure free course of treatment during the service period, we offer a free loan unit, whenever possible.



Complete Offers

20

Offer Compact

21

Offer Profi

22

Offer Hightech

23

Offer Podiatry

Offer Compact

Elegant functionality for professional treating

INSTITUTE EQUIPMENT

8533	Foot care chair Podo Xdream 2M	> P. 30
1184-12300	Cubus 400	> P. 54
8206	Magnifying lamp Standard	> P. 80
1946	Pedo Sprint 50 SP	> P. 114
8003-101200	Round chair with banana rest, color white, foot cross in chrome look	> P. 72



Offer Profi

Flexible treatment efficiency for
podiatry and foot care

INSTITUTE EQUIPMENT

8534	Foot care chair Podo Xdream 3M	> P. 34
1188-11230300	Cubus 800	> P. 56
8206.40	Magnifying lamp IQ LED	> P. 81
1997	SÜDA vac S	> P. 116
1259.21	SÜDA spray L	> P. 127
8202.4	Round chair FIT, white	> P. 77



Offer Hightech

Top features of top-class for demanding work - Perfect base for maximum treatment comfort

INSTITUTE EQUIPMENT

8540	Foot care chair Pioneer, 3 motors, color white	> P. 38
8206.72	Magnifying lamp Circle LED	> P. 82
1231-00	Treatment unit Spray Vac Unit	> P. 102
8202.5	Round chair Rodeo, color white	> P. 76



Offer Podiatry

Innovative treatment technology with the highest comfort for demanding applications

INSTITUTE EQUIPMENT

8542U14	Foot care chair Pioneer 5 motors, color Tomato	> P. 42
1143	Magnifying lamp Tevisio	> P. 83
1229-00	Compressed air treatment unit PodoAIRspeed, with Turbine handpiece, air motor handpiece, Three-function nozzle, UV compartment	> P. 106
1252.2	Compressor cabinet white, 30 liters	> P. 109
8002-1012U14	Round chair 360, banana rest, color tomato	> P. 72



Treatment chairs

28-35

Podo Xdream-
Line

36-47

Pioneer-Line

Podo Xdream 2M/3M



Podo Xdream 2M
> P. 30

Podo Xdream 3M
> P. 34

Classic-Contour
padding colors



(white)



U01 (cement)



U13 (terre)



U21 (prune)



U14 (tomato)



U23 (absinthe)

CATEGORY	Podo Xdream 2M (footrests manually)	Podo Xdream 3M (footrests manually)
MASSAGE SUITABILITY	no	yes
POWER UNIT	2 electric motors	3 electric motors
CONTROL	Foot switch	Foot switch
LENGTH (CM)	174 - 198cm	174 - 198 cm
WIDTH	56 cm without, 78 cm with armrest	56 cm without, 78 cm with armrest
SEAT DEPTH	47 cm	47 cm
WEIGHT	85 kg	85 kg
MAX. SEAT-HEIGHT	87 cm	87 cm
BACK ADJUSTMENT	27°	80°
SEAT ANGLE ADJUSTMENT	27°	27°
FOOTREST ADJUSTMENT	90°	90°
HEIGHT ADJUSTMENT (CM)	30 cm	30 cm
ANGLE OF ROTATION	0°	0°
LEG RESTS	segmented, extendable, height adjustable	segmented, extendable, height adjustable
MAX. WORK HEIGHT	125 cm	125 cm
MAX. LOAD ON SEAT	170 kg	170 kg
ENTRY LEVEL	57 cm	57 cm
HEADREST	extendable	extendable

Pioneer 3M/5M



Pioneer 3M
> P. 38

Pioneer 5M
> P. 42

CATEGORY	Pioneer 3M (footrests manually)	Pioneer 5M (footrests manually)
MASSAGE SUITABILITY	yes	yes
POWER UNIT	3 electric motors	5 electric motors
CONTROL	Foot switch	Foot switch
LENGTH (CM)	163 - 203 cm	163 - 203 cm
WIDTH	56 cm without, 78 cm with armrest	56 cm without, 78 cm with armrest
SEAT DEPTH	50 cm	50 cm
WEIGHT	110 kg	110 kg
MAX. SEAT-HEIGHT	100 cm	100 cm
BACK ADJUSTMENT	80°	80°
SEAT ANGLE ADJUSTMENT	27°	27°
FOOTREST ADJUSTMENT	90°	90°
HEIGHT ADJUSTMENT (CM)	48 cm	48 cm
ANGLE OF ROTATION	180°	180°
LEG RESTS	segmented, extendable, height adjustable laterally folded away	segmented, extendable, height adjustable laterally folded away
MAX. WORK HEIGHT	139 cm	139 cm
MAX. LOAD ON SEAT	200 kg	200 kg
ENTRY LEVEL	56 cm	56 cm
HEADREST	extendable	extendable

Cushion-Contour
padding colors



(white)



D1 (cement)



D2 (terre)



D3 (pearl)



Podò Xdream

Podo Xdream 2M

With an entry height of only 57 cm and a possible working height of up to 125 cm, it offers you and your customers comfort. In addition, the electric adjustment of the seat and the synchronous adjustment of the seat and back. The segmented, extendable and adjustable leg rests round off the perfect functionality.



Electric adjustment of the seat height

Synchronous adjustment of seat angle and backrest via footswitch

Armrest for lateral entry folded backwards

Extractable and removable headrest

Divided leg supports, extractable and height adjustable

Available in 6 padding colors



Podo Xdream 2M

The Podo Xdream with its new, ergonomic Classic Contour padding provides thanks to its 10 cm thick, breathable cushion with foam core and shell seat shape best seating comfort. The particularly durable, easy grained surface is resistant to heavy constant load at high treatment volume.

Resistant to disinfectants and detergents, it provides high hygienic safety.

TECHNICAL SPECIFICATIONS

Power consumption: 230 V / 50-60 Hz
 Power: 105 VA
 Length: 174 - 198 cm
 Headrest: 12 cm (extendable)
 Footrest: 12 cm (extendable)
 Width with armrests: 78 cm
 Backrest: 57 cm (length)
 Padding width: 55 cm
 Seat depth: 47 cm
 Seat height: 57 - 87 cm
 Back adjustment: 27°
 Seat angle adjustment: 27°
 Footrest: 90°
 Max. Working height: 125 cm
 Weight: 85 kg
 Max. Seat load: 170 kg

8533	Podo Xdream 2M
------	----------------

ACCESSORIES

8301.119	Terry cloth Podo Xdream, 8-pieces, white
8301.150	Soft Velvet cover Podo Xdream, 8 pieces, white
8301.1500	Soft Velvet cover for armrests, 1 pair
8530.101	Foot section protector Podo Xdream, 1 pair
8299.63	TEP mouss, upholstery Cleaner 150ml
8299.61	Microfiber cloth, 40 x 40 cm



Classic-Contour
 Polsterfarben



EQUIPMENT OPTIONS

8533-XXX*	Podo Xdream 2M, Classic-Contour-padding, color white, Standard colors*
-----------	--

* Please insert the corresponding color number at the end of the item number
 Classic-Contour padding colors: U01 (cement), U13 (terre), U21 (prune), U14 (tomato), U23 (absinthe); White is the standard color and do not require an additional specified color number

Podo Xdream 3M

3-motors model with electric adjustment of seat height, seat angle and separately adjustable backrest.

TECHNICAL SPECIFICATIONS

Power consumption: 230 V, AC, 50 - 60 Hz
 Power: 120 VA
 Length: 174-198 cm
 Headrest: 12 cm (extendable)
 Footrest: 12 cm (extendable)
 Width with armrests: 78 cm
 Backrest: 57 cm (length)
 Padding width: 55 cm
 Seat depth: 47 cm
 Seat height 57 - 87 cm
 Back adjustment: 80 °
 Seat angle adjustment: 27 °
 Footrest: 90 °
 Max. Working height: 125 cm
 Weight: 85 kg
 Max. Seat load: 170 kg

8534 Podo Xdream 3M

ACCESSORIES

8301.119 Terry cloth Podo Xdream, 8-piece, white
 8301.150 Soft Velvet cover Podo Xdream
 8301.1500 Soft Velvet cover for armrests, 1 pair
 8530.101 Foot section protector Podo Xdream, 2 pieces
 8299.63 TEP mousse, upholstery cleaner 150ml
 8299.61 Microfiber cloth, 40 x 40 cm
 8463.100 Massage headrest, white
 8463.100XXX* Massage headrest colored

EQUIPMENT OPTIONS

8534-XXX* Podo Xdream 3M, Classic-Contour-padding, white Standard colors*

* Please insert the corresponding color number at the end of the item number
 Classic-Contour padding colors: U01 (cement), U13 (terre), U21 (prune), U14 (tomato), U23 (absinthe); White is the standard color and do not require an additional specified color number

Separate electric adjustment of seat height, seat angle and backrest
 Retractable headrest and footrest
 Armrests in massage position foldable backwards
 Optional massage headrest





Pioneer

Pioneer 3M



The comfortable chiropractic chair offers a combination of high-quality stability criteria with technical highlights, successful design and a top price / performance ratio. This combination guarantees your customers the highest level of relaxed comfort and unique multifunctional and convenient treatment options.

- Separate electric adjustment of seat height, seat angle and backrest
- Padded armrests in massage folded position
- Fatigue-free seat comfort thanks to the body-contoured seat
- Rotating by 90 ° each left and right
- Retractable and removable headrest
- Leg supports can be swung outwards
- Multifunction foot control with 3-memory function
- Available in 6 padding colors

Vertical swing

Thanks to absolute vertical stroke, you do not need to adjust your working position even if you change the height adjustment.



Classic-Contour-pading

With its particularly light grained durable, surface, the Pioneer with Classic Contour padding is excellent for strong constant load at high treatment volume.

The 10 cm thick, breathable cushion with foam core in an ergonomic shell seat shape also provides perfect comfort. Best hygiene security guaranteed the water-proof padding with wet stop, which is resistant to disinfectants and cleaning agents.



Available in 6 colors



Cushion-Contour padding

The twelve centimeters thick, breathable and ergonomic shaped Cushion-Contour padding in a special mixed padding technology and in conjunction with the shell seat molding guarantees unique sitting and laying comfort.



Available in 4 upholstery colors



Pioneer 3M

Multifaceted and top treatment convenience through individual adjustment of the backrest, seat angle and height using multifunctional foot control with 3-fold memory.

TECHNICAL SPECIFICATIONS

Power consumption: 230 V, AC, 50 - 60 Hz
 Power: 105 VA
 Length: 163- 203 cm
 Headrest: 15 cm (extendable)
 Footrest: 25 cm (extendable)
 Leg rest: 34 - 59 cm (length)
 Width with armrests: 78 cm
 Back part: 52 cm (width)
 Padding width: 54 cm
 Seat depth: 50 cm
 Seat height: 56 - 100 cm
 Max. Working height: 139 cm
 Back adjustment: 80°
 Seat angle adjustment: 27 °
 Foot adjustment: 90 °
 Rotation angle: 180 °
 Weight: 85 kg
 Max. Seat load: 200 kg

8540 Pioneer 3M, base model



ACCESSORIES

8301.123	Terry cloth Classic Contour padding, 8 pieces, white
8301.124	Terry cloth Cushion-Contour padding, 8 pieces, white
8301.151	Soft Velvet cover Pioneer, 8 pieces, white
8301.1500	Soft Velvet cover for armrests, 1 pair
8540.100	Foot section protector, 4 pieces
8299.63	TEP mouss upholstery cleaner 150ml
8299.61	Microfiber cloth, 40 x 40 cm

ORDER SCHEME

		P	U/D
Model	Pioneer 3M, base model	854	
Padding/Motors = P	Classic-Contour-padding, 3 motors	0	
	Cushion-Contour-padding, 3 motors	1	
Padding color = U/D	Classic-Contour-padding white		U16
	Cushion-Contour padding white		
	Padding Standard colors*		XXX*

* Please insert the corresponding color number at the end of the item number (U / D):
 Classic-Contour padding colors: U01 (cement), U13 (terre), U21 (prune), U14 (tomato), U23 (absinthe);
 Cushion-Contour padding colors: D1 (cement), D2 (terre), D3 (pearl).
 White is the standard color and do not require an additional specified color number

Pioneer 5M



Perfect treatment efficiency thanks to adjustability of each individual treatment position by foot pressure. The legrests remain also with the fully automatic 5 motors model laterally rotatable. Each motor can be conveniently and individually adjusted by foot pressure or move thanks to memory function with only one foot print 5 motors simultaneously in the desired treatment position.

Separate electric adjustment of seat height, seat angle, backrest and leg rests

Padded armrests foldable in massage position

Fatigue-free seat comfort thanks to the body-contoured seat

Rotatable by 90 ° to the left and right

Retractable and removable headrest

Leg supports can be swung outwards

Multifunctional foot control with 3-way memory function

Available in 6 padding colors



5 motors for more treatment comfort. With the new, fully-automatic model the Pioneer offers even more comfortable conditions for concentrated and effective treatment. Completely controllable with the innovative multifunction foot switch the chair provides you perfect space for your hands and a professional treatment.

Pioneer 5M



The multifunction podiatry chair Pioneer now offers even more performance for maximum comfort and perfect efficiency during treatment.

With the 5 motors model, it is possible to adjust all positions automatically on the multifunction foot control: The convenient 3-way memory function drives your Pioneer with only a foot print fully in the treatment position you saved.

Multifunction foot switch



- Lifting and lowering of the chair and adjust the lying area by select key
- Display of the selected adjustment range
- Keys for moving up and down the selected area
- Memory buttons for direct selection of saved settings
- Memory storage key

ACCESSORIES

8301.123	Terry cloth Classic Contour padding, 8 pieces, white
8301.124	Terry cloth Cushion-Contour padding, 8 pieces, white
8301.151	Soft Velvet cover Pioneer, 8 pieces, white
8301.1500	Soft Velvet cover for armrests, 1 pair
8540.100	Foot section protector, 4 pieces
8299.63	TEP mouss upholstery cleaner 150ml
8299.61	Microfiber cloth, 40 x 40 cm
8542	Pioneer 5M, base model



TECHNICAL SPECIFICATIONS

Power consumption: 230 V, AC, 50 - 60 Hz
 Power: 105 VA
 Length: 163-203 cm
 Headrest: 15 cm (extendable)
 Footrest: 25 cm (extendable)
 Legrest: 34 - 59 cm (length)
 Backrest: 52 cm (length)
 Padding width: 54 cm
 Seat depth: 50 cm
 Seat height: 56 - 100 cm
 Max. Working height: 139 cm
 Back adjustment: 80°
 Seat angle adjustment: 27°
 Adjustable foot: 90°
 Angle of rotation: 180°
 Weight: 85 kg
 Max. Seat load: 200 kg

Cushion-Contour padding



Classic-Contour padding



ORDER SCHEME

		P	U/D
Model	Pioneer 5M, base model	854	
Padding/motors = P	Classic-Contour-padding, 5 motors	2	
	Cushion-Contour-padding, 5 motors	3	
Padding colors = U/D	Classic-Contour-padding white		U16
	Cushion-Contour padding white		
	Padding Standard colors*		XXX*

* Please insert the corresponding color number at the end of the item number (U / D):
 Classic-Contour padding colors: U01 (cement), U13 (terre), U21 (prune), U14 (tomato), U23 (absinthe);
 Cushion-Contour padding colors: D1 (cement), D2 (terre), D3 (pearl). White is the standard color and do not require an additional specified color number

SÜDA foot care chair AF1



On the base of the chair are 4 stable and swivel castors. With their help, the AF1 can be easily changed or moved. The castors can then be fixed easily at the fixed parking position.

The 10 cm thick, breathable smooth padding with foam core is particularly hard-wearing and has a slightly grained surface. The upholstery is waterproof with a wet stop and is resistant to disinfecting and cleaning agents.

8503 Foot care chair AF1, color white

8503.XXX Foot care chair AF1, Standard colors*

* Please insert the corresponding color number at the end of the item number (U / D): U01 (cement), U13 (terre), U21 (prune), U14 (tomato), U23 (absinthe). White is the standard color and do not require an additional specified color number.



TECHNICAL SPECIFICATIONS

Electrical connection: 230 V / 50-60 Hz
 Power: 105 VA
 Length: 166-190 cm
 Width with armrests: 76 cm
 Seat height: 51-85 cm
 Weight: 60 kg
 Max. Seat load: 150 kg

Height adjustable to only 51 cm

Swivel armrests for a comfortable seating entrance

Electric adjustment of the seat height by foot switch

The two-part leg support can be individually adjusted, including the inclination angle

Extendable and removable headrest and footrest

Available in 6 upholstery colors



The economically priced starter model is a versatile functional space and scores with its space-saving compactness.

Thanks to its horizontally adjustable backrest, the AF1 can also be used for massage applications in addition to foot care treatments. The adjustable headrest ensures an optimal lying position for the customers. For massages, it is re-inserted and forms a comfortable nose slit.

Institute equipment

50 - 69

Cabinets

70 - 77

Work chairs

78 - 85

Magnifying lamps

86 - 89

Foot baths

Cabinets



You can now find the perfect complement for your foot care studio. Whether your focus is on the space-saving accommodation of a device, or you need additional storage space - we offer numerous options for perfect interior design.



Cubus foot care 400
> P. 54



Pedo cabinet mini
> P. 62



Pedo cabinet maxi
> P. 63



Cubus foot care 600
> P. 55



Secretario
> P. 64



Roll car Studio Line
> P. 65



Cubus foot care 800
> P. 56



Boby III
> P. 66



Practica II
> P. 66



Wellness cabinet
> P. 67

Cubus

Cubus combines design and quality made in Germany for your beauty salon, spa or pedicure studio.



Brackets for water bowl, magnifying lamp, waste bins and other accessories

DIN-tested, certified mineral deposit surface is extra-rugged and easy to clean, offers germs no liability option and is resistant to chemicals

The pedicure handpieces can be attached to the side, so that the hoses do not hinder access to UV compartment and drawers

Controlled cable management into technical compartment

Equipment compartment for SÜDA pedicure equipment with wet / dry technology

Steamed, smooth-running wheels, flexible in positioning and stable, lockable front castors

accessories > P. 60

In clinical, cosmetic and podiatry areas high demands are placed on hygiene. Therefore for the storage shelf, a material was developed, which provides bacteria, fungi or microorganisms no contact.



The integrated technical compartment of the CUBUS foot care series holds the equipment power cords. About a main switch on the rear wall of the CUBUS all devices can easily be turned on and off.



All-over UV compartment

To ensure a seamless all-round irradiation of the stored instruments in the UV compartment as well as the corresponding hygienic safety, our UV compartments are equipped with all-over UV compartments.



Full extension

The spacious drawers of the CUBUS series are provided with full extension and automatic soft-closing function and are particularly easy to clean by using a special plastic in the interior. The matching drawer inserts are available in the accessories section.

Cubus foot care 400

The small all-rounder for the pedicure cabin



Width 40 cm

accessories > S. 60

FEATURES

- 1184-13300 Cubus foot care 400 white, with UV compartment
- 1184-13301 Cubus foot care 400 light grey, with UV compartment

- 40 cm wide, DIN-tested, certified mineral-shelf, developed for working areas with the most stringent hygiene requirements
- Holders for magnifying lamp, pedicure handpiece and further accessories
- Equipment compartment, with space for a pedicure device of wet or dry technology
- Integrated technical compartment which borrows and hides the power cord of the pedicure devices
- Main switch on the rear wall of the cubus over which all devices can be turned on and off
- A spacious drawer, incl. full extension and automatic soft-closing function
- Integrated all-over UV compartment with complete all-round irradiation
- Damped, smooth running castors, front castors lockable
- Other colors and decors optional available



Cubus foot care 600

The compact model with plenty of storage space for your whole equipment



Width 60 cm

accessories > S. 60

FEATURES

- 1186-12300 Cubus foot care 600, white
- 1186-12301 Cubus foot care 600, light grey

- 60 cm wide, DIN-tested, certified mineral-shelf, developed for working areas with the most stringent hygiene requirements
- Holders for magnifying lamp, pedicure handpiece and further accessories
- Equipment compartment, with space for a pedicure device of wet or dry technology and further treatment accessories
- Integrated technical compartment which borrows and hides the power cord of the pedicure devices
- Main switch on the rear wall of the cubus over which all devices can be turned on and off
- Three spacious drawers, incl. full extension and automatic soft-closing function
- Suitable drawer inserts (optional)
- Damped smooth running castors, front castors lockable
- Other colors and decors optional available



Cubus foot care 800

The all-rounder with integrated UV compartment and outstanding space



FEATURES

- 80 cm wide, DIN-tested, certified mineral-shelf, developed for working areas with the most stringent hygiene requirements
- Holders for magnifying lamp, pedicure handpiece and further accessories
- Equipment compartment, with space for a pedicure device of wet or dry technology
- Integrated technical compartment which borrows and hides the power cord of the pedicure devices
- Main switch on the rear wall of the cubus over which all devices can be turned on and off
- Four spacious drawers, incl. full extension and automatic soft-closing function
- Suitable drawer inserts (optional)
- Damped smooth running castors, front castors lockable
- Other colors and decors optional available

1188-11230300 Cubus foot care 800, white, with UV compartment


1188-11230301 Cubus foot care 800, light grey, with UV compartment



Your individual Cubus

Together we will design your individual CUBUS MADE IN GERMANY according to your requirements and matching to the needs and colors of your institute.

Choose freely from our numerous modification options:

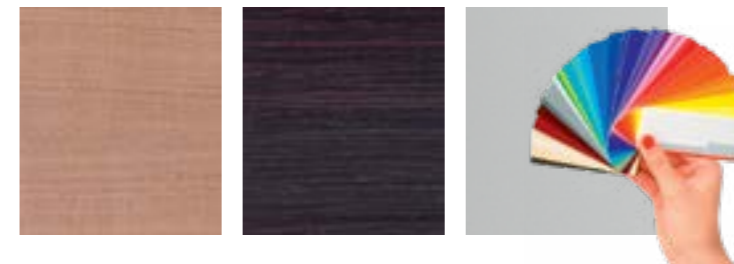

DECOR
Color design according to your requirements
Real wood veneer



STORAGE
Highboard with safety glass
equipment compartment



DRAWERS
with full extension and aluminum handle bar with Easy-Open for gentle and quiet closing



OPTIONAL
All-over UV compartment

Order scheme Cubus 400 / 600



The scope of delivery of the basic models includes:

- 1 body
- 2 drawers 24 cm high
- 2 handles
- 1 wooden cover plate
- 4 castors

1184	Cubus 400 Base model (W40 x H54 x D51 cm)
1186	Cubus 600 Base model (W60 x H54 x D51 cm)

CONFIGURATION OPTIONS	LOWER LEVEL	MEDIUM LEVEL	UPPER LEVEL	TYPE OF EXTRACT	BODY COLOR	1184-XXX* SURCHARGE	1186-XXX* SURCHARGE
Drawer 1 x 24 cm high	1					-	-
Drawer 2 x 12 cm high	2					surcharge	surcharge
Drawer 1 x 24 cm high		1				-	-
Drawer 2 x 12 cm high		2				surcharge	surcharge
UV-Drawer 1 x 24 cm high		3				surcharge	surcharge
Without equipment compartment and highboard			0			-	-
Drawer 1 x 24 cm high			1			surcharge	surcharge
Drawer 2 x 12 cm high			2			surcharge	surcharge
Equipment compartment and highboard			3			surcharge	surcharge
Handle				0		-	-
White					0	-	-
Light grey					1	surcharge	surcharge
Wenge decor					3	surcharge	surcharge
Wenge real wood veneer					6	surcharge	surcharge
Hickory real wood veneer					8	surcharge	surcharge

Please add the number code, which results from the equipment options to the item number in order to get the cubus in the configuration you have chosen.

If your Cubus is to be fitted with a suitable holder for Waldmann magnifying lamps, please attach a "W" to the order number. In the case of a later change of magnifying lamp, the holders can also be ordered separately.

Order scheme Cubus 800

The scope of delivery of the basic model includes:

- 1 body
- 4 drawers 24 cm high
- 4 handles
- 1 wooden cover plate
- 4 castors



1188	Cubus 800 Base model (W80 x H54 x D51 cm)
------	--

CONFIGURATION OPTIONS	LOWER LEVEL		MEDIUM LEVEL		UPPER LEVEL**		TYPE OF EXTRACT	BODY COLOR	1188-XXX* SURCHARGE
	left	right	left	right	left	right			
Drawer 1 x 24 cm high	1	1							-
Drawer 2 x 12 cm high	2	2							surcharge
Drawer 1 x 24 cm high			1	1					-
Drawer 2 x 12 cm high			2	2					surcharge
UV-Drawer 1 x 24 cm high			3	3					surcharge
Without equipment compartment and highboard					0	0			-
Drawer 1 x 24 cm high					1	1			surcharge
Drawer 2 x 12 cm high					2	2			surcharge
Equipment compartment and highboard					3	3			surcharge
Handle							0		-
White								0	-
Light grey								1	surcharge
Wenge decor								3	surcharge
Wenge real wood veneer								6	surcharge
Hickory real wood veneer								8	surcharge

* Please add the number code, which results from the equipment options to the item number in order to get the cubus in the configuration you have chosen. If your Cubus is to be fitted with a suitable holder for Waldmann magnifying lamps, please attach a "W" to the order number. In the case of a later change of magnifying lamp, the holders can also be ordered separately.

** In the upper level, the compartment and tray can only be selected for the entire width. However, a combination of the unit compartment on the one hand and drawers on the other side as well as different drawer heights are not possible.

Cubus Accessories

Does your CUBUS have an adapter sleeve for Waldmann magnifying lamps (Ø 16mm)? Please add a "W" at the end of the item number ordered to obtain an appropriate holder for your Waldmann adapter sleeves.



WASTE CONTAINER SET

- 1182.20 Waste container plastic
- 1182.21 Waste container stainless steel

Consisting of: Waste bin with tilting lid, Holder for waste, Adapter sleeve for magnifying lamp holder standard (For hanging in Magnifying Lamps holder Ø 13 mm)



FILING TRAY SET

- 1182.50 Filing tray set steel

Consisting of: Stainless steel tray, Cutlery tray 30 x 20 x 2.8 cm, Anti-slip mat, Stainless steel arm, Adapter sleeve for magnifying lamp holder Standard.



WATER BOWL SET

- 1182.30 Water bowl set, plastic

Consisting of: Water bowl (2.5 liter), Holder for drinking bowl, Adapter sleeve for magnifying lamp holder Standard



PREPARATIONS HOLDER SET

- 1182.40 Preparations holder set

Consisting of: Preparations holder 3 x Ø 9 cm, Holder for Preparations, Adapter sleeve for magnifying lamp holder Standard



ROD HOCK HOLDER UNIVERSAL

- 1182.15 Rod holder universal for accessories sets



DRAWER INSERT

- 1182.61 Cutlery tray and sliding upper shell, suitable for CUBUS 600

Drawer Division for deep thrusts 54.7 x 43 cm



DRAWER INSERT

- 1182.60 Cutlery tray and sliding front cover, suitable for CUBUS 400 and 800

Drawer subdivision for deep thrusts 34.7 x 43 cm



DRAWER INSERT

- 1182.63 2 x insert shell and sliding front cover, (transparent) for CUBUS 400 and 800

Drawer subdivision for deep thrusts 34.7 x 27.5 cm



HANPIECE HOLDER

- 1182.10 Handpiece holder M
For Medivac, Pedo Sprint 50, Sprint Eco, SÜDA vac S

- 1182.11 Handpiece holder D
For Pedo Sprint 100-D1

- 1182.13 Handpiece holder F
For für Pedo FS 3 und Air AS 1



HANPIECE HOLDER

- 1182.12 Handpiece holder
Universal for all wet technology devices



UV-COMPARTMENT

- 30132 UV-tube 8 Watt
- 1182.65 UV-compartment inner glass plate for Cubus 400 / 800
- 1182.66 UV-compartment inner glass plate for Cubus 600
- 32876 Filling tray for UV compartment

Pedo cabinet mini



The narrow version of the Maxi-cabinet offers all the benefits at only 43 centimeters wide enclosures. With particularly hard-wearing mineral worktop.



Width 43 cm

MODELS

8321	Pedo cabinet mini without UV compartment, light grey
8321.2	Pedo cabinet mini with UV compartment, light grey

FEATURES

- Dimensions: W 43 x H 80 x D 52 cm
- 2 drawers + 1 device tray slide or
- 1 drawer + 1 device tray slide + 1 UV compartment
- 1 magnifying lamp holder in the cover plate
- Cover plate pivotable with a few hand movement
- Large handle on the sides for pushing or can be used as a towel holder
- Easy care and moisture resistant
- 4 lockable castors
- 1 internal 3-way socket

ACCESSORIES

1175.002	Lupenleuchtenhalter, Aufnahme 16 mm für Waldmann-Leuchten
----------	---

Pedo cabinet maxi



With device tray slide and drawers as well as UV compartment. With particularly hard-wearing mineral worktop.



Width 81 cm

MODELLE

8320	Pedo cabinet maxi without UV compartment, light grey
8320.1	Pedo cabinet maxi 2 units extract trays without UV compartment, light grey
8320.6	Pedo cabinet maxi with UV compartment, light grey

FEATURES

- Size: W 81 x H 80 x D 52 cm
- Unit extract tray and drawers or
- Unit extract tray, drawers and UV compartment
- 2 Magnifying Lamps holder incorporated in the cover plate (right + left)
- Cover plate pivotable with a few hand movement
- Easy care and moisture resistant
- 4 lockable castors
- 1 internal 3-way socket

ACCESSORIES

1175.002	Magnifying lamps mounting sleeve for Waldmann lamps
----------	---

Secretario



MODELS

1170	Secretario without UV compartment, smooth white
1170.1	Secretario with UV compartment, smooth white



Width 58 cm

ACCESSORIES

8206.1	Magnifying Lamps wall and cabinet mounting
8206.1W	Magnifying Lamps wall and cabinet mounting for Waldmann lamps

FEATURES

- Size: W 58 × H 74 × D 42 cm
- Wide motors extraction tray
- 3 drawers or
- 2 drawers + 1 UV compartment
- Cover plate secured by raised side walls and the rear by a stainless steel railing
- 4 lockable castors
- Internal 3-way socket

HANDPIECE HOLDER

35527	Handpiece holder for Pedo Sprint ECO/50/100
1505.013	Handpiece holder for Pedospray NT 40/NT 40 light for self mounting on the cabinet



Roll car Studio Line



The exclusive furniture with tight radii adapts shapely in each cabin. With a large capacity extract and 4 drawers all accessories find their place. The comfortable soft feeders ensure that when closing nothing overturn or slips.



Upon request we provide you with a custom made roll car STUDIO LINE in your favorite color.

MODELS

1227	Roll car STUDIO LINE, white matt lacquer
1227.1	Roll car STUDIO LINE, light gray matt lacquer

FEATURES

- Dimensions: W 54.5 × H 86 × D 41 cm
- Cabinet with glass cover plate
- 4 large, smooth-running castors
- 1 × push large, bottom with handle shells
- 4 × push small, top with handle shells
- 2 × magnifying lamp holders standard
- Power supply in the technical compartment

ACCESSORIES

1175.005W	Magnifying lamp holder, 16 mm for Waldmann lamps
-----------	--

Boby III

Space miracle on 5 wheels:
The design classic made by Joe Colombo

8218	BOBY III white
8218.5	BOBY III silver



FEATURES

- Dimensions: W 41 × H 74 × D 43 cm
- Completely Made of plastic
- 2 rotary drawers



Practica II

Multi-purpose cart on 4 easy-running,
lockable castors

8317.2	PRACTICA II
--------	-------------



FEATURES

- Dimensions: W 57 × H 91 × D 42.5 cm
- 4 spacious drawers
- Compartment on the top can be pulled out to the front
- 2 storage boxes on the sides
- Sturdy and easy to clean plastic surfaces
- Color: white

Wellness cabinet

The elegant trolley in fine wood decor, with up to three
drawers, has plenty of space and storage space for your
working materials in the Institute/practice.

MADE IN
GER 2
WARRANTY

8322.5	Wellness cabinet, white
8322.9	Wellness cabinet, dark-faux wenge optic



ACCESSORIES

8206.1	Magnifying Lamps wall and cabinet mounting
8206.1W	Magnifying Lamps wall and cabinet mounting for Waldmann lamps
8322.51	Drawer, white
8322.91	Drawer, dark-faux wenge optic



FEATURES

- Dimensions: W 66.3 × H 80 × D 41.4 cm
- Devices and preparations cart on easy-rolling,
lockable castors for easy relocation
- Based classification with 3 deep shelves
- Without drawers, please order separately
- Sturdy and easy-care surfaces
- Pilot holes for a magnifying lamps mounting on
left or right already exists, the mounting must be
ordered separately

Cabinet



FEATURES

- Equipment and preparation trolleys on easy-running, lockable castors for convenient relocation
- Basic classification with 3 deep storage areas
- White lacquered metal frame
- Drawers are available optionally
- Robust and clean surfaces
- Dimensions: B 66.3 × H 80 × D 41.4 cm

ACCESSORIES

8406.6 Magnifying lamp holder, 13 mm

8406.6W Magnifying lamp holder, 16 mm for Waldmann lamps

8316.003 Bracket for water bowl holder

8322 Cabinet A, without drawer, White



8322.004 Cabinet B, with drawer, White



8322.005 Cabinet C, with 2 drawers, White

one4all

A simple and still comfortable cabinet for storing your treatment accessories.

This multitalent is suitable for foot care and cosmetic cabin. With an equipment opening extract and 3 drawers, the small, space-saving cabinet offers sufficient storage space for your treatment accessories.



FEATURES

- Dimensions: W 40 × H 76 × D 40 cm
- 3 drawers, 1 equipment opening extract
- 1 Magnifying lamp holder in the cover plate
- Large handle on both sides for pulling or as a towel holder
- 3-way socket inside

8315.10 Equipment cabinet one4all

Work chairs



A work chair must be perfectly matched to the users and their workplace. Especially in activities that require a high degree of concentration, it is crucial that chiropodists and podiatrists may practice in a pleasant position. SÜDA has the ideal work chair for each claim in the range of products.

2

WARRANTY

Combination diversity

Diversity of colors



Combi-chair program
> P. 72



Work stool Rocky
> P. 75



Work stool Syria
> P. 75



Swinging stool
> P. 76



Work stool Rodeo
> P. 76



Work stool FIT
> P. 77



00 (white)



U01 (cement)



U13 (terre)



U21 (prune)



U14 (tomato)



U23 (absinthe)

The combi-chair program

The combi-chair program allows you to put together your customized work chair. Here you can choose between different seat and backrest forms, gas springs, foot cross and castors respectively. For the upholstery color choose from our 6 standard colors. If you have a specific color request, please contact us.



ACCESSOIRES

8303.10	Terry-cloth cover for round chair, white
8303.11	Terry-cloth cover for saddle chair, white
8303.12	Terry-cloth cover for diamond rest, white
8303.13	Terry-cloth cover for banana-shaped rest, white
8303.14	Terry-cloth cover for trapeze rest, white
8303.15	Terry-cloth cover for shell chair, white

ORDER SCHEME	ITEM	L	G	F	R	U
Base model	Round chair 360	8002-	-	-	-	-
	Round chair	8003-				
	Saddle chair	8004-				
	Shell chair	8007-				
Backrest = L	No rest		0			
	Banana-shaped rest, incl. attachment		1			
	Diamond rest, incl. attachment		2			
	Trapeze rest		3			
Pneumatic spring = G	Standard pneumatic spring in chromium look		0			
	Long pneumatic spring in Chromium look		1			
Foot cross = F	Plastic foot cross			0		
	Foot cross in chromium look			1		
Castors = R	Plastic castors, 5 pieces				0	
	Castors in chromium look, 5 pieces				2	
Color = U	white					00
	Standard upholstery color*					XXX*

* Add color number at the end of the item number by: U01 (cement), U13 (terre), U21 (prune), U14 (tomato), U23 (absinthe)



BACKREST

Banana-shaped backrest
The banana-shaped backrest for supporting the arms can be combined with all chair models.



Diamond backrest
The diamond backrest is available for round and saddle chair only.



Trapeze backrest
The trapeze backrest can be combined with all chair models.

UPHOLSTERY



Round padding
Round seat cushion, seat diameter: 40 cm



Saddle padding
Ergonomically shaped seat seat, seat diameter: 40 cm



Shell padding
Ergonomically shaped molded wood shell, seat diameter: 43 cm

GAS SPRING



Standard gas spring in chrome look
min. seat height: 45.5 cm / max. seat height: 58 cm



Standard gas spring in chrome look
min. seat height: 45.5 cm / max. seat height: 58 cm



FOOT CROSS

Foot cross plastic
Ø 64,5 cm



Foot cross chrome look
Ø 67 cm



CASTORS

Castors plastic, for hard floors



Castors chrome look, for hard floors

Combinations

Round chair 360

The banana back of round chair 360 can be rotated independently from the seat. For longer or difficult work in a position to prevent the fatigue of the arm.



TECHNICAL SPECIFICATIONS

Seat diameter: 40 cm
Seat height, distance to the seat surface and backrest height adjustable

OPTIONS

- Standard gas spring or long gas spring
- Foot cross in plastic or in chrome look
- Castors on plastic or chrome look

Round chair

The flexible back support with ergonomic diamond rest remains in all movements in close contact with the back. Your back is always supported.

TECHNICAL SPECIFICATIONS

Seat diameter: 40 cm
Seat height, distance to the seat surface and backrest height adjustable



OPTIONS

- Banana-shaped backrest for supporting the arms, diamond backrest, trapeze backrest or without backrest
- Standard gas spring or long gas spring
- Foot cross in plastic or in chrome look
- Castors in plastic or chrome look

Saddle chair

Firmly in the saddle, so you sit on the saddle chair, with both feet firmly standing on the floor.



TECHNICAL SPECIFICATIONS


Seat diameter: 40 cm
Seat height, distance to the seat surface and backrest height adjustable

OPTIONS

- Banana-shaped backrest, diamond backrest, trapeze backrest or without backrest
- Standard gas spring or long gas spring
- Foot cross in plastic or in chrome look
- Castors in plastic or chrome look

Shell chair

The unique ergonomics of the chair provides you the perfect seating position, by giving your body the ideal stability.

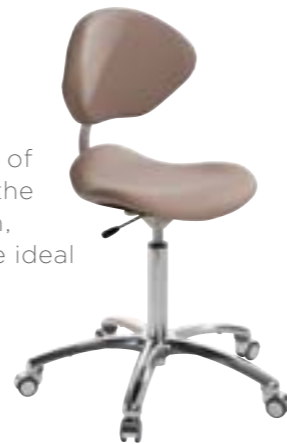

TECHNICAL SPECIFICATIONS

Seat diameter: 43 cm
Seat height, distance to the seat surface and height of the backrest adjustable

OPTIONS

- Banana-shaped backrest, diamond backrest, trapeze backrest or without backrest
- Standard gas spring or long gas spring
- Foot cross in plastic or in chrome look
- Castors in plastic or chrome look

Work stool Rocky

The saddle shaped supports the upright sitting posture and promotes blood circulation in the legs.



1268 Rocky, white

1268-XXX Rocky, colored*

* Add color number at the end of the item number by: U01 (cement), U13 (terre), U21 (prune), U14 (tomato), U23 (absinthe)

FEATURES

- Saddle cover made of easy-care, breathable, water steam-permeable leather
- Rotatable, ergonomic saddle seat
- Pneumatic spring „Lift-o-Mat“, protects spine and back
- High stability due to 5-foot bogie wheels with brakes
- Sturdy and yet light aluminum frame
- Available in 6 colors
- Seat circumference: 40 cm
- Seat height: 41 to 53 cm
- Weight: 5 kg

Work stool Syria

The backrest can serve as an armrest and adapts flexibly to your favorite sitting position in difficult or prolonged activities.



1301 Syria, white

1301-XXX Syria, colored*

* Add color number at the end of the item number by: U01 (cement), U13 (terre), U21 (prune), U14 (tomato), U23 (absinthe)

FEATURES

- Seat cover made from easy-care, breathable synthetic leather
- 360° rotatable back and armrests, separately rotating seat
- Safe stand by 5-foot bogie mounting with castors
- Gas spring in chrome look
- Pressure ring for height adjustment below on the underside of the seat cushion
- Available in 6 colors
- Seat circumference: 40 cm
- Seat height: 45 - 56 cm
- Weight: 10 kg

Swinging stool

360° free moving space, fully extended usable as a support aid.



8202.6 Swinging stool

FEATURES

- Ergonomically shaped seat
- Continuously height-adjustable
- Pendulum base foot made of beech, lacquered in stainless steel
- Anti-slip floor coating

Seat area: B 42 x D 32 cm
 Seat height: 53 - 79 cm
 Color white

Work stool Rodeo

This work stool gives you more flexibility in your daily work. Through his arm-holes on the back, you can during treatment, where a "steady hand" is needed, place the arms in the recesses. Use the backrest as work surface.



8202.5 Work stool Rodeo, white

8202.5-XXX Work stool Rodeo, colored*

Add color number at the end of the item number by: U01 (cement), U13 (terre), U21 (prune), U14 (tomato), U23 (absinthe)

FEATURES

- Back rest with practical arm recesses for treatment support
- Backrest tilt and height adjustable variable
- Adjustable height
- 5-foot-bogie safety braked castors
- Available in 6 colors

Seat surface: 42 cm
 Seat height: 41 - 50 cm
 Weight: 10 kg

Work stool FIT

For this work chair everything was designed for your convenience. While you can easily place your feet on the height-adjustable foot ring, the extra-large backrest with comfortable padding provides you everytime a comfortable seating position.



infinitely adjustable tilt angle

Height adjustable foot ring

FEATURES

- Ergonomically shaped seat and backrest upholstery
- Tilt angle and height of the seat adjustable
- 5-foot frame with safety braked castors
- Comfortable, height-adjustable foot ring for placing the feet
- Available in 6 colors

Seat surface: D 46 cm x W 44 cm
 Seat height: 41 to 53 cm
 Weight: approx. 11 kg

8202.4 Work stool FIT, white

8202.4-XXX Work stool FIT, colored*

* Add color number at the end of the item number by: U01 (cement), U13 (terre), U21 (prune), U14 (tomato), U23 (absinthe)

ACCESSOIRES

8303.38 Terry cloth, 2 pieces, white



white



U01 (cement)



U13 (terre)



U21 (prune)



U14 (tomato)



U23 (absinthe)

Magnifying lamps



Illuminate your workspace with our magnifying lamps for comfortable work even details. Select between the innovative LED models or you can rely on a donut concept and opt for the ergonomic features.

Magnifying lamp Basic LED



Magnifying lamp Standard LED



 5 years warranty

EQUIPMENT FEATURES

- Lightweight lamp with spring arm
- Metal, white lacquered
- Circular LED module
- Light emitting diodes with 50,000 hours of life expectancy
- 3-step dimming
- 13 mm spigot for insertion into magnifying lamp mounting

TECHNICAL SPECIFICATIONS

Electrical connection: 230 V/50 Hz
 Power: 9,5 watts
 Lens thickness: 3.0 diopters
 Magnifying glass diameter: 12.7 cm
 Magnifying head diameter: 16.8 cm
 Length jointed-arm: 100 cm
 Total projection: 123 cm
 Weight: 1,7 kg

8206.50 Magnifying lamp Basic LED

EQUIPMENT FEATURES

- Magnifying lamp with parallel triple arm
- Metal, white lacquered
- Two half-round LED modules (dimnable)
- Light emitting diodes with 50,000 hours of life expectancy
- Timer for automatic shutdown after 4 or 9 hours
- 13 mm spigot for insertion into magnifying lamp mounting

TECHNICAL SPECIFICATIONS

Electrical connection: 230 V/50 Hz
 Power: 11 watts
 Lens thickness: 3.0 diopters
 Magnifying glass diameter: 12.7 cm
 Magnifying head diameter: 19 cm
 Length of articulated arm: 105 cm
 Total projection: 138 cm
 Weight: 3 kg

8206.10 Magnifying lamp Standard LED

Magnifying lamp Magic / Magic LED

Available in two versions:
With ring fluorescent tube or with
LED ring lamp concept



EQUIPMENT FEATURES

- Magnifying glass strength 3.0 diopters each
- Provides comfortable warm-white light
- Very elegant, light magnifying lamp head
- Diameter of ring light tube only 15 mm
- LED: high light output by maximum energy efficiency
- 13 mm pin for insertion into magnifying lamps fixtures

TECHNICAL SPECIFICATIONS

Electrical connection: 230 - 240 V / 50 Hz
Power LED MAGIC: 56-pieces LED light 5 W
Power MAGIC: 22 W fluorescent light
Magnifying glass strength: 3 diopters
Diameter of magnifying glass: 12,7 cm
Length of lamp head: 26,0 cm
Length of lamp arm: 90 cm
Weight: approx. 2,5 kg

8206.34 Magnifying lamp Magic

8206.35 Magnifying lamp Magic LED

Magnifying lamp Standard



 5 years warranty

EQUIPMENT FEATURES

- Lamp with covered spring joint arm
- Metal, painted white
- Glass lens
- Visor
- 13 mm pin for insertion into magnifying lamps fixtures

TECHNICAL SPECIFICATIONS

Electrical connection: 230 V / 50 Hz
Power: 22 Watt
Magnifying glass strength: 3,0 diopters
Diameter of magnifying glass: 12,5 cm
Diameter of magnifying head: 26,5 cm
Length of lamp arm: approx. 104 cm
Weight: approx. 4,3 kg

8206 Magnifying lamp Standard

Magnifying lamp iQ LED

The magnifying lamp LED iQ with its modern and ergonomic design features has a crystal-clear lens with 1.75 times magnification and 60 ultrabright daylight LEDs for accurate detail in foot care.

The 60 LEDs of iQ have a life expectancy of up to 20,000 hours, which corresponds to a lifetime of over 20 years. An exchange of the lamps is not required. The iQ LED Magnifying Lamp embodies the ideal support for optimal working in details .



EQUIPMENT FEATURES

- 4-steps dimming for optimum light intensity
- Flexible joint between head and spring balanced arm allows effortless movement with a long range and a precise and stable positioning
- Optional 5-diopter lens
- "Quick Swap Lens" system allows easy replacement of the lens
- Luminous intensity of about 6,000 lux for area illumination with fine detail

TECHNICAL SPECIFICATIONS

Voltage: 100 - 240V / 50 - 60Hz
Power: 12 Watt
Magnifying glass strength: 3/5
Diameter of magnifying glass: 17.5 cm
Diameter of magnifying head: 23 cm
Length of lamp arm: 118 cm
Weight: 3,6 kg

8206.40 Magnifying lamp iQ LED

ACCESSOIRES

8206.401 5-diopter lens

Magnifying lamp Circle / Circle LED

Unusual design combined with optimum functionality - the ring magnifying lamp Circle is not just a piece of jewelry at work, it will open up a new quality when looking at your work area.



This cold light magnifying lamp can be changed over to a daylight lamp by exchanging the light bulb with the daylight tube (8407.1).

 **5 years warranty**

EQUIPMENT FEATURES

- Cold light ring tube with low heat development and flicker-free light or
- Round LED module, dimmable 10-100%
- Precision glass lens
- 13 mm spigot for insertion into magnifying lamp mounting
- Delivery including table clamp

TECHNICAL SPECIFICATIONS

Electrical connection: 230 V / 50 - 60 Hz
 Power: 24 Watt cold light
 Magnifying glass strength: 3,5 diopters
 Diameter of magnifying glass: 16,5 cm
 Diameter of magnifying head: 27,0 cm
 Length of lamp arm: approx. 104 cm
 Weight: approx. 3,6 kg

8206.71 Magnifying lamp Circle,
white / grey



TECHNICAL SPECIFICATIONS

Electrical connection: 100-240 V / 50-60 Hz
 Power: 10 Watt (24 LED)
 Magnifying power: 3.5 diopters
 Diameter of magnifying glass: 16,5 cm
 Diameter of magnifying head 24.5 cm
 Length of lamp arm: 100 cm
 Weight: 4,9 kg

8206.72 Magnifying lamp Circle LED
white / grey

Magnifying lamp head made of aluminum

Safety cover

Rugged aluminum housing

High mobility thanks to newly designed lamp arm

Stable technology for continuous use in professional use

2 handles for easy attachment of the light head to the right working position, suitable for right and left handers.

Magnifying lamp Tevisio

LED lighting technology for high performance. The new TEVISIO Waldmann offers maximum work output illumination paired with an exceptional design. Equipped with cutting edge LED technology it features not only the highest light output by its unique flexibility.

The innovative hinge technology allows perfect positioning in each treatment situation. Convince yourself of functionality and aesthetics of the new TEVISIO Waldmann, of course, MADE IN GERMANY.



Light intensity dimmable to 10% of maximum strength

Hygienic on/off button

Portable, lockable joints

Scratch resistant glass lens

EQUIPMENT FEATURES

- Scratch resistant glass lens with 3.5 diopters
- LED lighting technology
- Visor
- 16 mm pin for insertion into magnifying lamp mounting

TECHNICAL SPECIFICATIONS

Electrical connection: Wall wart 100-240 V / 50-60 Hz
 Power: 14 W
 Magnifying glass strength: 3,5 Dioptrien
 Diameter of magnifying glass: 16,5 cm
 Diameter of magnifying head: 26,5 cm
 Length of lamp arm: approx. 99 cm
 Luminous intensity: 6.000 Lux

1143 Magnifying lamp Tevisio



Magnifying lamp accessoires



HALOGEN LAMP

8915.001 Halogen lamp 20 Watt

For suitcase light 1117.1



COLD LIGHT TUBE 9 WATT

30466 Cold light tube 9 Watt

For magnifying lamp VISTA 1141 / SNT.



REPLACEMENT RING TUBE MAGNIFYING LAMPS LUXO AND WALDMANN

8407.20 Replacement ring tube 22 Watt

For ring magnifying lamp LUXO and Waldmann



REPLACEMENT RING TUBE MAGNIFYING LAMP CIRCUS

8206.700 Replacement ring tube 22 Watt

Cold light tube for magnifying lamp CIRCUS



REPLACEMENT RING TUBE MAGNIFYING LAMP WALDMANN LRE 122

1142.003 Replacement ring tube T R16

For magnifying lamp Waldmann LRE 122



TABLE MOUNT

8406.3 Table mounting fixture 13 mm

For all magnifying lamps except Waldmann



TABLE MOUNT

8206.21 Table mounting fixture 13 mm

For all magnifying lamps except Waldmann



TABLE MOUNT

8206.2 Table mounting fixture 13 mm

For all magnifying lamps VISTA and LUXO



HOLDER FOR WATER BOWL HOLDER

8316.003 Holder for water bowl holder for cabinets A / B / C

MAGNIFYING LAMPS HOLDER CABINET MOUNTING

8406.6 Mounting 13 mm

8406.6W Mounting 16 mm for Waldmann

MAGNIFYING LAMP WALL / CABINET MOUNTING

8206.1 Mounting 13 mm

8206.1W Mounting 16 mm for Waldmann

ADAPTER SLEEVE (POM) FOR MAGNIFYING LAMP HOLDER

1182.001 Mounting 13 mm

1182.001W Mounting 16 mm for Waldmann

MAGNIFYING LAMP HOLDER CUBUS CABINETS

1175.003 Mounting 13 mm

1175.003W Mounting 16 mm for Waldmann

MAGNIFYING LAMP HOLDER UNOVAG / PODOAIRSPEED

1175.005 Mounting 13 mm

1175.005W Mounting 16 mm for Waldmann

MAGNIFYING LAMP HOLDER SLEEVE WALDMANN-LAMPS (16 MM)

1175.004 Mounting 13 mm

1175.002 Mounting 16 mm for Waldmann

For PEDO cabinets with integrated magnifying lamp holder in the cover plate

STATIV

8406.1 Tripod for magnifying lamps

5-foot tripod. With static weights. ø 60 cm. Adjustable height. Weight: 4,3 kg. For all magnifying lamps except Waldmann. Mounting ø 13 mm for the insertion of magnifying lamps (except Waldmann).

Foot baths



Foot bath Hygienic
> P. 88



Mobile Foot bath
> P. 89



Foot bath foursquare
> P. 89

2
WARRANTY

Modern foot baths allow you maximum hygiene and safety, even at rapid customer acceptance. Due to exchangeable plastic inserts the spread of bacteria and fungi is stopped consistently.

Foot bath Hygienic



This modern foot bath allows you maximum hygiene and safety, even at rapid customer acceptance. By interchangeable plastic inserts, you spared cleaning the bath. Before use, place one of the accurately fitting plastic inserts in the tray and fill it with water. After the foot bath, drain the tray and discard the plastic insert. The propagation of bacteria, fungi and the infections is hereby prevented.



TECHNICAL SPECIFICATIONS

Dimensions: W 30 x H 11,5 x L 39 cm
 Filling capacity: 5 liters
 Weight: 420 g

8153.40 Foot bath Hygienic

ACCESSOIRES

8153.41 Plastic inserts, 80 pieces

Mobile Foot bath



stable, white lacquered tube stell frame
 4 smooth-running castors with brakes
 Clipboard grating

TECHNICAL SPECIFICATIONS

Dimensions: W 32 x H 35 x L 45 cm
 Tank capacity: 10 Liter
 Optical filling capacity: 5 liters

8153.3 Mobile foot bath

ACCESSOIRES

8153.31 Foot bath bags, 100 pieces

Foot bath Foursquare

Comfortable handle bar is used as a towel holder, while allowing a displacement of the tray, without having to bend down. If necessary, the hanger is simply to dismantle by loosening one screw.

4 smooth-running castors with brakes
 Handle bracket and towel holder



TECHNICAL SPECIFICATIONS

Dimensions: W40 x H 15 x T 40 cm
 Tank capacity: 15 Liter
 Optical filling capacity: approximately 8 liters

8153.6 Foot bath foursquare

ACCESSOIRES

8153.61 Foot bath bags, 100 pieces

Treatment systems

97 - 109

Studio Line

100

Cabinet systems with wet and suction technology

104

Cabinet systems with compressed air technology

111 - 128

Flexi Line

112

Foot care equipment with suction technology

118

Foot care equipment with wet technology

129

Accessories

Studio Line cabinet systems

MADE IN
GER

2
WARRANTY



Spray Vac Unit
>P. 102



podoAIRspeed
>P. 106

CATEGORY	Wet technology	Suction technology	Compressed air technology
HANDPIECE			
SPEED (RPM/ MIN)	6.000 - 40.000	6.000 - 40.000	Handpiece = max. 27.000 Turbine = max. 400.000
LEFT/RIGHT-HANDED ROTATION	Yes	Yes	Yes
DIAMETER	20 mm	23,5 mm	Turbine = 16,4 mm Handpiece = 20 mm
WEIGHT	163 g	112 g	Turbine = 65 g Handpiece = 207 g
LENGTH	140 mm	134 mm	Turbine = 120 mm Handpiece = 165 mm
LIGHT FUNCTION	Yes	No	Optional
CONTROL UNIT / CABINET			
MAX SUCTION POWER	-	210 l/min	-
NOISE LEVEL			
AT MIN. SUCTION OUTPUT	-	45 db (A)	-
AT MAX. SUCTION OUTPUT	-	49 db(A)	-
DIGITAL DISPLAY	Yes		No
MEMORY FUNCTION	3 memory locations		No
FOOT SWITCH	Optional		Yes
THREE FUNCTIONAL NOZZLE	Yes		Yes
UV COMPARTMENT	Yes		Yes
WIDTH	100 cm		100 cm
DEPTH	48 cm		48 cm
HEIGHT	81 cm		81 cm
WEIGHT	80 kg		76 kg
ELECTRICAL CONNECTION	230 V / 50 - 60 Hz		230 V / 50 - 60 Hz

Flexi Line – Suction technology

MADE IN
GER

2
WARRANTY



Pedo Sprint
ECO SP
>P. 113



Pedo Sprint
50 SP
>P. 114



Pedo Sprint
100 SP
>P. 115



SÜDA vac S
>P. 116

CATEGORY	Suction technology	Suction technology	Suction technology	Suction technology
HANDPIECE				
SPEED (RPM/ MIN)	6.000 - 16.000	6.000 - 25.000	6.000 - 40.000	6.000 - 40.000
LEFT/RIGHT-HANDED ROTATION	No	Yes	Yes	Yes
DIAMETER	23,5 mm	23,5 mm	23,5 mm	23,5 mm
WEIGHT	119 g	119 g	112 g	112 g
LENGTH	134 mm	134 mm	134 mm	134 mm
LIGHT FUNCTION	No	No	No	No
CONTROL UNIT				
MAX SUCTION POWER	180 l/min	180 l/min	180 l/min	210 l/min
NOISE LEVEL				
AT MIN. SUCTION OUTPUT	53 db(A)	53 db(A)	53 db(A)	45 db(A)
AT MAX. SUCTION OUTPUT	60 db(A)	60 db(A)	60 db(A)	49 db(A)
DIGITAL DISPLAY	No	Partially	Yes	Yes
MEMORY FUNCTION	No	No	3 memory locations	3 memory locations
OPTIONAL FOOT SWITCH	No	Yes	Yes	Yes
WIDTH	220 mm	220 mm	220 mm	240 mm
DEPTH	180 mm	180 mm	180 mm	226 mm
HEIGHT	140 mm	140 mm	140 mm	126 mm
WEIGHT	3 kg	3 kg	3 kg	3,6 kg
ELECTRICAL CONNECTION	220-240 V/ 50 - 60 Hz	220-240 V/ 50 - 60 Hz	220-240 V/ 50 - 60 Hz	110-240 V/ 50-60 Hz

Flexi Line – Wet technology



Asteria 1000 > P. 119
 Pedospray NT 40 > P. 120
 Pedospray NT 40 Light > P. 121
 Südatronic Komfort > P. 122

CATEGORY	Wet technology	Wet technology	Wet technology	Wet technology
----------	----------------	----------------	----------------	----------------

HANDPIECE

SPEED (RPM/ MIN)	6.000 – 30.000	5.000 – 40.000	8.000 – 40.000	5.000 – 40.000
LEFT/RIGHT-HANDED ROTATION	Yes	Yes	Yes	Yes
DIAMETER	22 mm	22 mm	22 mm	22 mm
WEIGHT	125 g	120 g	115 g	148 g
LENGTH	150 mm	135 mm	135 mm	149 mm
LIGHT FUNCTION	No	No	Yes	No

CONTROL UNIT

DIGITAL DISPLAY	Partially	Yes	Yes	Yes
MEMORY FUNCTION	No	3 memory locations	3 memory locations	3 memory locations
OPTIONAL FOOT SWITCH	Yes	Yes	Yes	Yes
WIDTH	220 mm	250 mm	250 mm	269 mm
DEPTH	180 mm	250 mm	250 mm	195 mm
HEIGHT	140 mm	80 mm	80 mm	135 mm
WEIGHT	3,2 kg	3,2 kg	3,2 kg	4 kg
ELECTRICAL CONNECTION	230 V/50 - 60 Hz	230 V/50 - 60 Hz	230 V/50 - 60 Hz	230 V/50 Hz



SÜDA spray S > P. 126
 SÜDA spray L > P. 127
 Südatronic Starlight > P. 123
 Südatronic Starlight Power > P. 124

CATEGORY	Wet technology	Wet technology	Wet technology	Wet technology
----------	----------------	----------------	----------------	----------------

HANDPIECE

SPEED (RPM/ MIN)	6.000 – 40.000	6.000 – 40.000	5.000 – 40.000	8.000 – 40.000
LEFT/RIGHT-HANDED ROTATION	Yes	Yes	Yes	Yes
DIAMETER	20 mm	20 mm	22 mm	22 mm
WEIGHT	163 g	163 g	151 g	129 g
LENGTH	140 mm	140 mm	149 mm	137 mm
LIGHT FUNCTION	No	Yes	Yes	Yes

CONTROL UNIT

DIGITAL DISPLAY	Yes	Yes	Yes	Yes
MEMORY FUNCTION	3 memory locations	2 memory locations	3 memory locations	3 memory locations
OPTIONAL FOOT SWITCH	Yes	Yes	Yes	Yes
WIDTH	240 mm	240 mm	269 mm	269 mm
DEPTH	225 mm	225 mm	195 mm	195 mm
HEIGHT	122 mm	122 mm	135 mm	135 mm
WEIGHT	2,8 kg	2,8 kg	4 kg	4 kg
ELECTRICAL CONNECTION	24 V/50 - 60 Hz	24 V/50 - 60 Hz	230 V/50 Hz	230 V/50 Hz



Studio Line

Compact treatment units for pedicure and podology

The high-end cabinet systems of the studio line from SÜDA combine high-performance technology with an optimal storage space offer. The special design makes it easier to comply with essential hygiene standards.

Thanks to a special, constructive implementation, all units are characterized by an average 25% lower noise level compared to conventional treatment devices, which makes the work much more comfortable for the customer as well as for you as an operator.

Studio Line

The perfect solution for your equipment



Highest hygiene standards

Only during the treatment, the side flaps are open and the handpieces are readily available. After cleaning they are protected against pollution and contamination by closing the doors. All the accessories are stored as well, nothing is open. More Hygiene: the easily cleaned and disinfected foil display.

Incredibly quiet

During the design stage of StudioLine much emphasis was placed on minimizing noise. So the system technology works within the treatment units almost noiseless. So the burden of operating noise could be reduced by 20 to 25% compared to conventional foot care devices.

Compact combination of treatment system + storage space

Thanks to ergonomically placed handpieces in the side flaps, there is plenty of storage space for preparations and working utensils in the drawers, compartments and extracts.



High quality durable work place with integrated brackets for magnifying lamps and accessories.



Drawers with plenty of storage space



Optional with UV compartment



Design in elegantly rounded lines



Wet and suction technology

Studio Line



In the new Studio Line cabinet system, wet and suction technology is used in combination and offers the greatest possible flexibility for the user. Benefit from the advantages of a compact treatment unit with perfectly adjustable performance parameters. Optionally, two special attachments are available for the wet handpiece, which also make your work easier.

Special essay handpieces

Omnicut

Suitable for micromotors NBX, NBX LED, NLX Nano and SÜDA MM-WS/ WS L

The OMNICUT attachment is a useful addition to the highquality water technology appliances. Simply remove the hand piece from the Pododont micromotor and attach OMNICUT; the rotating mini blades remove even stubborn, thick calluses in a matter of seconds, without leaving any skin dust behind! The blades are replaceable and can be sterilised. Optimal cutting quality is guaranteed even after several treatments.



1503.8 OMNICUT-attachment with 2 Blades

1503.800 OMNICUT-Blades, 20 pieces (10 set)

SÜDA transmission handpiece

Suitable for micromotors NBX, NBX LED, NLX Nano and SÜDA MM-WS/WS L

This extremely handy and lightweight hand piece reaches up to 200,000 rpm. This high rotational speed combined with small diamond cutters enables exact, precise and vibration-free operation on corns, calluses, nail fold area and nails. This leads to a shorter operation time because of the short 'touch', no pressure and precision which is a big advantage for therapists as well as customers/patients.



1398 SÜDA transmission handpiece

Simply remove the handpiece from the micromotor.



...and attach the OMNICUT attachment or transmission handpiece.



Spray Vac Unit

Fine even spray and powerful suction

TECHNICAL SPECIFICATIONS

Control unit/cabinet:

- Electrical connection: 230 V / 50 Hz
- Suction capacity max.: 210 l / min
- Noise level: 49 db (A) at 100% and 45 db (A) at 20% suction power
- Dimensions: W 100 x H 81 x D 48 cm
- Weight: 80 kg

Aspiration handpiece:

- Speed: 6.000 - 40.000 rpm
- Dimensions: Ø 23.5 x L 134 mm
- Weight: 112 g

Spray handpiece SÜDA:

- Speed: 6.000 - 25.000 rpm
- Dimensions: Ø 20 x L 140 mm
- Weight: 163 g

Spray handpiece NSK:

- Speed: 6.000 - 40.000 rpm
- Dimensions: Ø 22 x L 149 mm
- Weight: 151 g

Please note the speed of the spray handpieces:

SÜDA HS-WS M4 L

Recommended up to max. 25,000 rpm
NSK PX65L

Recommended up to max. 40,000 rpm

1231-00	Spray Vac Unit With spray handpiece SÜDA HS-WS M4 L
---------	---

1231-100	Spray Vac Unit With spray handpiece NSK PX65L
----------	---

If, in contrast to the standard version, you need a holder for Waldmann magnifying lamps, please order the article 1175.005W and replace it on the cabinet.



SÜDA wet and suction technology

The new Studio Line cabinet system combines wet and suction technology. Benefit from the advantages of a compact treatment unit with perfectly adjustable performance parameters.

In-house development - "made in Germany"

Branded products from SÜDA are "made in Germany" and come almost exclusively from the idea through the development right up to production and logistics directly from the production site and logistics center near Leipzig.



Fine uniform spray mist thanks to special technology

Extremely powerful and maintenance-free suction turbine for high traction

The large, clear control panel provides intuitive user-friendliness.

A UV compartment is standard equipment.

Integrated three-function nozzle

3 drawers and a shelf in the right side-chamber create optimal storage space.

The large work surface offers enough space with the greatest possible working space economy.

The lamp holder integrated into the worktable ensures perfect compatibility with all common magnifying lamps models.

Resistant materials as on the cabinet surface and the work surface, make the unit particularly insensitive and hygienic.

Thanks to the robust castors, the cabinet system can be moved effortlessly.

The service flap ensures uncomplicated, easily accessible maintenance.

MADE IN GER

Compressed air technology

Studio Line



With a very interesting price / performance ratio, air operated dental units offer a real alternative to micro-motor-driven treatment systems. Since the so-called pneumatic or compressed air systems are not driven by an electric / micro motor, but with compressed air, these systems are almost maintenance free and are characterized by extremely low operating costs.

The SÜDA compressed air treatment unit can optionally be configured with various spray and turbine handpieces. The three-function nozzle is included as standard.



SPRAY HANDPIECES

The air-powered motor enables vibration-free cutting in all rotational speeds. The motor is firmly secured in the handpiece by locking the cutter, independently of the working direction or speed.



TURBINE HANDPIECES

With a turbine you can perform finest work, fast and extremely precise. The distance for e.g. a Clavus (corn) does not need any preparation. The comfortable work without any vibration and also - due to the high speed - without pressure, makes the treatment for the patient absolutely painless and therefore extremely pleasant.



THREE-WAY NOZZLE

The three-way nozzle is used for a quick preparation of the work area: rinsing, cleaning or drying. Secondly, the three functions are essential during treatment, for example, when working on a unguis (Nail), the area in the narrow furrow can be reliably rinsed and dried in order to then treat quickly further.

podoAIRspeed

Pneumatic treatment system - convincing in price, performance and design



MADE IN
GER

SÜDA Compressed air system

The compressed air powered treatment unit podo AirSpeed opened both in its technical components, the ergonomics and efficiency as well as functionality and design, new dimensions in the compressed air-powered treatment equipment.

With its extremely favorable price / performance ratio, it provides for small foot care institutes a real alternative to micro-motor operated suction / or wet engineering systems. Air motors with spray handpiece and turbine: The treatment functions of compressed air units can be regulate during treatment comfortable by use of the foot switch, so that you can work intensively, accurately and highly efficiently.

Good accessibility of the control valves and filters for a pleasant workplace economy.

Easy filling of the water tank.

Optimized switch and control functions allow intuitive operation of the device.

The large, clear control panel ensures the best ease of use.

A UV compartment is standard equipment.

Integrated three-function nozzle

The cabinet unit offers plenty of storage space and also has a large storage area.

The lamp holder integrated into the large work surface ensures perfect compatibility with all common magnifying lamp models.

Resistant materials as on the cabinet surface and the work surface, make the unit particularly insensitive and hygienic.

Thanks to the robust castors the cabinet system can be moved effortlessly.

High maintenance-friendliness due to easy access to all relevant functional parts.



Speed control via foot switch

Rocker lever for switching the fine spray function on/off

Blow-Up function for drying and cleaning the working field



TECHNICAL SPECIFICATIONS

Control unit/cabinet:

- Electrical connection: 230 V / 50 Hz
- Dimensions: W 100 x H 81 x D 48 cm
- Weight: 76 kg

Spray handpiece (with air motor):

- Speed: up to 27.000 rpm
- Dimensions: Ø 20 x L 165 mm
- Weight: 207 g

Turbine handpiece:

- Speed: up to 400.000 rpm
- Dimensions: Ø 16.4 x L 120 mm
- Weight: 65 g

1229-00 podoAIRspeed

ACCESSORIES

1389.9	STERIA Spray, service spray for air motors, handpieces and turbines, 500 ml
1389.2	Nozzle standard
1389.6	Nozzle for turbines
1256.002	Connection set with 3 m hose for extension
1256.001	Hose per additional meter
1229-00	podoAIRspeed Base model

The scope of delivery of the treatment unit as basic model includes:

- Cabinet body white with rail
- Integrated UV compartment
- 2 x magnifying lamp holder standard
- Air motor with spray handpiece SÜDA HS-WS M4
- Turbine handpiece Triple Spray / SÜDA TT-WS M4
- Three functional nozzle
- 4 castors, thereof 2 can be locked

CONFIGURATION POSSIBILITIES HANDPIECES

	TURBINE HANDPIECE	SPRAY HANDPIECE	1229-XX* SURCHARGE
Turbines handpiece Triple Spray / SÜDA TT-WS M4	0		-
Turbine handpiece LED / SÜDA TLT-WS M4	1		surcharge
Turbine handpiece NSK Pana Max	2		surcharge
Spray handpiece SÜDA HS-WS M4		0	-
Spray handpiece NSK PX 65		1	surcharge

* Please add the number code that results from the handpiece configuration to the item number. This gives you the order number for the podoAIRspeed in the equipment variant you have chosen. If, in contrast to the standard version, you need a holder for Waldmann magnifying lamps, please order the item 1175.005W and replace it on the cabinet.

Handpiece SÜDA HS-WS M4/
Air motor SÜDA AM-WS M4



Handpiece PX65/
Air motor SÜDA AM-WS M4



Turbine handpiece triple
Spray / SÜDA TT-WS M4



Turbine handpiece with LED/
SÜDA TLT-WS M4 with light function



Turbine handpiece / NSK Pana Max



Three functional nozzle



Compressor podoAIRspeed

At SÜDA a dry-running air compressor is used which delivers oil-free compressed air directly to the control module of the treatment unit by means of a thin hose. The technological differences between the SÜDA air compressor and the commercially available oil-air compressor should be taken into account when making your purchase decision. The advantage of dry-type dryers is the supply of oil-free air, which prevents the penetration of oil-laden air into the system and significantly less contaminates the device. This makes the system almost maintenance-free, which makes our SÜDA air compressor an excellent choice.

1) Pressure regulator / condensate separator / shut-off valve (red)

2) Pressure gauge / safety valve / pressure switch



EQUIPMENT FEATURES

- Front door for compressor access
- Robust work surface
- 4 castors, of which 2 can be locked

TECHNICAL SPECIFICATIONS

Electrical connection: 230 V / 50 Hz
Power consumption: 1150 W
Suction capacity: 230 l / min
Dimensions: W 69 x D 77 x H 52 cm

Weight:

- Cabinet with compressor: 82 kg
- Cabinet without compressor: 38 kg
- Compressor: 44 kg

1252.2 Compressor cabinet white, 30 liters

1252.3 Cabinet without compressor

1252.10 Compressor 30 liters



Flexi Line

Flexible treatment equipment
for foot care and podology

SUCTION TECHNOLOGY

Flexi Line



Since developing the first pedicure device with integrated suction SÜDA has optimized these technologies, such as the successful Pedo Sprint series, in close collaboration with chiropodists and podiatrists. At the same time new technologies are revolutionizing the market such as the SÜDA vac.

Pedo Sprint ECO silentpower

With the ECO PEDO Sprint in the Sprint series, we meet the demand for a simple and powerful foot care system. This system is characterized by the ingeniously simple in operation and a sensational price-performance ratio. The Pedo Sprint ECO did not compromise the quality of the components used.

FEATURES

- Compact, easy to clean housing with disinfectable surfaces
- Micro motor handpiece with flexible suction hose and dust cover
- Proven centrifugal force for fast milling change
- Analog rotary control for suction power and milling speed
- Side handpiece holder
- Wired footswitch (optional)



Cases >P. 140

1936	Pedo Sprint ECO SP
------	--------------------

ACCESSORIES

8921	Podiatry case „PSmobil“
8926	Podiatry case “Mobil”
1630	Podiatry case “Smartcase”
8913.12	Podiatry case „MINIPRAX“
8211.024	Foot switch
8336	Handpiece holder
8415.2	Micropore dust bags (pack of 5 pcs)
1997.001	Change collet Set
35600	Gear box silentpower

TECHNICAL SPECIFICATIONS

Control unit:

- Electrical connection: 220-240 V, 50-60 Hz
- Suction capacity max.: 180 l / min
- Noise level: 60 db (A) at 100% and 53 db (A) at 20% suction power
- Dimensions: B 250 x H 140 x D 185 mm (including handpiece holder)
- Weight: 3 kg

Micro motor handpiece:

- Speed: 6.000 - 16.000 rpm
- Dimensions: Ø 23.5 x L 134 mm
- Weight: 119 g, with hose connection 340 g

SCOPE OF DELIVERY

- Foot care device Pedo Sprint ECO SP incl. suction handpiece
- 1 Micropore dust bag (8415.2)
- Power cord

Pedo Sprint 50 silentpower

The main focus of the Pedo Sprint 50 SP with a maximum micromotor speed of 25,000 rpm is the high torque. Especially in the low speed range, the device is characterized by a particularly strong pulling force. For this reason, it is also ideally suited for the use of larger instruments, such as a cap grinder.

FEATURES

- Compact, easy to clean housing with disinfectable surfaces
- Micromotor handpiece with flexible suction hose and dust cover
- Proven centrifugal force for fast milling change
- Operation via keypad to control suction power and milling speed
- Display for digital speed display
- Right / left rotation
- Display of dust bag change
- Handpiece holder on the side
- Wired foot switch (optional)



Cases >P. 140

1946	Pedo Sprint 50 SP
------	-------------------

ACCESSORIES

8921	Podiatry case „PS Mobil“
------	--------------------------

8926	Podiatry case „Mobil“
------	-----------------------

1630	Podiatry case „Smartcase“
------	---------------------------

8913.12	Podiatry case „Miniprax“
---------	--------------------------

8211.024	Foot switch
----------	-------------

8336	Handpiece holder
------	------------------

8415.2	Micropore dust bags (pack of 5 pcs)
--------	-------------------------------------

1997.001	Change collet Set
----------	-------------------

35600	Gear box silentpower
-------	----------------------

TECHNICAL SPECIFICATIONS

Control unit:

- Electrical connection: 220-240 V AC, 50-60 Hz
- Suction capacity max.: 180 l / min
- Noise level: 60 db (A) at 100% and 53 db (A) at 20% suction power
- Dimensions: B 250 x H 140 x D 185 mm (including handpiece holder)
- Weight: 3 kg

Micromotor handpiece:

- Speed: 6.000 - 25.000 rpm
- Dimensions: Ø 23.5 x L 134 mm
- Weight: 119 g, with hose connection 340 g

SCOPE OF DELIVERY

- Foot care device Pedo Sprint 50 SP incl. suction handpiece
- Micropore dust bag (8415.2)
- Power cord

Pedo Sprint 100 silentpower

The large clear LCD display of the Pedo Sprint 100 SP shows all important operating parameters such as suction power, speed, direction of rotation as well as maintenance and repair instructions. The "thinking" electronics reminds us when a dust bag and filter change becomes due and also points in time to maintenance dates.

FEATURES

- Compact, easy to clean housing with disinfectable surfaces
- Micromotor handpiece with flexible suction hose and dust cover
- Proven centrifugal force for fast milling change
- ON / OFF button directly on the handpiece
- Operation via keyboard to control suction power and milling speed
- Memory function for storing settings
- Right / left rotation
- Display for digital display of suction power, speed, direction of rotation, dust bag change, maintenance / repair
- Handpiece holder on the side
- Wired foot switch (optional)



Cases >P. 140

1956	Pedo Sprint 100 SP
------	--------------------

ACCESSORIES

8921	Podiatry case „PS Mobil“
------	--------------------------

8926	Podiatry case „Mobil“
------	-----------------------

1630	Podiatry case „Smartcase“
------	---------------------------

8913.12	Podiatry case „Miniprax“
---------	--------------------------

8211.024	Foot switch
----------	-------------

8336	Handpiece holder
------	------------------

8415.2	Micropore dust bags (pack of 5 pcs)
--------	-------------------------------------

1997.001	Change collet Set
----------	-------------------

35600	Gear box silentpower
-------	----------------------

TECHNICAL SPECIFICATIONS

Control unit:

- Electrical connection: 220-240 V AC, 50-60 Hz
- Suction capacity max.: 180 l / min
- Noise level: 60 db (A) at 100% and 53 db (A) at 20% suction power
- Dimensions: B 250 x H 140 x D 185 mm (including handpiece holder)
- Weight: 3 kg

Micromotor handpiece:

- Speed: 6.000 - 40.000 rpm
- Dimensions: Ø 23.5 x L 134 mm
- Weight: 141 g, with hose connection 365 g

SCOPE OF DELIVERY

- Foot care device Pedo Sprint 100 SP incl. suction handpiece
- Micropore dust bag (8415.2)
- Power cord

SÜDA vac S

The quietest unit with integrated suction

Quiet Revolution

The SÜDA Silent Air technology makes this unit with 49 dB(A) the quietest in the market.

Operation

User-friendly front panel with intelligent display easily and intuitively in all operations.

Guaranteed no maintenance

The innovative suction turbine works largely without wear. A guarantee of 3 years will save costly maintenance.

Handpiece lightweight

Effortless work thanks to perfect ergonomics – without compromising on high engine power and suction power.

Hygiene

The valuable aluminum surface can be cleaned perfectly. Filters are changed easily and hygienically.

Mobility

Das SÜDA vac is robust, lightweight, compact in dimensions and



Since developing the first pedicure device with integrated suction, SÜDA has continuously developed through its own research. Particular attention was now again placed on maximizing the performance of the suction in conjunction with the clearly noticeable reduction in noise.

FEATURES

- Compact, easy to clean housing with disinfectable surfaces
- Micromotor handpiece with flexible suction hose and dust cover
- Proven centrifugal force for fast milling change
- ON / OFF button directly on the handpiece
- Operation via keypad to control suction power and cutter speed
- Memory function for storing settings
- Right / left rotation
- Display for digital display of suction power, speed, direction of rotation, dust bag change, maintenance / repair
- Side handpiece holder
- Wired footswitch (optional)



Cases > P. 140

LISTEN BY YOURSELF

Convince yourself now of the "quiet revolution" and listen to the new dimension of suction technology online.

SCOPE OF DELIVERY

- Foot care device SÜDA vac S including suction handpiece
- 1 Micropore dust bag (8415.2)
- Power cord

TECHNICAL SPECIFICATIONS

Control unit:

- Electrical connection: 110-240 V / 50-60 Hz
- Power: 160 W
- Suction capacity max.: 210 l / min
- Noise level: 49 db (A) at 100% and 45 db (A) at 20% suction power
- Dimensions: W 24 x H 12.5 x D 22.5 cm
- Weight: 3.5 kg

Micro motor handpiece:

- Speed: 6.000 - 40.000 rpm
- Dimensions: Ø 23.5 x L 134 mm
- Weight: 112 g, with hose connection 317 g

1997	SÜDA vac S
1997.1	SÜDA vac S with wireless foot switch

ACCESSORIES

8930.3	SÜDA mobility Working unit
8926	Foot care case "Mobil"
8211.024	Foot switch, wired
8336	Handpiece support
8415.2	Micropore dust bag, 1 pack of 5 pieces
1997.001	Collet chuck set
36077	Replacement filter SÜDA vac S

Wet technology

Flexi Line



For decades SÜDA offers a wide range of wet technology devices. Special emphasis is placed on the continuous refinement of the spray image. This in turn is considered to be crucial for treatment efficiency and quality.

ASTERIA 1000

Lightweight with mini-handpiece for maximum comfort

The ASTERIA 1000 wet technology foot care device offers a lot of performance at an incredibly low price. This is mainly due to our innovative optimized hand piece technology. Thus the ASTERIA 1000 represents the ideal device for the entry into the wet technology or the perfect second device for your practice or mobile use.

FEATURES

- Easy to clean housing with disinfectable surfaces
- Micromotor handpiece with flexible spray hose
- Fine fogging, even when using distilled water
- Clean head system
- Easy removable water tank for easy refilling
- Proven centrifugal force for fast milling change
- Clear control panel for control of spray power and cutter speed
- Right / left rotation
- Digital speed display
- Side handpiece holder
- Wired footswitch (optional)

SCOPE OF DELIVERY

- Foot care device ASTERIA 1000 incl.
- Spray handpiece PD 30
- Power cord

TECHNICAL SPECIFICATIONS

- Control unit:
- Electrical connection: 100-230 V / 50 Hz
 - Power: 100 W
 - Dimensions: W 22 x H 14 x D 18 cm
 - Weight: 2.3 kg

- Micro motor handpiece:
- Speed: 6.000 - 30.000 rpm
 - Dimensions: Ø 22 x L 150 mm
 - Weight: 125 g



Cases > P. 140

1480	ASTERIA 1000
------	--------------

ACCESSORIES

8924	Podiatry case „AS Mobil“
------	--------------------------

8926	Podiatry case „Mobil“
------	-----------------------

1630	Podiatry case „Smartcase“
------	---------------------------

8913.12	Podiatry case „Miniprax“
---------	--------------------------

8211.024	Foot switch, wired
----------	--------------------

8336	Handpiece support
------	-------------------

9320.600	CARE & CLEAN NT Special, 100 ml
----------	---------------------------------

8334.001	Podo Liquid, 5l canister
----------	--------------------------

8334.007	Aqua distillate, 5l canister
----------	------------------------------

PEDOSPRAY NT 40

A pedicure system which combines many years of experience in engineering and modern technology with the latest findings of modern foot care. New standards were set in the pedicure equipment with the development of PEDOSPRAY NT 40



Cases > P. 140

SCOPE OF DELIVERY

- Foot care device PEDOSPRAY NT 40 incl. Spray handle Pododent
- Power cord

TECHNICAL SPECIFICATIONS

Control unit:

- Electrical connection: 230 V / 50 Hz
- Power: 140 W
- Dimensions: W 24 x H 9,5 x D 25 cm
- Weight: 3.5 kg

Micro motor handpiece:

- Speed: 5.000 - 40.000 rpm
- Dimensions: Ø 22 x L 135 mm
- Weight: 120 g

1507	PEDOSPRAY NT 40
------	-----------------

ACCESSORIES

1616.5	Podiatry case „de Luxe“
8211.024	Foot switch, wired
8336	Handpiece support
9320.600	CARE & CLEAN NT Special, 100 ml
8334.001	Podo Liquid, 5l canister
8334.007	Aqua distillate, 5l canister

FEATURES

- Easy to clean housing with disinfected surfaces
- Micromotor handpiece with flexible spray hose
- Fine fogging, even when using distilled water
- Clean head system
- Easy removable water tank for easy refilling
- Proven centrifugal force for fast milling change
- Right / left rotation
- Large control panel with hygienic keyboard keypad and digital display for controlling milling speed, direction of rotation and spray switching
- Individual spray regulation for optimal dust binding and cooling
- Memory function for storing settings
- Side handpiece holder
- Wired footswitch (optional)

PEDOSPRAY NT 40 light

Powerful wet technology foot care device based on the PEDOSPRAY NT 40 with a particularly light handpiece. 10 LED illuminate the field of treatment at the level of natural daylight.



Cases > P. 140

SCOPE OF DELIVERY

- Foot care device PEDOSPRAY NT 40 light incl. Spray handpiece PDL 40 with light function
- Power cord

TECHNICAL SPECIFICATIONS

Control unit:

- Electrical connection: 230 V / 50 Hz
- Power: 140 W
- Dimensions: W 24 x H 9,5 x D 25 cm
- Weight: 3.5 kg

Micro motor handpiece:

- Speed: 8.000 - 40.000 rpm
- Dimensions: Ø 22 x L 135 mm
- Weight: 115 g

1507.2	PEDOSPRAY NT 40 light
--------	-----------------------

ACCESSORIES

1616.5	Podiatry case „de Luxe“
8211.024	Foot switch, wired
8336	Handpiece support
9320.600	CARE & CLEAN NT Special, 100 ml
8334.001	Podo Liquid, 5l canister
8334.007	Aqua distillate, 5l canister

FEATURES

- Easy to clean housing with disinfected surfaces
- Micromotor handpiece with flexible spray hose
- Handpiece tip with switchable LED light ring for illuminating the working field
- Fine fogging, even when using distilled water
- Clean head system
- Easy removable water tank for easy refilling
- Proven centrifugal force for fast milling change
- Right / left rotation
- Large control panel with hygienic keyboard keypad and digital display for controlling milling speed, direction of rotation, spray and light switching
- Individual spray regulation for optimal dust binding and cooling
- Memory function for storing settings
- Side handpiece holder
- Wired footswitch (optional)

SÜDATRONIC Komfort

The shorter and comfortable Profi-Line handpiece PX65 in conjunction with the NBX micro motor distinguished itself by its carrying capacity. Both components are ideally suited for continuous operation.



Cases > P. 140

SCOPE OF DELIVERY

- Foot care device SÜDATRONIC Komfort including professional line handpiece PX65 and micromotor NBX
- Power cord
- Steria Spray 500 ml

TECHNICAL SPECIFICATIONS

Control unit:

- Electrical connection: 230 V / 50 Hz
- Power: 140 W
- Dimensions: W 26,9 x H 13,5 x D 19,5 cm
- Weight: 4 kg

Micro motor handpiece:

- Speed: 5.000 - 40.000 rpm
- Dimensions: Ø 22 x L 149 mm
- Weight: 148 g

1540.6	SÜDATRONIC Komfort
--------	--------------------

ACCESSORIES

8923	Foot care case „WS mobil“
8211.024	Foot switch, wired
8336	Handpiece support
1398	Dental attachment handpiece
1503.8	OMNICUT attachment
9320.600	CARE & CLEAN NT Special, 100 ml
8334.001	Podo Liquid, 5l canister
8334.007	Aqua distillate, 5l canister

FEATURES

- Easy to clean housing with disinfectable surfaces
- Micromotor dental handpiece with flexible spray hose
- Micromotor and handpiece with light titanium housing and scratch-resistant DURACOAT surface
- Fine fogging, even when using distilled water
- Clean head system
- Easy removable water tank for easy refilling
- Torsion system for secure holding of the milling cutters
- Right / left rotation
- Large control panel with hygienic keypad and digital display for controlling milling speed, direction of rotation and spray switching
- Individual spray regulation for optimal dust binding and cooling
- Memory function for storing settings
- Side handpiece holder
- Wired footswitch (optional)
- Option for Omnicut and translation handpiece (optional)
- Handpiece attachment without micromotor autoclavable

SÜDATRONIC Starlight

The Profi-Line handpiece PX65L offers in connection with the micromotor NBX LED an integrated light system for an optimal illumination of the treatment field.



Cases > P. 140

SCOPE OF DELIVERY

- Foot care device SÜDATRONIC Starlight incl. Profi-Line handpiece PX65L and Micromotor NBX LED
- Power cord
- Steria Spray 500 ml

TECHNICAL SPECIFICATIONS

Control unit:

- Electrical connection: 230 V / 50 Hz
- Power: 140 W
- Dimensions: W 26,9 x H 13,5 x D 19,5 cm
- Weight: 4 kg

Micro motor handpiece:

- Speed: 5.000 - 40.000 rpm
- Measurements: Ø 22 x L 149 mm
- Weight: 151 g
- Gewicht: 151 g

1540.8	SÜDATRONIC Starlight
--------	----------------------

ACCESSORIES

8923	Foot care case „WS mobil“
8211.024	Foot switch, wired
8336	Handpiece support
1398	Dental attachment handpiece
1503.8	OMNICUT attachment
9320.600	CARE & CLEAN NT Special, 100 ml
8334.001	Podo Liquid, 5l canister
8334.007	Aqua Destillata, 5l-Kanister

FEATURES

- Easy to clean housing with disinfectable surfaces
- Micromotor dental handpiece with flexible spray hose
- Micromotor and handpiece with light titanium housing and scratch-resistant DURACOAT surface
- Fine fogging, even when using distilled water
- Clean head system
- Easy removable water tank for easy refilling
- Torsion system for secure holding of the milling cutters
- Right / left rotation
- Large control panel with hygienic keypad and digital display for controlling milling speed, direction of rotation and spray switching
- Individual spray regulation for optimal dust binding and cooling
- Memory function for storing settings
- Side handpiece holder
- Wired footswitch (optional)
- Option for Omnicut and translation handpiece (optional)
- Handpiece attachment without micromotor autoclavable

SÜDATRONIC Starlight Power

In combination with the extremely light and compact micro-motor NLX nano, the Profi-Line handpiece PX65L has been designed in such a way that precise, vibration-free and quiet working combined with maximum comfort is possible under high stresses. The light, ergonomically optimized titanium housing also ensures low-fatigue work, the integrated light system for optimum illumination of the treatment surface.



SCOPE OF DELIVERY

- Foot care device SÜDATRONIC Starlight Power incl. Profi-Line handpiece PX65L and micromotor NLX nano
- Power cord
- Steria Spray 500 ml

TECHNICAL SPECIFICATIONS

Control unit:

- Electrical connection: 230 V / 50 Hz
- Power: 140 W
- Dimensions: W 26,9 x H 13,5 x W 19,5 cm
- Weight: 4 kg

Micro motor handpiece:

- Speed: 8.000 - 40.000 rpm
- Measures: Ø 22 x L 137 mm
- Weight: 129 g

1540.10	SÜDATRONIC Starlight Power
---------	----------------------------

ACCESSORIES

8923	Foot care case "WS mobil"
8211.024	Foot switch, wired
8336	Handpiece support
1398	Dental attachment handpiece
1503.8	OMNICUT attachment
9320.600	CARE & CLEAN NT Special, 100 ml
8334.001	Podo Liquid, 5l canister
8334.007	Aqua distillate, 5l canister

Cases > P. 140

FEATURES

- Easy to clean housing with disinfectable surfaces
- Micromotor dental handpiece with flexible spray hose
- Micromotor and handpiece with light titanium housing and scratch-resistant DURACOAT surface
- Fine fogging, even when using distilled water
- Clean head system
- Easy removable water tank for easy refilling
- Torsion system for secure holding of the milling cutters
- Right / left rotation
- Large control panel with hygienic keypad and digital display for controlling milling speed, direction of rotation and spray switching
- Individual spray regulation for optimal dust binding and cooling
- Memory function for storing settings
- Side handpiece holder
- Wired footswitch (optional)
- Option for Omnicut and translation handpiece (optional)
- Handpiece attachment without micromotor autoclavable

SÜDA spray-Line

Finest atomized spray, low noise, easy transport

Refined

Through our optimized hand piece, we achieved a fine mist for ideal dust suppression.

Quieter

The operating noise has been further reduced by the considerably lowered sound pressure.

Lighter

To facilitate mobile pedicure, the weight was further reduced by the innovative pump technology. A comfortable transport is secured.

Easy filling

The device operates with a unique new, unpressurized liquid system; therefore the water tank is not pressurized and thus allows an easy filling, which is also possible during operation.

Stronger hand piece

The powerful, brushless micro motors in the hand pieces allow work reliably.

Working safely

Due to the use of low-voltage-driven, electrical components and the use of an external 24 V power supply, the device is particularly safe.



SÜDA spray S

The spray has been adjusted in the light of medical, professional foot care and further optimized. Even with existing models of SÜDA wet technique you could rely on an outstanding dust control in which the new SÜDA spray with a further optimized binding of dust just allowed an ideal treatment of allergic persons.

FEATURES

- Compact, easy to clean housing with disinfectable surfaces
- Micromotor dental handpiece with flexible spray hose
- Fine fogging, even when using distilled water
- Pressureless filling system with easily removable water tank
- Torsion system for secure holding of the milling cutters
- Right / left rotation
- Large control panel with hygienic keyboard keypad and digital display for controlling milling speed, direction of rotation and spray switching
- Maintenance and repair indicator
- Individual spray regulation for optimal dust binding and cooling
- Memory function for storing settings
- Side handpiece holder
- Wired footswitch (optional)
- Wireless foot switch (optional depending on model)
- Option for Omnicut and dental attachment handpiece (optional)
- Handpiece attachment autoclavable



Cases >P. 140

Please note the speed of the spray handpieces:

SÜDA HS-WS M4
recommended to max. 25.000 rpm
BienAir PM 1:1
recommended to max. 40.000 rpm

1259.20	Control unit SÜDA spray S
1259.201	Control unit SÜDA spray S with wireless foot switch
1258.3	Handpiece SÜDA HS-WS M4
1258.9	Handpiece BienAir PM 1:1

ACCESSORIES

8930.3	SÜDA mobility Working unit
8926	Foot care case "Mobil"
8211.024	Foot switch, wired
8336	Handpiece support
9320.600	CARE&CLEAN NT-Spezial, 100 ml
8334.001	Podo-Liquid, 5 l-Can
8334.007	Aqua Destillata, 5 l-Can

SCOPE OF DELIVERY

- Control unit SÜDA spray S including spray hose and micromotor SÜDA MM-WS
- Handpiece SÜDA HS-WS M4 or BienAir PM 1:1
- Power cord
- Steria Spray 500 ml

TECHNICAL SPECIFICATIONS

Control unit

- Electrical connection: 100-240 V / 50 Hz (input), 24 V (output)
- Power: 120 W
- Dimensions: W 24 x H 12.2 x D 22.5 cm
- Weight: 3.3 kg

Micromotor handpiece SÜDA HS-WS M4

- Speed: 6.000 - 40.000 rpm
- Dimensions: Ø 20 x L 140 mm
- Weight: 163 g

Micromotor handpiece BienAir PM 1:1

- Speed: 6.000 - 40.000 rpm
- Dimensions: Ø 20 x L 144 mm
- Weight: 145 g

SÜDA spray L

Powerful wet technology foot care device based on the SÜDA spray S with a light handpiece for optimal illumination of the treatment field.



Cases >P. 140

Please note the speed of the spray handpieces:

SÜDA HS-WS M4
recommended to max. 25.000 rpm
NSK PX65L
recommended to max. 40.000 rpm

1259.21	Control unit SÜDA spray L
1259.211	Control unit SÜDA spray L mit Funkfußtaster
1258.6	Handpiece SÜDA HS-WS M4 L
1530.81	Handpiece NSK PX65L

ACCESSORIES

8930.3	SÜDA mobility Working unit
8926	Foot care case "Mobil"
8211.024	Foot switch, wired
8336	Handpiece support
9320.600	CARE&CLEAN NT-Spezial, 100 ml
8334.001	Podo-Liquid, 5 l-Can
8334.007	Aqua Destillata, 5 l-Can

FEATURES

- Compact, easy to clean housing with disinfectable surfaces
- Micromotor dental handpiece with flexible spray hose
- Fine fogging, even when using distilled water
- Pressureless filling system with easily removable water tank
- Torsion system for secure holding of the milling cutters
- Right / left rotation
- Large control panel with hygienic keyboard keypad and digital display for controlling milling speed, direction of rotation and spray switching
- Maintenance and repair indicator
- Individual spray regulation for optimal dust binding and cooling
- Memory function for storing settings
- Side handpiece holder
- Wired footswitch (optional)
- Wireless foot switch (optional depending on model)
- Option for Omnicut and dental attachment handpiece (optional)
- Handpiece attachment autoclavable

SCOPE OF DELIVERY

- Control unit SÜDA spray S including spray hose and micromotor SÜDA MM-WS L
- Handpiece SÜDA HS-WS M4 or NSK PX65L
- Power cord
- Steria Spray 500 ml

TECHNICAL SPECIFICATIONS

Control unit:

- Electrical connection: 100-240 V / 50 Hz (input), 24 V (output)
- Power: 120 W
- Dimensions: W 24 x H 12.2 x D 22.5 cm
- Weight: 3.3 kg

Micromotor handpiece SÜDA HS-WS M4:

- Speed: 6.000 - 40.000 rpm
- Dimensions: Ø 20 x L 140 mm
- Weight: 163 g

Micromotor handpiece BienAir PM 1:1:

- Speed: 6.000 - 40.000 rpm
- Dimensions: Ø 20 x L 144 mm
- Weight: 145 g

Nail M110

Nailer without suction for professional nail design



The Nail M110 offers technical innovations and modern design. So the device is to be mounted directly on a waist belt on the body and can be used by rechargeable batteries even without cable. For professional natural and artificial nail treatment in continuous use.

SCOPE OF DELIVERY

- Nail polisher M110
- Micromotor handpiece
- Footswitch
- Equipment and handpiece support
- Waist belt
- Charger

TECHNICAL SPECIFICATIONS

Control unit:

- Electrical connection: 100-240 V / 50-60 Hz (input), 13 V (output)
- Recharge time: 2.5 hours
- Operating time: 5 hours (with maximum rotation) 9 hours (with medium rotation)
- Dimensions: W 6,7 x H 3 x D 13,2 cm
- Weight: 185 kg

Micromotor handpiece:

- Speed: max. 25.000 rpm
- Dimensions: Ø 24 x L 136 mm
- Weight: 120 g

1436 Nail M110

Suction technology accessories

MICRO-DUST BAG FOR PEDO SPRINT / SÜDA VAC

8415.2 1 pack of 5 pieces

Dust bag suitable for all devices PEDO Sprint 100/50 / ECO and SÜDA vac S. By the 5-layers special fabric with ultrasonically welded seams, best filter characteristics are guaranteed.



MICROPORE DUST BAG FOR AIR AS1 / PEDO FS3

8412.1 1 pack of 5 pieces

Bag that fit the pedicure devices AIR AS1 / PEDO FS3. Through the five-layered special fabric with ultrasound welding seams, best filter characteristics are guaranteed.



MEDIVAC DUST BAG

1995.001 Micropore dust bags Medivac, 5 pieces

Textured, multi-layer depth filter according to IFA. Micro particulate filter with the highest dust separation quality (in D) to dust class M.



Handpieces for wet technology

Aufsatz- und Turbinenhandstücke für Spray Vac Unit, podoAIRspeed, SÜDATRONIC und SÜDA spray



ESSAY HANDPIECES (WITHOUT MICRO MOTOR)

1258.3	Handpiece SÜDA HS-WS M4
1258.9	Handpiece BienAir PM 1:1
1530.61	Handpiece NSK PX65

ESSAY HANDPIECES WITH LIGHT FUNCTION (WITHOUT MICRO MOTOR)

1258.6	Handpiece SÜDA HS-WS M4 L
1530.81	Handpiece NSK PX65L

TURBINE HANDPIECE



1258.1	Turbine handpiece Triple Spray / SÜDA TT-WS M4
1258.2	Turbine handpiece LED / SÜDA TLT-WS M4
1396.1	Turbine handpiece NSK Pana Max

OMNICUT

1503.8	Omnicut-essay
1503.800	blades 20 pieces (10 sets)



The rotating mini-blades carry even stubborn, thickened cornea in seconds. No milling dust is generated! The blades are interchangeable and sterilizable. 2 blades (= 1 set) included in delivery. Suitable for NBX, NBX LED, NLX nano and SÜDA MM-WS / WS L micro motors

SÜDA DENTAL ATTACHMENT HANDPIECE

1398	Dental attachment handpiece
------	-----------------------------



Handy and light handpiece with up to 200,000 rpm. Precise, fast and precise work. Pleasant, painless treatment. Suitable for NBX, NBX LED, NLX nano and SÜDA MM-WS / WS L micro motors

Accessories for wet technology



AQUA DESTILLATA

8334.007	Can, 5 l
Distilled water. Demineralized water gem. VDE 0510, an essential prerequisite for trouble-free work with wet technology devices.	



PODO-LIQUID

8334.001	Can, 5 l
Distilled water. Demineralized water gem. VDE 0510, an essential prerequisite for trouble-free work with wet technology devices.	



CARE&CLEAN NT-SPEZIAL

9320.600	100 ml
Alcoholic, odorless, germicidal concentrate to improve the atomization, by breaking the surface tension, and for longer shelf life of distilled water. Mix contents of the bottle with 5 liters of distilled water.	



WL-CLEAN / WL-DRY

1389.10	Combination pack 500 ml/300 ml
Universal, ready-to-use, aldehyde- and alcohol-free rinsing and cleaning solution for non-protein-fixing internal cleaning and decontamination before disinfection / sterilization of hollow body instruments such as turbines, handpieces and elbows. With drying spray for rapid and dry cleaning of the rinsed instruments.	



WL-ADAPTER

1389.101	WL adapter 02 made of stainless steel
Adapter for spray water and spray air ducts of all handpieces and angle pieces with ISOS interface	

Accessories for wet technology



STERIA SPRAY

1389.9 500 ml-Tin

1389.2 Nozzle

1389.6 Nozzle for turbines

Your high-performance hand pieces need regular maintenance. Under normal loading daily 1 x, with heavy loads several times a day, flush with the STERIA Spray - the and pieces will thank you through a long life. Nozzle ordered separately if necessary.



W&H SERVICE OIL F1

1390 400 ml can

For the care of turbine handpieces, spray handpieces and air motors. For a longer life of the handpieces. Sterilization-resistant and extremely economical.



PANA SPRAY PLUS

1389.1 500 ml-can

Specially suitable for turbine handpieces, spray handpieces and air motors. For effective cleaning and lubrication before each autoclaving process.



BIOTEST WATER QUALITY

9320.810 Biotest, 10 tests / pack

Water quality - germ-free and safe: With the Biotest you will find a perfect way to complement your hygiene management to a simple, sensible option for measuring your water quality.



BILPRON

9320.800 Bilpron, 1 l

With Bilpron you get the ideal hygienic solution for cleansing the waterways. Removes both organic and inorganic residues safely and reliably. Has a broad spectrum of germs. Prevents new formation of biofilm. Minimum of 12 hours in the unit.

Accessories



HANDPIECE SUPPORT

8336

Accessory handpiece support made of robust soft plastic. Absolutely non-slip and washable - for safe working on all work surfaces. Dimensions: W 130 x H 30 x D 48 mm



HANDPIECE SUPPORT UNIVERSAL

8920.102

Accessory handpiece support made of robust soft plastic. Absolutely non-slip and washable - for safe working on all work surfaces. With slots for accommodating six rotating instruments. Dimensions: B 150 x H 12 x D 60 mm



FOOT SWITCH

8211.024 Foot switch, wired

Fits to systems of Pedo Sprint series, Asteria 1000, Pedospray NT 40 / NT 40 light, SÜDATRONIC and SÜDA spray

Mobile Foot care

136 - 139

Trolley system
SÜDA mobility

140 - 147

Foot care cases

148

Case lamps

149

Legrests

150 - 151

Accessories

Trolley system



SÜDA, as your long-term partner in the mobile pedicure, proves again its innovative strength in developing the pedicure suitcase on wheels practically new. The following points were crucial in the development:

Total mobility

The innovative trolley system makes transport easy not only on flat surfaces but also on slopes. You can choose according to your needs between the ball bearing and stair climber trolley.

Functionality is the key

The two-piece, intelligent trolley system was designed directly on the transport with the trolley system but also works well independently. In this case, the system is divided into two separate units, both of them can be fixed on the trolley. For simple home visits, it is possible to transport only the work unit (AE) manually to the customer. When treating several

customers at a same place, such as in a nursing home, it may be necessary, to have the supply unit as well on site. In this combination, the work unit is placed on the supply unit, whereby the ideal working height for the foot care device is achieved and you can work directly from the trolley.

Sophisticated look

The increased focus on the subject "hygiene" was assured that the case system is simple and almost germ-free clean. Materials and colors were ideally adapted to the requirements of mobile pedicure. Let you be convinced by this thoughtful, innovative case concept.

SÜDA mobility

Discover the unlimited possibilities of our trolley system and fill out the new dimension of mobile pedicure.

MADE IN
GER

Additional storage space is always very helpful especially when you have to work on different customers on a same place. The supply unit offers you the possibility to store larger quantities of consumables or sales products.

Whether you use the work unit in combination with the trolley or rely only on the work unit - you'll easily be able to accommodate pedicure device magnifying lamp and tray.

Both the resistant outer material as well as the entire interior are almost germ-free clean. Do not allow dirt and moisture while using a suitable cleaning solution from SÜDA.

The sturdy wheels of the stair climber trolley could be retrofitted by the ball bearing model. So you are really prepared for any type of ground. The retrofitting of the wheels is done quickly and easily via a click mechanism.



SÜDA mobility-Set



SÜDA MOBILITY WORK UNIT

8930.3 Work unit (AE)

Dimensions W 35 x H 25 x D 35 cm

Features: Removable and separate use of AE Placed on the trolley VE the AE reaches the ideal working height without that the pedicure device must be taken out of the AE Storage space for all current SÜDA devices Strap available (included) Support for magnifying lamp and tray Storage compartments (for eg footswitch) Foldable outside pockets for extra storage.



SÜDA MOBILITY CARRYING AID

8930.5 Carrying aid

Dimensions Ø 15, L 65 cm

Features: Transport aid with shoulder strap for footrest and/or stool by single use of the AE.



SÜDA MOBILITY SUPPLY UNIT

8930.4 Supply Unit (VE)

Dimensions: W 35 x H 48 x D 35 cm / Weight: VE (empty) 5.8 kg

Features: Removable, used separately VE. Case for instrument boxes. Additional storage space (e.g. for sales materials). Strap available (included). Footrest and / or stool side fixable to VE. Includes a drawer unit of 3 drawers.

Optional: A second drawer unit (8930.7) can be ordered if needed or instead of a second drawer unit, a cool bag (8930.6) can be ordered.



SÜDA MOBILITY COOL BAG

8930.6 Cool bag

Dimensions: W 19 x H 24 cm x 27 cm T

Features: Cooler bag for cooling consumables.

Can be integrated into VE.



SÜDA MOBILITY DRAWER BOX

8930.7 Drawer box

Dimensions: W 19 x H 24 x D 27 cm

Features: drawer unit to store materials. Can be integrated into VE.



SÜDA MOBILITY TROLLEYS

8930.1 Trolley stair climber

Dimensions: H 108 cm / Weight: Trolley stair climber (empty) 2.7 kg

Features: Ideal to climb stairs

SÜDA Recommendation: Perfect especially for urban areas suitable.

8930.2 Trolley ball bearing

Dimensions: H 108 cm / Weight: Trolley Ball Bearing (empty) 2.5 kg

Features: Rugged trolley with unicycle system

SÜDA Recommendation: Perfect especially for rural areas, and uneven ground suitable.



SÜDA MOBILITY TROLLEY WHEELS

Separate stair climber wheels as replacement or as supplement for your trolley.

Easy to attach or to replace on the trolley.

8930.8 Stair climber, 1 piece

8930.9 Ball bearings, 1 piece



SÜDA MOBILITY SETS

8930.100 Complete set of stair climber

Consisting of work unit, supply unit and trolley stair climber

8930.200 Complete set of ball bearings

Consisting of working unit, supply unit and trolley ball bearing

ACCESSORIES

8919 Suitcase light Flexilight

8206.36 Suitcase magnifying light
Vario LED

8206.361 Adapter for varioLEDflex

8206.38 Suitcase magnifying lamp
Mobil LED

8206.381 Adapter for Mobil LED

8204.3 Workstool PROFI

1104.1 Footrest MINI

7251 Instrument case stainless steel
(B 205 x H 40 x D 100 mm)

8311.11 Hygiene Box, 180 ml

8311.16 Hygiene-Box, 280 ml

8311.15 Hygiene-Box, 700 ml

4862 Melamine tray small
(W 190 x H 15 x D 150 mm)

4863 Melamine tray large
(W 235 x H 15 x D 180 mm)

9224.90 Combined dispenser

2172 Collection box for blades

SÜDA cases

With a function-optimized inner dividing, the comfort suitcase SÜDA mobile series provide a comprehensive and comfortable facility to go. Equipment and accessories are perfectly fitted into the hygienic, easy to clean hard-shell form, so that no wobbling or slippage is possible. The central positioning of the foot care device ensures optimal weight distribution, which has a positive effect on the carrying comfort.

Incl. comfortable, easily stowable tablet, used as a work surface

Compartment with holder for case light Flexi-Light

2 recessed sockets for connecting equipment and light

Lateral power connection for long power cord

Controlled hose guide with handpiece / fixation

Compartments for 2 hygienic boxes for cutters or fluids


Storage compartment for foot

Incl. 2 large boxes for instruments and dressing material

Very handy: you can easily change the filter without taking the device out of the case.



MADE IN GER

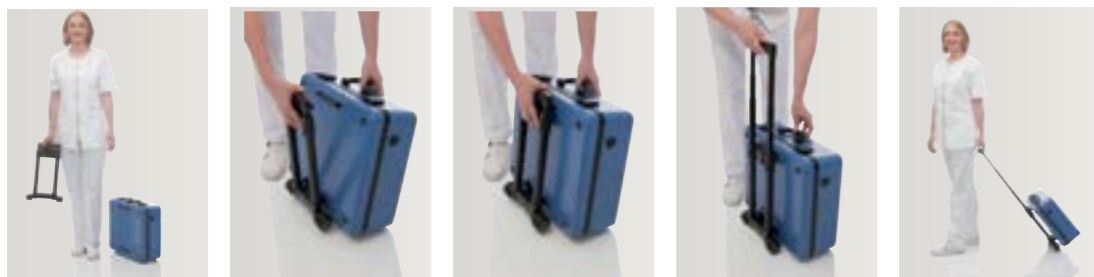

Transport roller is available as an option and can be easily installed.

PS mobil

Pedicure Cases "PS mobile" for all suction technology systems of Pedo Sprint

8921 Fußpflegekoffer »PS mobil«



ACCESSORIES

8920.100	Transport roller
8920.102	Handpiece holder
8919	Suitcase light Flexi-Light
8311.11	Hygiene box small (empty)
1614.200	Universal case cable, 5 m
8920.101	Carrying strap

SCOPE OF DELIVERY

- 2 × instrument boxes large
- Double socket,
- Tray
- Holder for case lamp 8919
- Holder for 2 × hygiene boxes 8311.11
- Holder for carrying strap 8920.101
- Holder for transport roller 8920.100
- Power cable angled
- Storage compartment for foot switch

TECHNICAL SPECIFICATIONS

Case Dimensions: W 48 × H 39.5 × D 16.5 cm
Weight: 5.5 kg

Equipment, case lamps, instruments and other decorative items are not included.

AS mobil

Foot care case "AS mobil" for wet technology system Asteria 1000.

8924 Foot care case »AS mobil«



ACCESSORIES

8920.100 Transport roller

8920.102 Handpiece holder

8919 Suitcase light Flexi-Light

8311.11 Hygiene box small (empty)

1614.200 Universal case cable, 5 m

8920.101 Carrying strap

SCOPE OF DELIVERY

· 2 × instrument boxes large

· Double socket,

· Tray

· Slots for bottles

· Holder for case lamp 8919

· Holder for 2 × hygiene boxes 8311.11

· Holder for carrying strap 8920.101

· Holder for transport roller 8920.100

· Power cable angled

· Storage compartment for foot switch

TECHNICAL SPECIFICATIONS

Equipment, case lamps, instruments and other decorative items are not included.

Case Dimensions: W 48 × H 39.5 × D 16.5 cm

Weight: 5.5 kg

WS mobil

Foot care case for all wet technology devices of the Südatronic and Waterspray series.

8923 Foot care case »WS mobil«



ACCESSORIES

8920.100 Transport roller 64,90 €

8920.102 Handpiece holder 14,90 €

8919 Suitcase light Flexi-Light 99,00 €

8311.11 Hygiene box small (empty) 2,90 €

1614.200 Universal case cable, 5 m 9,00 €

8920.101 Carrying strap 20,30 €

SCOPE OF DELIVERY

· 2 × instrument boxes large

· Double socket,

· Tray

· Slots for bottles

· Holder for case lamp 8919

· Holder for 2 × hygiene boxes 8311.11

· Holder for carrying strap 8920.101

· Holder for transport roller 8920.100

· Power cable angled

· Storage compartment for foot switch

TECHNICAL SPECIFICATIONS

Equipment, case lamps, instruments and other decorative items are not included.

Case Dimensions: W 48 × H 39.5 × D 16.5 cm

Weight: 5.5 kg

Mobil

Foot care case universally suitable for all suction devices of the series Pedo Sprint and SÜDA vac as well as for the wet technology devices SÜDA spray and Asteria 1000.

8926 Foot care case »Mobil«



ACCESSORIES

8920.100	Transport roller
8920.102	Handpiece holder
8206.38	Crankcase light Mobil LED
8206.381	Adapter for mobile LED
8919	Suitcase light Flexi-Light
8311.11	Hygiene box, 180 ml
8311.15	Hygiene box, 700 ml
1614.200	Universal case cable, 5 m
8920.101	Carrying strap

Equipment, case lamps, instruments and other decorative items are not included.

SCOPE OF DELIVERY

- Double socket
- Slots for hygienic boxes 8311.11 and 8311.15
- Fixing lugs for carrying strap 8920.101
- Holding device for transport roller 8920.100
- Power cable angled

TECHNICAL SPECIFICATIONS

Dimensions: B 56 x H 41 x D 16 cm
Weight: 4.8 kg

Miniprax

Foot care cases "Miniprax" for suction technology systems of Pedo-Sprint Series as well as the wet technology systems Asteria 1000, Waterspray and Südatronic.

8913.12 Foot care case »Miniprax«



ACCESSORIES

8919	Suitcase light Flexi-Light
2179.1	Drilling operation for instruments
8913.20	Carrying strap
1614.200	Universal case cable, 5 m

Equipment, case lamps, instruments and other decorative items are not included.

Sturdy case for daily use

Surfaces made of scratch-resistant plastic and therefore very resistant and durable

Attractive design: black with rounded abutting edges of aluminum. Compact format

Well-organized, functional „inner life“



SCOPE OF DELIVERY

- Double socket
- Mounting lugs for carrying strap
- 2 keys for the snap locks
- Drill essay for instruments
- Three-part hard shell for accessories
- Outside pocket
- Milling stand aluminum small 28800
- 2 x clear glass bottle 50 ml
- Power cable angled

TECHNICAL SPECIFICATIONS

Dimensions: B 42 x H 38 x D 24 cm
Weight: 4.5 kg

Smartcase

Foot care case "SmartCase" for all suction technology systems of Pedo-Sprint series and for wet technology system Asteria 1000.

1630 Foot care case »Smartcase«



ACCESSORIES

8919	Case lamp Flexi-Light
8913.20	Carrying strap
1614.200	Universal case cable, 5 m

SCOPE OF DELIVERY

- Accessory box for small parts
- Mounting lugs for carrying strap
- Simple socket
- Outside pocket
- Cutter case (empty) 8311.5
- Power cable angled

TECHNICAL SPECIFICATIONS

Equipment, case lamps, instruments and other decorative items are not included.

Case dimensions: W 49 x H 20 x D 33.5 cm
Weight: 5.2 kg

De Luxe

Foot care case "De Luxe" for wet technology systems Pedospray NT 40 and NT 40 light.

1616.5 Foot care case »De Luxe«



ACCESSORIES

8919	Case lamp Flexi-Light
8913.20	Carrying strap
1614.200	Universal case cable, 5 m

SCOPE OF DELIVERY

- Double socket
- Cutter case (empty) 8311.5
- Mounting lugs for carrying strap
- Outside pocket
- Clear glass bottle 50 ml
- Power cable angled

TECHNICAL SPECIFICATIONS

Case Dimensions: W 44 x H 14 x T 36 cm
Weight: 5 kg

Case lamps

varioLEDflex

The varioLEDflex stands out due to the exceptional illumination. In this case, the light is optimally focused, and then formed into a uniform light cone. The 63 durable SMD LEDs also ensure the high light output. An extremely comfortable, intuitive handling is achieved as well by the low weight, compact dimensions, ergonomic magnifying handle and the flexible, smooth gooseneck stand. Best cleaning properties allow easy cleaning and disinfection of the magnifying lamp. The low energy consumption gives the lamps a sensational life of about 50,000 hours.



FEATURES

- Large field of view with 76 mm lens Ø at 2.5 x magnification / 6 diopters
- Illumination: approx. 6,000 lx (in 150 mm working distance)
- Fastening via adapter in the case or on the table clamp
- Cable length: approx. 2.5 m

SUITCASE LIGHT varioLEDflex

8206.36	Suitcase magnifying lamp varioLEDflex
8206.361	Adapter varioLEDflex (for case)
8406.21	Table mounting
8206.363	Adapter varioLEDflex (for table mounting)



Mobil LED

CASE MAGNIFYING LAMP MOBILE LED

- 3 diopters glass lens
- Magnifier diameter 95mm
- Light module with 30 LEDs
- Mounting via adapter in case or via table clamp (included)

8206.38	Case magnifying lamp Mobil LED
8206.381	Adapter Mobil LED (for case)



Flexi-Light

CASE LAMP FLEXI-LIGHT

8919	Case lamp Flexi-Light
------	-----------------------

The LED case lamp with its very small power supply and the extremely flat lamp head can be placed space-saving in the case and is characterized by its high luminosity and low weight. Length: 70 cm. Lighting head: B 20 x H x D 1 3.5 cm



Legrests

Legrest MINI



FEATURES

With small legrest padding. Super light and handy with integrated carrying handle. Easy to use. Ideal for outpatient treatments. Collapsible light-weight tripod with height-adjustable center column and telescopic legs. Surface with easy-care plastic coating.

DIMENSIONS

Padding dimensions:
W 12 x H 3 x D 16 cm
Height adjustment: 54-115 cm
Weight: 1.1 kg

1104.1 Legrest MINI

Legrest MAXI



FEATURES

With large backrest padding. Super light and handy with integrated carrying handle. Easy to use. Ideal for ambulant treatments. Collapsible light-weight tripod with height-adjustable center column and telescopic legs. Surface with easy-care plastic coating.

DIMENSIONS

Padding dimensions:
W 14.5 x H 3 x D 21 cm
Height adjustment: 54-115 cm
Weight: 1.2 kg

1104.2 Legrest MAXI

Legrest KOMFORT



FEATURES

With ergonomically shaped leg support padding and adjustable tilt angle. More comfort and better support due to special design with lateral support surfaces. Extremely robust tripod for high stability and good ground adhesion. With height-adjustable center column and pull-out feet. Surface with easy-care plastic coating. Easy to fold.

DIMENSIONS

Padding dimensions:
W 13 x H 5 x D 21 cm
Height adjustment: 49-73 cm
Weight: 1.6 kg

1104.3 Legrest KOMFORT

Available from September 2017

Accessories



HANDPIECE SUPPORT

8336 Handpiece support

Additional handpiece support made of robust soft plastic. Absolutely non-slip and washable - for safe working on all work surfaces. Dimensions: W 130 x H 30 x D 48 mm



HANDPIECE SUPPORT UNIVERSAL

8920.102 Handpiece support universal

Accessory handrail made of robust soft plastic. Absolutely non-slip and washable - for safe working on all work surfaces. With slots for accommodating six rotating instruments. Dimensions: B 150 x H 12 x D 60 mm



FOOT SWITCH

8211.024 Foot switch, wired

This product is compatible with Pedo Sprint ECO /50/100, SÜDA vac, Asteria 1000, Pedospray NT 40 / NT 40 light, SÜDATRONIC and SÜDA spray.



UNIVERSAL CASE CABEL

1614.200 Universal case cable, 5 m

With the 5 m long universal cable case you are guaranteed connection! Power cord with Schuko plug and matching suitcase coupling.



ACCESSORY BAG HELPER

8913.9 Accessory bag helper

Spacious transportation help for your utensils for mobile foot care: Balsam and disinfection, the gang cutter, paving and much more - in the "Helper" you can place everything clearly and accessible. In that way you expand the space in the treatment suitcase and distribute the weight of your "mobile pedicure cabin" on both shoulders.

Features: Made from durable, black nylon, water-repellent. Robust zipper. Plastic tray as a flexible room divider; so you can organize the "inner life" of the helper on your own view. Wide shoulder strap for optimal weight distribution. Interior dimensions: W 39.5 x H 11.5 x X 29.5 cm. Weight: 2,0 kg



WORK STOOL

8204.3 Work stool Profi

8405.800 Shoulder bag for work stool

For optimal sitting in the outpatient pedicure. Excellent stability through stable tripod. Seat made of breathable natural leather. Is set in seconds! Seat height: 37 cm. Weight: about 1 kg.

Instruments

Big impact, small tool

154 - 173

Cutters and
grinders

178 - 197

Nippers and
scissors

174 / 198

Accessories

Cutters and grinders

A Strong Partner for Top-Quality Made in Germany

In order to be able to comprehensively satisfy the extremely sophisticated demands on rotating instruments for foot care and podiatry, the SÜDA cutter product range was developed together with a competent and very experienced partner – the company, Gebrüder Brasseler. They are known with the brand, KOMET, as one of the worldwide leading manufacturers of rotating dental instruments. Internationally, the brand sets standards for the highest quality, precision as well as reliability, and for almost 90 years now. Almost 1,000 employees work in Lemgo to manufacture the most modern dental instruments. The KOMET brand originated in a dental bur factory, which was founded in 1923. In the meantime, in addition to their products in the area of dentistry, the company also sells rotating instruments and saw blades for use in hip and knee endoprosthetics, bone, ENT and neurosurgery. The international jewellery industry, jewellers and goldsmiths are also supplied with the finest tools for working on rings, earrings and bracelets. In order to be able to implement the immense experience and know-how of the KOMET brand to meet the sophisticated demands of podiatry and chiropody perfectly, SÜDA took on the task of developing the SÜDA cutter product range with the united powers and support of experienced podiatrists and pedicurists. The result of this intense development process are now the instruments for all problems and application demands available at SÜDA.

Impressive Variety – Unbeatable in Price and Performance

The advantage of the SÜDA cutters lies not only in the superior quality, but also in a top price versus performance relationship. The variety offered is also convincing: In addition to the well-known favourites and standard forms of the previous cutters, you will now find several new shapes and materials for even more comfortable use. A clear plus are the professionally sealed plastic foil packaging, which guarantees you the first use of these medical products.

All essential information, for example, article numbers for re-orders as well as maximum speed and optimal speed for treatment can be found on the label of the packaging.

The appropriate cutters for your needs

Whether for grinding callused skin, cutting nails, treating fungus or polishing artificial nails – for every type of work, there are rotating instruments in countless forms, various granulations and made out of diverse materials available today. The number of providers is currently limitless, which makes it all the more difficult for you as a user to decide which are the ideal solutions in this area. As an experienced specialist for all issues regarding podiatry and chiropody, SÜDA sees it as their responsibility to offer unique, perfect quality – also where filing machines are concerned. Therefore, we have added our own cutter program from the SÜDA brand. See the cutters on the following pages, clearly classified and the recommended application areas are explained.



Ceramic cutters

The perfect innovation for first class treatment



The SÜDA ceramic cutters are ideal for handling highly thickened nails, very hard skin and calluses and artificial nail materials under prosthetic or cosmetic corrections. In nail salon they are the ultimate instrument for a speedy removal of plastics, such as UV gel or acrylic. Both, users and their customers and patients report alike repeatedly by a gentle way of working that is felt to be particularly pleasant. The white, aesthetic appearance also emphasize the impression of a hygienic working environment.

The work of these instruments consists of a high-performance all-ceramic also used in orthopedics on ball heads for hip joints or in the dental field for implant abutments. The biocompatible material is highly resistant to pressure, wear and chemical disinfectants and thus extremely durable.



- Extremely good cutting properties
- Impressively smooth surfaces
- No smearing and clogging of the gear teeth
- Perfect smoothness due to reduced vibration
- Low heat than metal drills
- High resistance to disinfectants
- Sterilizable by autoclaving in hot-air sterilizer

ACR-GEAR TEETH



28302.060	Ceramic cutter (ACR) pear	↻ opt. 20.000 rpm
28304.040	Ceramic cutter (ACR) pyramid	↻ opt. 20.000 rpm

Excellent for treatment of thickened nails of wood or very hard mycotic nails. Particularly profitable staggered toothting, also for artificial nails materials such as gel or acrylic. Leaves a very smooth surface which does not have to be reworked.

GSQ-GEAR TEETH



28307.060	Ceramic cutter (GSQ) pear	↻ opt. 20.000 rpm
28301.040	Ceramic cutter (GSQ) pyramid	↻ opt. 20.000 rpm
28305.023	Ceramic cutter (GSQ) pyramid	↻ opt. 20.000 rpm

Universally suitable for processing heavily callused skin, calluses, thickened nails or nail art materials. Very cut joyful coarse cross-cut teeth. Depending on the shape ideal for a speedy nail removal, the selective callus removal or calluses in the furrow.

EQ-GEAR TEETH



28303.060	Ceramic cutter (EQ) pear	↻ opt. 20.000 rpm
-----------	--------------------------	-------------------

Combined toothting for strong removal and finishing with only one instrument. Medium staggered tooth system and fine cross-cut toothting, suitable for removing and reworking thickened natural nails and nail art materials (UV gel).

A-GEAR TEETH



28306.023	Ceramic cutter (A) ball	↻ opt. 6.000 rpm
-----------	-------------------------	------------------

Particularly sharp teeth for selective processing of rough hornification on the skin as well as for milling corns and fungal infections at the nail.



TIP

28309	Ceramic cutter set, 7 pieces
-------	------------------------------

Ceramic mills should after disinfection if required only be cleaned with brushes with metal-free bristles (for example, item no. 2643). Otherwise, metal debris may resign in the pores of the ceramic material, which leads to a blackening of the ceramic.

For best results, we ask you to note the recommended rotational speed.

Tungsten carbide cutters

No chance for mycoses - the "fungal nail cutters" of SÜDA

Provided with a special coating, these cutters are specifically designed for an effective but gentle removal of nail mycosis by rough and brittle nails. In addition, they are outstandingly suitable for selective removal softer corneal areas (eg. after a foot bath) or gentle milling of fissures. Depending on the application we recommend for these instruments a rotation speed of 15.000-20.000 rpm / min.



- Specially developed cutter construction for the area of chiropody / podiatry
- Less heat development through better surface cooling
- Can be very well used even at low to medium rotational speed in the extraction technology.
- Extremely low clogging of the chip space, for better removal and less cleaning.
- The material carbide correspondingly long service life due to low wear.
- This distinct advantage compared to diamond cutters lies in application as well as economically.



TIP	
28204	Cutter set carbide metal gold, 4 pieces in cutter stand aluminium small
Are you still working with the traditional RF-steel cutters? Then we recommend you to switch to carbide tools for more work comfort, better performance and longevity. Our cutter in the gear tooth system GP, P, and GSQ FSQ are well appropriate.	

GP-GEAR TEETH



28201.040	Tungsten carbide cutters (GP) pyramid	☺ opt. 15.000 rpm
Coarse igear tooth forming for the removal of rough and brittle, mycotic nails.		

P-GEAR TEETH



28202.031	Tungsten carbide cutters (P) pyramid	☺ opt. 15.000 rpm
28200.040	Tungsten carbide cutters (P) pyramid	☺ opt. 15.000 rpm
28203.070	Tungsten carbide cutters (P) pyramid	☺ opt. 15.000 rpm
Coarse toothing for the removal of rough and brittle, fungal nails.		

GSQ-GEAR TEETH



28350.060	Tungsten carbide cutters (GSQ) pyramid	☺ opt. 20.000 rpm
28351.060	Tungsten carbide cutters (GSQ) pear	☺ opt. 20.000 rpm
28352.040	Tungsten carbide cutters (GSQ) pyramid	☺ opt. 20.000 rpm
28352.070	Tungsten carbide cutters (GSQ) pyramid	☺ opt. 20.000 rpm
28354.023	Tungsten carbide cutters (GSQ) Cone	☺ opt. 25.000 rpm
28356.023	Tungsten carbide cutters (GSQ) Cone	☺ opt. 25.000 rpm
28357.016	Tungsten carbide cutters (GSQ) Cone	☺ opt. 25.000 rpm

Rough cut joyful cross-cut toothing (GSQ) for removing heavy skin and calluses. Removal of thickened nails and nail art material. Treat of plastics orthoses. Benefits: Effective cutting without smearing. High material removal rate. Gentle or vigorous milling depending on pressure.

FSQ-GEAR TEETH



28340.060	Tungsten carbide cutters (FSQ) pyramid	☺ opt. 20.000 rpm
28341.060	Tungsten carbide cutters (FSQ) pear	☺ opt. 20.000 rpm
28342.040	Tungsten carbide cutters (FSQ) pyramid	☺ opt. 20.000 rpm
28342.070	Tungsten carbide cutters (FSQ) pyramid	☺ opt. 20.000 rpm
28343.023	Tungsten carbide cutters (FSQ) chuck cone	☺ opt. 25.000 rpm
28344.023	Tungsten carbide cutters (FSQ) chuck cone	☺ opt. 25.000 rpm
28345.014	Tungsten carbide cutters (FSQ) drop	☺ opt. 25.000 rpm

Fine cut joyful cross-cut toothing (FSQ) for the final nail and skin smoothing. Fine work on nails and nail art material. Modeling of orthoses plastics. Benefits: Effective cutting without smearing. Particularly smooth milling at low pressure. Little heat generation by small cutting segments. Achieving optimum surface quality.

For best results, we ask you to note the recommended rotational speed.

SGE-GEAR TEETH



2820.060	Tungsten carbide cutters (SGE) pear	↻ opt. 15.000 rpm
2821.050	Tungsten carbide cutters (SGE) pyramid	↻ opt. 15.000 rpm

Super Coarse staggered toothting (SGE) for processing greatly thickened nails (woode nails, grypotische nails). Advantages: High cutting performance for particularly efficient removal of nail substance. Quiet operation for ergonomic working.

GE-GEAR TEETH



2835.060	Tungsten carbide cutters (GE) pear	↻ opt. 20.000 rpm
2836.050	Tungsten carbide cutters (GE) pyramide	↻ opt. 20.000 rpm

Coarse staggered toothting (GE) for the removal of thickened nails and gross shaping. Processing of artificial nail mass. Advantages: High cutting performance for a speedy removal of nail substance. Quiet operation for ergonomic working.

E-GEAR TEETH



2822.060	Tungsten carbide cutters (E) pyramid	↻ opt. 20.000 rpm
2823.040	Tungsten carbide cutters (E) pear	↻ opt. 20.000 rpm
2824.060	Tungsten carbide cutters (E) pear	↻ opt. 20.000 rpm
2825.031	Tungsten carbide cutters (E) pyramid	↻ opt. 20.000 rpm
2826.060	Tungsten carbide cutters (E) pyramid	↻ opt. 20.000 rpm
2827.060	Tungsten carbide cutters (E) roll	↻ opt. 20.000 rpm

Medium Cross cut (E) for universal nail processing and shaping. Processing of artificial nail mass. Advantages: Can be used universally. Quiet operation for ergonomic working.


EF-GEAR TEETH



2830.040	Tungsten carbide cutters (EF) pyramid	↻ opt. 20.000 rpm
2830.060	Tungsten carbide cutters (EF) pyramid	↻ opt. 20.000 rpm
2831.060	Tungsten carbide cutters (EF) pear	↻ opt. 20.000 rpm
2832.031	Tungsten carbide cutters (EF) pyramid	↻ opt. 20.000 rpm
2837.018	Tungsten carbide cutters (EF) chuck cone	↻ opt. 25.000 rpm
2837.023	Tungsten carbide cutters (EF) chuck cone	↻ opt. 25.000 rpm

Fine staggered toothting (EF) for final processing of nail surface. Nail smoothing to facilitate the subsequent polishing. Processing of artificial nail mass.

Advantages: Quiet operation for ergonomic working. Soft, pressure poor milling.

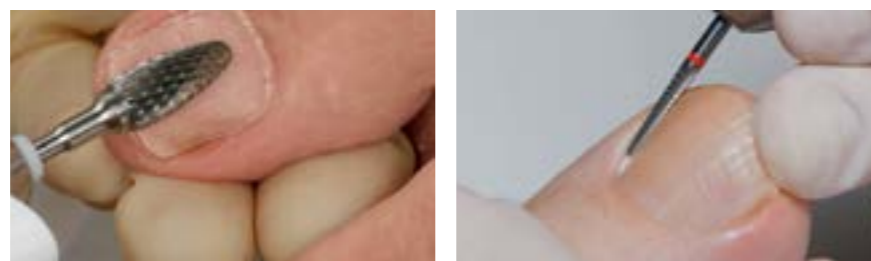

S-GEAR TEETH



2838.010	Tungsten carbide cutters (S) ball	↻ opt. 25.000 rpm
2838.014	Tungsten carbide cutters (S) ball	↻ opt. 25.000 rpm
2838.018	Tungsten carbide cutters (S) ball	↻ opt. 25.000 rpm
2838.023	Tungsten carbide cutters (S) ball	↻ opt. 25.000 rpm

Fast cutting teeth (S) for milling hardened skin and nail mycoses.

Removal of corns under the nail plate. Advantages: High cutting output with quiet operation. Long-lasting cutting sharpness.


STITCH CUTTER



2839	Stitch cutter with triangular working part	↻ opt. 6.000 rpm
------	--	------------------


For cutting and for processing of plastic orthoses.

Standard cutters used anti-clockwise.

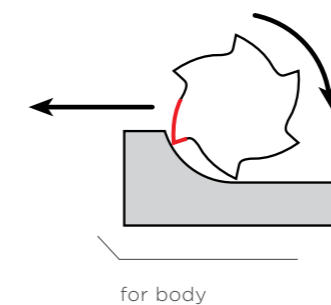

for body

Left-handed cutters used anti-clockwise.

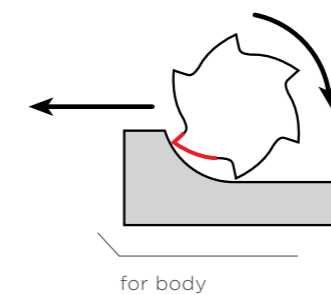

for body

We will not let you lie to the left - SÜDA cutter especially for lefthanded
For all cutters with a defined cutting it unfortunately is not enough to provide the pedicure device to anticlockwise. In order to work as a left-handed with these instruments properly, therefore, require special left-handed cutter having a corresponding toothting. Only by left-handed rotation movement of the device an effective removal by conventional milling in body direction is made possible.

L-GEAR TEETH



2836L014	Tungsten carbide cutters (L) ball	↻ opt. 25.000 rpm
2836L040	Tungsten carbide cutters (L) pyramid	↻ opt. 20.000 rpm
2836L060	Tungsten carbide cutters (L) pear	↻ opt. 20.000 rpm
2836L070	Tungsten carbide cutters (L) pyramid	↻ opt. 15.000 rpm

Cutters with blade geometry specifically for reverse operation. Conventional milling in body direction. Labelling with purple color ring. Shaft labelling "left/ right".

For best results, we ask you to note the recommended rotational speed.

Diamond grinders

The cylindrical shape is particularly suitable for smoothing and trimming of nails and for deburring the cut edges after cutting with a nail clipper. Also the processing of nail surfaces in the area of the nail wall and removing cuticles done with these tools easily.

In addition to the standard form, cylindrical grinder in special shapes are available which offer a lower risk of injury and therefore additional security by rounded edges or bevels. Therefore they are also recommended especially for the treatment of diabetic foot.

- 2876.055
Instrument with rounded edge
- Smoothing nails and cuticles
 - Minimized risk of injury
 - More safety by high risk patient



- 2874.051
Instrument is only top diamond coated
- Ideal for shortening nails
 - Adjacent skin is not ground down
 - Safety bevel on the front side edge



- 2874.055
Instrument is only laterally diamond coated
- Efficient grinding possible to the edge of the nail
 - Significantly minimized risk of injury
 - Safety bevel on the front side edge



TIP

When using diamond cutters for the rough removal of nails it is possible that it comes to unpleasant high heat generation due to the strong abrasion. Therefore we recommend you to use specific carbide cutters for this work, because less heat is generated by the cutting removal. In addition, you will benefit from the much longer life of these instruments.

MITTLERE KÖRNUNG 2ER-SET

2865.010	Diamond grinder medium ball	☺ opt. 40.000 rpm
2865.012	Diamond grinder medium ball	☺ opt. 40.000 rpm
2865.014	Diamond grinder medium ball	☺ opt. 40.000 rpm
2865.016	Diamond grinder medium ball	☺ opt. 40.000 rpm
2865.023	Diamond grinder medium ball	☺ opt. 40.000 rpm






MITTLERE KÖRNUNG

28661.033	Diamond grinder medium chuck cone	☺ opt. 30.000 rpm
2866.040	Diamond grinder medium chuck cone	☺ opt. 30.000 rpm
2867.033	Diamond grinder medium chuck cone	☺ opt. 30.000 rpm
2868.040	Diamond grinder medium chuck cone	☺ opt. 30.000 rpm
2870.016	Diamond grinder medium chuck cone	☺ opt. 35.000 rpm
2870.023	Diamond grinder medium chuck cone	☺ opt. 35.000 rpm
2871.016	Diamond grinder medium drop	☺ opt. 35.000 rpm
2872.023	Diamond grinder medium drop	☺ opt. 35.000 rpm
2873.060	Diamond grinder medium cone	☺ opt. 20.000 rpm
2879.047	Diamond grinder medium cone	☺ opt. 25.000 rpm
2879.060	Diamond grinder medium cone	☺ opt. 20.000 rpm
2874.051	Diamond grinder medium cylinder only top diamond coated	☺ opt. 20.000 rpm
2874.055	Diamond grinder medium cylinder only top diamond coated	☺ opt. 20.000 rpm
2875.037	Diamond grinder medium cylinder	☺ opt. 30.000 rpm
2876.055	Diamond grinder medium cylinder with rounded edge	☺ opt. 20.000 rpm
2877.055	Diamond grinder medium cylinder	☺ opt. 20.000 rpm
2878.100	Diamond grinder medium cylinder	☺ opt. 10.000 rpm

For the universal processing of corneal and nails. Smoothing after corneal ablation by rough cutter and grinder. Deburring the nail edges after cutting with the nipper. Processing of nail fold and cuticle.

For best results, we ask you to note the recommended rotational speed.

FINE GRAIN

	28531.012	Diamond grinder fine flame	☺ opt. 35.000 rpm
	28551.016	Diamond grinder fine chuck cone	☺ opt. 35.000 rpm
	28551.021	Diamond grinder fine chuck cone	☺ opt. 35.000 rpm
	2856.016	Diamond grinder fine drop	☺ opt. 35.000 rpm
	2857.023	Diamond grinder fine drop	☺ opt. 35.000 rpm
	2858.033	Diamond grinder fine chuck cone	☺ opt. 30.000 rpm
	2859.041	Diamond grinder fine cone	☺ opt. 25.000 rpm
	2860.060	Diamond grinder fine cone	☺ opt. 20.000 rpm
	2862.034	Diamond grinder fine cylinder	☺ opt. 30.000 rpm
	2862.051	Diamond grinder fine cylinder	☺ opt. 20.000 rpm






For final smoothing of nails before polishing.
Fine processing of nail edges and nail fold. Removing cuticles.

COARSE GRAIN

	28811.016	Diamond grinder coarse flame	☺ opt. 35.000 rpm
	28821.023	Diamond grinder coarse chuck cone	☺ opt. 35.000 rpm
	2883.042	Diamond grinder coarse chuck cone	☺ opt. 25.000 rpm
	2884.044	Diamond grinder coarse cone	☺ opt. 25.000 rpm
	2885.060	Diamond grinder coarse cone	☺ opt. 20.000 rpm
	2886.037	Diamond grinder coarse cylinder	☺ opt. 30.000 rpm
	2886.055	Diamond grinder coarse cylinder	☺ opt. 20.000 rpm

To remove calluses and for rough nail processing.

SUPER COARSE GRAIN

	28360.090	Diamond grinder super coarse pear	☺ opt. 15.000 rpm
	2887.060	Diamond grinder super coarse cone	☺ opt. 20.000 rpm
	2888.047	Diamond grinder super coarse cone	☺ opt. 25.000 rpm
	2888.060	Diamond grinder super coarse cone	☺ opt. 20.000 rpm
	2889.058	Diamond grinder super coarse cylinder	☺ opt. 20.000 rpm

For removal of rough skin. Grinding thickened nails.



EXTRA COARSE GRAIN



2890.060	Diamond grinder extra corase cone	☺ opt. 20.000 rpm
2890.080	Diamond grinder extra corase pear	☺ opt. 15.000 rpm

This extra coarse diamond grinders are ideal for removing very thick skin and calluses. The helical surface interruption ensures a reduced heat generation and an improved removal of the ablated cornea particles. Benefits: maintaining a high removal rate. Longer grinding greatly by lower tendency to smudging and clogging of the surface. Suitable for suction technology due to reduced heat generation. Optimal use by spray feed in wet technology devices.

HOLE DIAMOND



28375.100	Hole diamond medium	☺ opt. 15.000 rpm
28375.130	Hole diamond medium	☺ opt. 10.000 rpm
28376.100	Hole diamond coarse	☺ opt. 15.000 rpm
28376.130	Hole diamond coarse	☺ opt. 10.000 rpm
28377.100	Hole diamond super coarse	☺ opt. 15.000 rpm
28377.130	Hole diamond super coarse	☺ opt. 10.000 rpm

The perforated diamond grinders are ideal for the flat corneal processing and are especially popular because of their cap mold for many users. They are available in various sizes and grain sizes and therefore universally applicable depending on the intended removal. Advantages: Time-saving processing of larger corneal surface. Better heat dissipation through the interrupted with holes hollow body. Support the suction in the suction technology. Optimize instrument cooling in the wet technique.

HOLLOW DIAMOND



28370.090	Hollow diamond medium	☺ opt. 15.000 rpm
28370.115	Hollow diamond medium	☺ opt. 10.000 rpm
28371.090	Hollow diamond coarse	☺ opt. 15.000 rpm
28371.115	Hollow diamond coarse	☺ opt. 10.000 rpm
28372.090	Hollow diamond super coarse	☺ opt. 15.000 rpm
28372.115	Hollow diamond super coarse	☺ opt. 10.000 rpm

These diamond cutter are recommended universally for a flat corneal processing. The hollow body carries here again in a heat reduction on the instrument surface. Because of the hand piece towards closed mold these instruments, however, are primarily used for the suction technology as an optimum striking of the spray on the working part but in the wet technique can not be ensured in any case.

For best results, we ask you to note the recommended rotational speed.

RF-Steel cutters

RF-steel cutters are one-piece stainless steel construction and are counted among the classical used instruments that come in the field of foot care for many years. They are depending on the shape and cutting construction for all current work on the cornea, corns, skin and nails.

ROSE-HEAD BURS



2801.010	Rose-head burs RF, set of 2	U opt. 25.000 rpm
2801.014	Rose-head burs RF, set of 2	U opt. 20.000 rpm
2801.018	Rose-head burs RF, set of 2	U opt. 20.000 rpm
2801.023	Rose-head burs RF, set of 2	U opt. 15.000 rpm

Simply toothed, for milling corns and fungal infections at the nail.

ROSE-HEAD BURS



2802.010	Rose-head burs RF cross-cut, set of 2	U opt. 25.000 rpm
2802.014	Rose-head burs RF cross-cut, set of 2	U opt. 20.000 rpm
2802.018	Rose-head burs RF cross-cut, set of 2	U opt. 20.000 rpm
2802.023	Rose-head burs RF cross-cut, set of 2	U opt. 15.000 rpm

Cross-cut toothings, for milling corns and calluses on the skin.

FISSURE CUTTERS



2805.010	Fissure cutters RF, set of 2	U opt. 25.000 rpm
2805.012	Fissure cutters RF, set of 2	U opt. 20.000 rpm
2805.014	Fissure cutters RF, set of 2	U opt. 20.000 rpm
2805.016	Fissure cutters RF, set of 2	U opt. 15.000 rpm

Cross-cut toothings, for work on the nail fold and for milling hardened skin on skin.

TIP

Particularly when processing of hard and thick nails, subject to cutting RF steel mills due to the material to greater wear. We recommend for these applications, therefore the use of our much more durable carbide cutters.



CALLUS CUTTER

2806.050	Steel cutters RF roll, set of 2	U opt. 10.000 rpm
2806.060	Steel cutters RF roll, set of 2	U opt. 10.000 rpm
2807.040	Steel cutters RF oval, set of 2	U opt. 10.000 rpm
2809.050	Steel cutters RF cone, set of 2	U opt. 10.000 rpm
2809.060	Steel cutters RF cone, set of 2	U opt. 10.000 rpm
2810.100	Steel cutters RF roll	U opt. 8.000 rpm

Cross-cut toothings, for universal processing of hard skin and calluses.



NAIL CUTTER

2811.100	Steel cutters RF roll	U opt. 8.000 rpm
----------	-----------------------	------------------

Simply toothed, for the removal of thickened nails.



NAIL FOLD CLEANER

2808.012	Nail fold cleaner RF, set of 2	U opt. 8.000 rpm
----------	--------------------------------	------------------

Hexagon working part without sharp edges for gentle removal of deposits in the nail fold and excess cuticle.



HOLLOW CUTTERS WITH TOOTHED CUTTING

2815.018	Hollow cutter RF toothed	U opt. 6.000 rpm
2815.023	Hollow cutter RF toothed	U opt. 6.000 rpm
2815.027	Hollow cutter RF toothed	U opt. 6.000 rpm

For cutting of calluses and hard skin grooves.



HOLLOW MILLS WITH GROUND CUTTING

2816.018	Hallow cutter RF ground	U opt. 6.000 rpm
2816.023	Hallow cutter RF ground	U opt. 6.000 rpm
2816.027	Hallow cutter RF ground	U opt. 6.000 rpm

To remove deep-seated corns.

For best results, we ask you to note the recommended rotational speed.

Corundum grinder, Polishers and brushes

TIP
Corundum grinder are among the instruments of the lower price segment, are not medical devices and are subject to heavy wear. We therefore recommend alternatively to use diamond cutters.

CORUNDUM PINK



2892.1	Corundum grinder pink cylinder	☺ opt. 30.000 rpm
2892.2	Corundum grinder pink roll	☺ opt. 30.000 rpm
2892.3	Corundum grinder pink chuck cone	☺ opt. 30.000 rpm
2892.4	Corundum grinder pink cone	☺ opt. 30.000 rpm
2892.5	Corundum grinder pink drop	☺ opt. 30.000 rpm
2892.6	Corundum grinder pink chuck cone	☺ opt. 30.000 rpm
2892.7	Corundum grinder pink drop	☺ opt. 30.000 rpm

To smooth glass surfaces and deburring nail edges as well as for fine reworking of skin and nails.

CORUNDUM WHITE



2653	Corundum white cylinder, set of 2	☺ opt. 15.000 rpm
2653.1	Corundum white cylinder, set of 2	☺ opt. 30.000 rpm

For roughening of the nail plate in the B/S brace system and for grinding work in the nail prosthetics.

POLISHERS ANTHRACIT



28101.055	Prepolishers anthracite, set of 2	☺ opt. 10.000 rpm
28101.100	Prepolishers anthracite, set of 2	☺ opt. 8.000 rpm

For smoothing and pre-polishing of natural and artificial nails.

POLISHERS KHAKI



28102.055	Gloss polishers khaki, Set of 2	☺ opt. 10.000 rpm
28102.100	Gloss polishers khaki, Set of 2	☺ opt. 8.000 rpm

To gloss polishing of natural and artificial nails.

GOAT HAIR BRUSH



7345	Goat hair brush, set of 2	☺ opt. 6.000 rpm
------	---------------------------	------------------

For the final nail cleaning after treatment.

For best results, we ask you to note the recommended rotational speed.

SÜDA cap grinders

Quality at reasonable prices

Grinder caps are for decades the classic instrument for corneal treatment. They are still used in many foot care practices, as good removal results can be achieved with them. To ensure a hygienic work, pay attention to the intended use of the grinder caps as a disposable item as a disinfection of grinder caps is not possible. The grinders themselves can be disinfected in an ultrasonic bath and sterilizable in an autoclave at max. 121 ° C (gentle program for heat-sensitive instruments).



CAP GRINDERS

2591.10	Grinder 10 mm	1 pcs.
2591.110	Caps round 10 mm	Package of 10 pcs.
2591.1100	Caps round 10 mm	Package of 100 pcs.
2591.111	1 cap grinder round 10 mm made of rubber + 50 Caps	1 grinder + 50 Caps
2591.20	Grinder 13 mm	1 pcs.
2591.210	Caps round 13 mm	Package of 10 pcs.
2591.2100	Caps round 13 mm	Package of 100 pcs.
2591.211	1 cap grinder round 13 mm made of rubber + 50 Caps	1 grinder + 50 Caps
2591.30	Grinder 13 mm	1 pcs.
2591.310	Caps sharp 13 mm	Package of 10 pcs.
2591.3100	Caps sharp 13 mm	Package of 100 pcs.
2591.311	1 cap grinder sharp 13 mm made of rubber+ 50 caps	1 grinder + 50 Caps
2591.40	Grinder 5 mm 1 pcs.	1 pcs.
2591.410	Caps round 5 mm	Package of 10 pcs.
2591.4100	Caps round 5 mm	Package of 100 pcs.
2591.411	1 cap grinder round 5 mm made of rubber + 50 caps	1 grinder + 50 Caps
2591.50	Grinder 7 mm	1 pcs.
2591.510	Caps round 7 mm	Package of 10 pcs.
2591.5100	Caps round 7 mm	Package of 100 pcs.
2591.511	1 cap grinder round 7 mm made of rubber + 50 caps	1 grinder + 50 Caps

TIP
Cap grinders are not suitable due to the material for use in wet technique hand pieces. We recommend for the use of diamond cutters, z. B. hole diamond in cap shape.

Due to the occurrence of heat and to protect your hand pieces from excessive wear we recommend to operate the grinder cap with the following speeds:

- Max. 6.000 U / min. with caps with diameter 10 and 13mm
- Max. 8.000 U / min. with caps with diameter 5 and 7 mm

Standard-Grinder caps

In proven quality



GRINDER CAPS ROUND

2618	Grinder 13 R, Ø 13 mm	↻ opt. 6.000 rpm
7307.001	Grinder caps 13 R, 10 pcs.	
7307.0010	Grinder caps 13 R, 100 pcs.	
2612	Grinder 10 R, Ø 10 mm	↻ opt. 8.000 rpm
7305.001	Grinder caps 10 R, 10 pcs.	
7305.0010	Grinder caps 10 R, 100 pcs.	
2621	Grinder 7 R, Ø 7 mm	↻ opt. 8.000 rpm
7308.001	Grinder caps 7 R, 10 pcs.	
7308.0010	Grinder caps 7 R, 100 pcs.	
2620	Grinder 5 R, Ø 5 mm	↻ opt. 8.000 rpm
7309.001	Grinder caps 5 R, 10 pcs.	
7309.0010	Grinder caps 5 R, 100 pcs.	



GRINDER CAPS SHARP

2619	Grinder 13 S, Ø 13 mm	↻ opt. 6.000 rpm
2616.1	Grinder caps 13 S, 10 pcs.	
2616.2	Grinder caps 13 S, 100 pcs.	
2608	Grinder 7 S, Ø 7 mm	↻ opt. 8.000 rpm
2608.1	Grinder caps 7 S, 10 pcs.	
2608.2	Grinder caps 7 S, 100 pcs.	

BOX OF CAP GRINDERS

26010	Box of cap grinders, 14 pcs.
-------	------------------------------

Content: 1 × 2618 grinder 13 R, 1 × 2612 grinder 10 R, 1 × 2621 grinder 7 R, 1 × 2620 grinder 5 R, 1 × 2619 grinder 13, 1 × 2608 grinder 7 S, 2 × 7307,001 grinder cap 13 R, 1 × 7305,001 grinder cap 10 R, 1 × 7308,001 grinder cap 7 R, 1 × 7309,001 grinder cap 5 R, 1 × 2616.1 grinder cap 13 S, 1 × 2608.1 grinder cap 7 S



For best results, we ask you to note the recommended rotational speed.

Cutter sets



28137 Cutter set suction technique

SCOPE OF DELIVERY

- 2802.010 Round bur RF cross-cut, set of 2
- 2802.014 Round bur RF cross-cut, set of 2
- 2805.010 Fissures cutter RF, Set of 2
- 2806.050 Steel cutter RF roll, set of 2
- 2815.023 Hollow cutter RF toothed
- 2808.012 Furrow cleaner RF, Set of 2
- 2876.055 Diamond grinder medium cylinder
- 2859.041 Diamond grinder fine cone
- 2591.20 Grinder round 13 mm
- 2591.210 Grinder caps round 13 mm, 10 pieces
- 2591.40 Grinder round 5 mm
- 2591.410 Grinder caps round 5 mm, 10 pieces
- 28801 Cutter box made of plastic



28138 Cutter set callused skin

SCOPE OF DELIVERY

- 2830.060 Carbide cutters (EF) pyramid
- 2888.047 Diamond grinder super coarse cone
- 2888.060 Diamond grinder super coarse cone
- 28376.100 Hole diamond rough
- 28377.100 Hole diamond super coarse
- 28801 Cutter box made of plastic



28135 Cutter set wet technique

SCOPE OF DELIVERY

- 2865.010 Diamond grinder medium ball, set of 2
- 2865.012 Diamond grinder medium ball, set of 2
- 28551.016 Diamond grinder fine cone
- 2888.047 Diamond grinder super coarse cone
- 2888.060 Diamond grinder super coarse cone
- 2859.041 Diamond grinder fine cone
- 2815.023 Hollow cutter RF toothed
- 28301.040 Ceramic cutter (GSQ) pyramid
- 28377.100 Hole diamond super coarse
- 28801 Cutter box made of plastic



28136 Cutter set nails

SCOPE OF DELIVERY

- 28301.040 Ceramic cutter (GSQ) pyramid
- 2865.010 Diamond grinder medium ball, set of 2
- 2876.055 Diamond grinder medium cylinder
- 2870.016 Diamond grinder medium cone
- 28551.016 Diamond grinder fine cone
- 2859.041 Diamond grinder fine cone
- 2856.016 Diamond grinder fine drop
- 28801 Cutter box made of plastic

Turbine instruments

We bring you up to speed SÜDA instruments for high performance turbine

SÜDA offers a selection of diamond cutters in very conventional forms which are suitable for all common work.

The FG-instruments (= friction grip) with a shaft diameter of 1.6 mm are provided specifically for use in turbine handpieces of dental units with compressed air technology.



ATTENTION!
Shaft diameter 1.6 mm suitable for turbine handpieces.



28124.009	Diamond grinder FG ball	☺ opt. 400.000 rpm
28125.012	Diamond grinder FG ball	☺ opt. 400.000 rpm
28120.014	Diamond grinder FG ball	☺ opt. 400.000 rpm
28128.016	Diamond grinder FG ball	☺ opt. 400.000 rpm
28123.014	Diamond grinder FG cone	☺ opt. 400.000 rpm
28126.014	Diamond grinder FG trapeze	☺ opt. 400.000 rpm
28122.014	Diamond grinder FG trapeze	☺ opt. 400.000 rpm
28127.033	Diamond grinder FG buoy	☺ opt. 100.000 rpm

For best results, we ask you to note the recommended rotational speed.

Cutter and grinder accessories



CUTTER STAND SMALL

28800 Cutter stand small

Made of anodized aluminum for admission of 15 instruments. Can be disinfected and sterilized up to 180 ° C. Cutters not included.



CUTTER STAND MEDIUM

28803 Cutter stand medium

Made of anodized aluminum for admission of 23 instruments. Can be disinfected and sterilized up to 180 ° C. Cutters not included.



CUTTER STAND LARGE

28805 Cutter stand large

Made of anodized aluminum for admission of 23 instruments. Can be disinfected and sterilized up to 180 ° C. Cutters not included.



CUTTER STAND MADE OF PLASTIC

28801 Cutter stand made of plastic

With transparent cover for admission of 10 instruments.

28802 Cutter stand made of plastic

With transparent cover for admission of 35 instruments.

8311.5 Cutter stand made of plastic

With holder for brass brush in the lid for admission of 30 instruments. Without content.



CUTTER STAND TURBINE

28810 Cutter stands turbine

Made of anodized aluminum for admission of 10 FG instruments with 1.6 mm shaft diameter. Can be disinfected and sterilized up to 180° C



CLEANING

2640 Brass disc brush

For cleaning cutters.
Attention! Not recommended for ceramic cutters or diamond grinders!



2643 Cleaning brush



2643.1 Replacement brushes

With nylon bristles, sterilizable. Suitable for ceramic cutter and diamond grinder.



INSTRUMENT DISINFECTION

2644 DC-1 concentrate, 1 Liter

For disinfection and cleaning of rotary instruments.
1 liter makes 100 liters of solution (1%), in an ultrasonic bath 50 liters (2%).



CUTTER CROSS TWEEZERS

2176 Cutter cross tweezers

For better removal of the cutter from the handpiece or the disinfecting solution.



FRÄSER-PRÜFLEHRE

2695 Cutter check gauge

To check deformations of cutter shafts.
With this tool you can check the concentricity of your cutter.
Can the cutter shank insert smoothly into the bore of the check gauge, a concentricity is ensured. If this is not or only partially the case, you should sort the affected cutter, in order to avoid damage to your handpiece.

General use and safety for the application of rotary pedicure instruments

AREA OF APPLICATION

The general operating and safety instructions listed here apply to all SÜDA cutters and attention should be paid.

PROPERLY APPLIED

- It is to be sure to use only technically and hygienically perfect, maintained and cleaned handpieces.
- The instruments must be inserted as far as possible into the handpiece.
- The instruments are to bring up to speed before applying to the object.
- Tilting or levering the instruments increases the fracture and injury and should be avoided.
- Because of dust due to milling and grinding it is recommended to work with a suction or wet technology device. Also a mouthguard should be worn at work.
- To avoid thermal damage caused by rotary instruments at low speed or with adequate cooling.
- Excessive pressing force must be avoided to minimize risk of injury and high heat generation.
- These on the labels recommended maximum speed must not be exceeded.
- Improper use leads to poor work results and may cause injuries.

SPEED RECOMMENDATIONS

- In case of grinding instruments (diamond grinder, corundum grinder, cap grinder) following rule can be applied: The greater the grinding head the lower the speed to be selected.
- Usually a range between 6.000-40.000 rpm / min. is provided.- In cutting instruments (ceramic cutters, carbide cutters, carbide cutters RF) the size of the milling head plays a subordinate role. Here, it de-

pends on the material and the teeth of the cutter.

- For ceramic and carbide cutters, revolutions per minute between 15.000 and 25.000 are required, for RF steel cutters usually somewhat less.
- Recommendations for optimal speed, see the catalog respectively at the corresponding item number. This is intended to serve only as indication and may vary slightly depending on the user and treatment case.
- The maximum speed must not be exceeded for safety reasons. Please take this information from the label on the packaging of the instrument.

SELECTION OF USED ROTATING INSTRUMENTS

- Broken off and shapeless cutting grinders cause vibrations and lead to rough surfaces and increase temperature development.
- Blank area on diamond cutters indicate lack of abrasive grain and can be evidence of blunt instruments. These cutters must be sorted out because you can cause increased heat generation on the skin or the nail.
- Deformed or not round burs to be sorted out immediately. They destroy the chucks of the handpieces and can lead to a violation of the device.
- If a cutter fall down so it shall be tested on concentricity before the next use. For this purpose there is a cutter test gauge (Item no. 2695).

INDICATIONS FOR STORING AND HANDLING

- All SÜDA cutters are non sterile packaged. Therefore, it is advised to disinfect prior to first use the cutter properly and sterilize. Thus is prevented that any harmful germs are brought to the skin.
- To disinfect the cutters are fully covered in a disinfectant bath with VAH-listing (e.g. care & clean instrument disinfection 9320.100).

- The exposure of each agent must be observed. But the time also should not be exceeded significantly, since otherwise there is a risk that the cutters corrode.
- After disinfection remove cutters from the bath and rinse thoroughly under running water. Remove any adhering particles with a cleaning brush while continuously rotating. Check instruments with cavity carefully for particles. Repeat this procedure until no more contamination visible.
- Then dry completely the cutters.
- Undergo a visual examination of purified instruments and sort out if damaged.
- Optionally place dry cutters in dry heat sterilization or in the autoclave. When autoclave sterilization it is recommended to weld to designated sterile bags (item no. 8235.206). To ensure that the instrument remains sterile until the next use.
- Even with the heat sterilization, the bagging with temperature hose (Item no. 1412.1) is recommended.
- Ideally, store the instruments in a designated cutter stand (Item 28800/28805). This stand then put in a plastic box for storage. Ideal for mobile use.
- Storage should be protected against dust-, moisture- and recontamination.
- Don't store nail polisher in a UV-protected tray and store against UV.

INDICATIONS FOR MATERIAL COMPATIBILITY

- Basically, only disinfectants with VAH listing should be used.
- Also polishers and corundum grinders should be disinfected with suitable non-alcoholic products such as Carborundol (Article 4636).
- Ceramic cutters, carbide cutters, diamond cutters and RF-steel cutters are suitable for sterilization up to 180 ° C.
- Corundum grinder brushes, polishers and cap carriers CANNOT be sterilized.

PRODUCT LIABILITY

- Any damage arising from non-compliance with the instruction manual leads to the exclusion of liability by the company IONTO Health & Beauty GmbH. This is particularly the case with incorrect use of rotary instruments by the user.
- Complaints of cutters must be made in writing form at least one week after receipt of goods to the customer. Any subsequent complaints can not be recognized. The IONTO Health & Beauty GmbH reserves to examine the complaining material and where appropriate to provide compensation.

Explanation of symbols on the packaging

- optimum revolutions per min (rpm)
- max. allowed revolutions per min (rpm)
- Packaging unit
- Purchase order number
- CE mark for medical devices
- Attention should be paid of enclosed documents
- ISO-number
- Lot number
- For single use only

Nippers and scissors



Naturally stainless

Durable and easy-care

Resistant to disinfectants

Heat resistant and can be sterilized up to 180° C

Excellent cutting quality

Consistently accurate (no joint play) even after many applications

Besides the right equipment for the practice or mobile use well-functioning tools is a prerequisite for successful treatment. After all, you get the customers and patients very close with your instruments.

So that you always have the right thing at hand, you will find at SÜDA a wide range of instruments, such as corner clamps in five different cutting variants from strong to extra fine. Our two product lines "Premium" and "Excalibur" leave it for different needs nothing to be desired.

Medi-blades and nail splitters



2210 MEDI-blade holder

Made of stainless steel. Easy change of blade, sterilizable.



MEDI-blades careBASIC, 100 pieces

- 2166.1 Size 1
- 2166.2 Size 2
- 2166.3 Size 3
- 2166.4 Size 4



MEDI-blades careBASIC, 100 pieces

- 2166.5 Size 5
- 2166.6 Size 6
- 2166.8 Size 8
- 2166.10 Size 10

Precisely ground hollow-gouge blades made of high-quality carbon steel. For precise work with horns on toes, fore-foot and heels as well as in the nail fold and intermediate toe area. Individually sterile packed.



2211.6 Holder for nail splitter blades

For all heavy work around the nail fold.



2212.6 Nail splitter blade 06, 1 piece

Nail splitter blades 06, 10 pieces

Non-sterile packed.

Scalpels



2169 Scalpel holder No. 3

Made of stainless steel, easy change of blade, sterilizable.



Aesculap-blades, 100 pieces

- 2164.10 Size 10
- 2164.11 Size 11
- 2164.13 Size 13
- 2164.15 Size 15

Individually sterile packed.
Suitable for scalpel holder No. 3 (item 2169)



Scalpel blades careBASIC, 100 pieces

- 2167.10 Size 10
- 2167.11 Size 11
- 2167.13 Size 13
- 2167.15 Size 15

Individually sterile packed.
Suitable for scalpel holder No. 3 (item 2169)



2171 Scalpel holder No. 4

Made of stainless steel, easy change of blade, sterilizable.



Aesculap-blades, 100 pieces

- 2165.18 Size 18
- 2165.19 Size 19
- 2165.20 Size 20
- 2165.21 Size 21

Individually sterile packed.
Suitable for scalpel holder No. 4 (item 2171)



Scalpel blades careBASIC, 100 pieces

- 2167.18 Size 18
- 2167.19 Size 19
- 2167.20 Size 20
- 2167.21 Size 21

Individually sterile packed.
Suitable for scalpel holder No. 4 (item 2171)



2178 OR-scalpel holder

Made of stainless steel, simple blade replacement by OR technique, sterilized.



OR-blades, 10 pieces

OR-blades, 100 pieces

- 2178.00 Size 0
- 2178.01 Size 1
- 2178.02 Size 2
- 2178.03 Size 3
- 2178.10 Size 00

Non-sterile packed



2172.1 Blade puller Scalpel-Ex

Made of stainless steel. Disinfecting and sterilizing. For your safety when changing blade! Essential accessory for preventing accidents.
Note: Only for blade holder No. 3 and 4.

APPLICATION

Insert scalpel in the blade puller until it stops, puller squeezing until the blade bends upward, then pull out the scalpel holder - ready!



ACCESSORIES

2172 Collection box for blades

For safe, hygienic storage and disposal of used scalpel blades and needles.

7110.5 Multi-Safe sani

Handy waste box for scalpels and other disposable blades. With universal puller / dresser. Filling volume: 0.2 liters. Similar to the picture.

Wave pedicure nippers

Innovative forms in top quality made in Solingen

Surely you set a high value when selecting your instruments not only on impeccable quality or high-quality appearance but also on functionality and a comfortable grip. Then you our WAVE nippers are perfectly adapted. The Solinger quality guarantees an excellent finish and impresses with its pronounced handiness. The nippers are sometimes significantly lighter than comparable tools and allow you by its slender shape a less tiring work. In order to guarantee you to receive the original quality, the WAVE instruments were provided with a design protection of the European Office for Harmonisation in Alicante. This prevents that you are disappointed by cheap imitations.



WAVE HEAD CUTTER STANDARD

2300 WAVE Head cutter standard, 14 cm

Classical form for sided section. Easier work through the lean WAVE form.



WAVE CUTTER HALF CUTTER

2300.1 WAVE Cutter half cutter, 14 cm

Half-page head cutter for unilateral work on the nail. By WAVE form ideal handling.



WAVE CORNER NIPPERS STANDARD

2303.11 WAVE Corner nippers standard, 11,5 cm

2303.13 WAVE Corner nippers standard, 13 cm

Corner nippers with stable standard cuts with rounded cutting jaws.



WAVE CORNER NIPPERS FLAME POINTED

2304.11 WAVE Corner nippers flame pointed, 11,5 cm

2304.13 WAVE Corner nippers flame pointed, 13 cm

Corner nippers with side cutting jaws for deeper nails.



WAVE NAIL NIPPER CURVED

2301 WAVE Nail nippers curved, 14 cm

Nail nipper conceived as a side cutter. Perfect ergonomics by Wave form.



WAVE NAIL NIPPER CURVED

2301.1 WAVE Nail nippers curved, 13 cm

Light version without lock. Suitable for small and medium-sized nails.



WAVE SKIN NIPPER

2302 WAVE Skin nipper, 10 cm

For precise removal of cuticles. Cutting area of 6 mm.

Wave Safty Pedicure nippers

Especially for the treatment of patients at risk, we recommend the WAVE safety nippers. By blunt rounded tips a possible risk of injury is significantly reduced. This advantage is of course otherwise estimated in daily work.



Suitable for the treatment of diabetics.



WAVE SAFETY CORNER NIPPERS FLAME POINTED

2306.11 WAVE safety corner nippers flame pointed, 11,5 cm

2306.13 WAVE safety corner nippers flame pointed, 13 cm

Corner pliers with side tapered cutting and blunt tips. Suitable for deeper nails especially for patients at risk.

WAVE SAFETY HEAD CUTTER STANDARD

2305 WAVE safety Head cutter standard, 14 cm

Classic shape with blunt tips sided injury-free cutting.

Premium-Instruments




High grade stainless steel in Solingen quality, manufactured on precision machines by toolmakers who know their trade: So once the legendary reputation of German workmanship was founded. And so the instruments of PREMIUM class are still manufactured in the house in Solingen EKS®. Solid instruments, which have grown in materials and workmanship highest demands, long stay functionally reliable and therefore rarely need to be replaced.



PREMIUM HEAD CUTTER STANDARD

2111 PREMIUM Head cutter standard, 14 cm

Instrument for precise machining of nails. Classical form sided section. Sectional area of 17 mm.

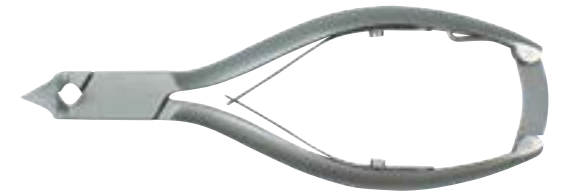

PREMIUM HEAD CUTTER FLAME POINTED

2111.2 PREMIUM Head cutter flame pointed, 14 cm

Bilateral head cutter with a tapered, pointed cutting. Also suitable for deeper nails. Sectional area of 20 mm.



PREMIUM HEAD CUTTER HALF CUT

2111.6 PREMIUM Head cutter half cut, 14 cm

Instrument for precise machining of nails. Half-sections for single-sided section. Sectional area of 15 mm.



PREMIUM NAIL NIPPERS CURVED

2110 PREMIUM nail nippers curved, 13 cm

Stable nail nippers with buffer spring. Also suitable for very thick nails. Sectional area of 19 mm.



PREMIUM NAIL NIPPERS CURVED

2110.1 PREMIUM nail nippers curved, 13 cm

Handy, small nail nippers with curved edge. Lightweight. Sectional area of 16 mm.



PREMIUM CORNER NIPPERS SLIM

2103 PREMIUM corner nippers slim, 10 cm

2103.1 PREMIUM corner nippers slim, 11,5 cm

2103.2 PREMIUM corner nippers slim, 13 cm

Slim corner nippers with relief grinding cutting for removing small nail corners.



PREMIUM CORNER NIPPERS SHARP POINTED

2115.2 PREMIUM corner nippers sharp pointed, 10 cm

2115 PREMIUM corner nippers sharp pointed, 11,5 cm

2115.1 PREMIUM corner nippers sharp pointed, 13 cm

Stable corner nippers with tapered jaws.



PREMIUM NAIL SPLITTER

2118.1 PREMIUM Nail splitter, 13 cm

Especially for processing of ingrown toenails. Safe-section through the two differently shaped jaw. Sectional area of 19 mm.



PREMIUM CORNER NIPPERS STANDARD

2113.2 PREMIUM corner nippers standard, 10 cm

2113 PREMIUM corner nippers standard, 11,5 cm

2113.1 PREMIUM corner nippers standard, 13 cm

Corner nippers with stable standard cuts with rounded cutting jaws.



PREMIUM CORNER NIPPERS FLAME POINTED

2116.2 PREMIUM corner nippers flame pointed, 10 cm

2116 PREMIUM corner nippers flame pointed, 11,5 cm

2116.1 PREMIUM corner nippers flame pointed, 13 cm

Corner nippers with side cutting jaws tapered for deeper nails.



PREMIUM CORNER NIPPERS SUPERFINE

2108 PREMIUM corner nippers superfine, 10 cm

2108.1 PREMIUM corner nippers superfine, 11,5 cm

2108.2 PREMIUM corner nippers superfine, 13 cm

Particularly narrow corner nippers with relief grinding cutting for removing nail smallest corners.



PREMIUM SKIN NIPPERS

2114 PREMIUM Skin nippers, 10 cm


Skin nipper for precise removal of cuticles. Sectional area of 7 mm.



PREMIUM SKIN NIPPERS ERGONOMIC

2090 PREMIUM skin nipper ergonomic, 10 cm

Skin nipper with angled cutting head. Relieves the wrist and makes work easier. Sectional area of 7 mm.

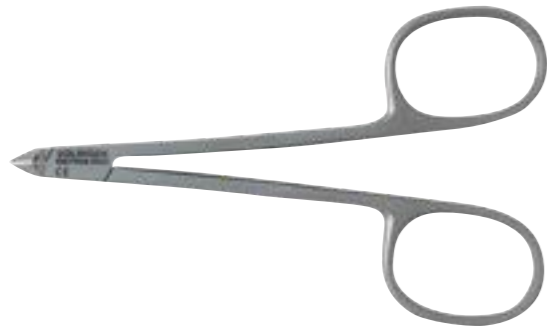

PREMIUM SKIN PINCH NIPPERS

2112 PREMIUM SKIN pinch nippers, 10 cm

Fine cuticle nipper with scissor grip. Ideal for removing overhanging cuticles and the processing of corns. Sectional area of 5 mm



PREMIUM NAIL SCISSORS CURVED

2119 PREMIUM skin scissors curved, 10 cm

Robust nail scissors with curved edge for trimming nails.



PREMIUM SKIN SCISSORS CURVED

2117 PREMIUM nail scissors curved, 10 cm

Classic cuticle scissors with curved blade for removing cuticles.



PREMIUM SKIN SCISSORS SPIRE POINTED

2117.1 PREMIUM skin scissors spire pointed, 10 cm

Extra fine cuticle scissors with spire pointed for difficult work on the cuticles.



PREMIUM BANDAGE SCISSORS

7030.3 PREMIUM bandage scissors pointed/ blunt, 14 cm

Scissors for bandages and dressings.



2118.8 PREMIUM double probe extra slim, 17 cm

At the suggestion of our customers: Special design with slim spatulas.



2118.5 PREMIUM double probe, 17 cm

Classic double probe in excellent quality.



2118.6 PREMIUM double probe with spoon, 17 cm

Double probe with a flat spoon.



2118.7 PREMIUM Nail corner file twice, 14 cm

Classical corner file straight / curved.



2098 PREMIUM Spatula blunt, 17 cm

Especially suitable for gentle scraping of the cornea, which for example was pretreated with SÜDAcare Kera special.



2105.3 PREMIUM splitter tweezer "Feilchenfeld", 8 cm

2105.4 PREMIUM splitter tweezer "Feilchenfeld", 11,5 cm

Splitter tweezer in mini and standard format, for precisely work. From two parts precision welded, therefore no lateral displacement possible.



2187.4 Skin curette Ø 4 mm

2187.7 Skin curette Ø 7 mm

Disposable instrument with sharply polished stainless steel loop. For the easy removal of keratinized or macerated skin, For example in the intermediate toe region. Individually sterile packed.

Excalibur-Instruments



The EXCALIBUR series instruments are distinguished by an exceptionally good price-performance ratio and therefore are also recommended especially for students or young professionals, where the question of cost often comes first. They are manufactured of stainless steel, very durable and naturally resistant to chemical disinfectants and sterilizing heat to 180° C.



EXCALIBUR HEAD CUTTER STANDARD

2128 EXCALIBUR head cutter standard, 14 cm

The classic nail nipper with blades on both sides for precise machining of nails. Sectional area 24 mm.



EXCALIBUR HEAD CUTTER HALF CUT

2126 EXCALIBUR head cutter half cut, 14 cm

Nail nipper with one side cutting head for unilateral cutting. Sectional area 12 mm.



EXCALIBUR NAIL NIPPERS STRAIGHT

2129.11 EXCALIBUR nail nippers straight, 13 cm

Nail nippers with straight cut. Good for strong nails, wooden nails, etc. sectional area 20 mm.



EXCALIBUR NAIL NIPPERS CURVED

2129 EXCALIBUR Nail nippers curved, 10 cm

Table nail nippers for stronger nails. Equipped with bow spring. Small and portable, ideal for home visit suitcase. Sectional area of 20 mm.



EXCALIBUR NAIL NIPPERS WITH TRANSMISSION

2134 EXCALIBUR nail nippers with transmission, 17 cm

Strong nail nippers with transmission. Cuts even thickest deformed nails effortlessly with little effort. Sectional area of 20 mm.



EXCALIBUR NAIL NIPPERS CURVED

2132 EXCALIBUR nail nippers curved, 12 cm

2133 EXCALIBUR nail nippers curved, 13 cm

Stable nail nippers with buffer spring. Also suitable for very thick nails. Sectional area of 16 mm.



EXCALIBUR NAIL NIPPERS CURVED

2132.1 EXCALIBUR nail nippers curved, 12 cm

Excellent small nail nippers with curved cutting. Fits comfortably in the hand, almost as a classic corner nippers. Sectional area of 14 mm.



EXCALIBUR NAIL SPLITTER

2118 EXCALIBUR Nail splitter, 13 cm

Special nippers for reliable processing of ingrown toenails. Safe-section through the two differently shaped jaws. Sectional area 18 mm.



EXCALIBUR CORNER NIPPERS STANDARD

2130 EXCALIBUR corner nippers standard, 10 cm

2131.11 EXCALIBUR corner nippers standard, 11,5 cm

2131.13 EXCALIBUR corner nippers standard, 13 cm

Corner nippers with stable standard cuts with rounded cutting jaws.
Sectional area in 2130 14 mm. Sectional area 2131.11 16 mm. Sectional area 2131.13 20 mm.



EXCALIBUR CORNER NIPPERS POINTED

2130.11 EXCALIBUR corner nippers pointed, 11,5 cm

2130.13 EXCALIBUR corner nippers pointed, 13 cm

Stable corner pliers with tapered jaws. Sectional area 11 cm 16 mm. Sectional area of 13 cm 20 mm.



EXCALIBUR SKIN NIPPERS

2120 EXCALIBUR Skin nippers, 10 cm

To remove the cuticles. The short beveled cutting also allows for difficult work at the nail edges.
Sectional area: 5 mm.



EXCALIBUR SKIN NIPPERS

2120.10 EXCALIBUR skin nippers, 10,5 cm

2120.15 EXCALIBUR skin nippers, 12 cm

Skin nippers for precise removal of cuticles. Model 2120.15 also suitable for small calluses. Sectional area 2120.10 6 mm. Sectional area 2120.15 10 mm.



EXCALIBUR CORNER NIPPERS FLAME POINTED

2135.11 EXCALIBUR corner nippers flame pointed, 11,5 cm

2135.13 EXCALIBUR corner nippers flame pointed, 13 cm

Corner nippers with side cutting jaws tapered for deeper nails. Sectional area 2135.11 18 mm. Sectional area 2135.13 20 mm.



EXCALIBUR CORNER NIPPERS FLAME POINTED MINI

2140.11 EXCALIBUR corner nippers flame pointed mini, 11,5 cm

2140.13 EXCALIBUR corner nippers flame pointed mini, 13 cm

Corner nippers with extra narrow top and slim cut surfaces. Optimal for ingrown, curled or hard to reach corners nail.



EXCALIBUR SKIN PINCH NIPPERS SHORT

2124 EXCALIBUR Skin pinch nippers short, 8 cm

Cuticle nippers with scissor grip for removing overhanging cuticles. With short legs.
This better handling with difficult work. Sectional area of 5 mm.



EXCALIBUR SKIN PINCH NIPPERS SHORT

2125 EXCALIBUR Skin pinch nippers short, 10 cm

Cuticle nippers with scissor grip for removing overhanging cuticles. Sectional area of 8 mm.



EXCALIBUR SKIN PINCH NIPPERS

2122 EXCALIBUR skin pinch nippers,
10 cm

Cuticle nippers with scissor grip for removing over-
hanging cuticles. Sectional area of 8 mm.



EXCALIBUR SKIN PINCH NIPPERS STRAIGHT

2123 EXCALIBUR skin pinch nippers
straight, 10 cm

Cuticle nippers with scissor grip for removing over-
hanging cuticles. With straight
blade. Also suitable for smaller nail corners.
Sectional area of 10 mm.



EXCALIBUR NAIL SCISSORS CURVED

2139 EXCALIBUR Nail scissors curved,
9 cm

Robust nail scissors with curved edge and sharp cut.
Also suitable for harder nails.



EXCALIBUR BANDAGE SCISSORS

2138 EXCALIBUR bandage scissors,
13 cm

Scissors for bandages and dressings.
Cutting sharp / blunt.



EXCALIBUR SPLITTER TWEEZERS »FEILCHENFELD«

2175 EXCALIBUR Splitter tweezers,
11 cm

2175.1 EXCALIBUR Splitter tweezers,
7,5 cm



2144.1 EXCALIBUR double probe,
17 cm

Classic double probe right / left with spatulas.



2189 EXCALIBUR double probe
with spoon, 17 cm

Double probe with a flat spoon.



2188 EXCALIBUR nail corner file
twice, 14.5 cm



2146 EXCALIBUR nail corner file
simple 15 cm



2142 EXCALIBUR double instrument
straight / angled, 16 cm



2143 EXCALIBUR double instrument
straight / curved, 14 cm



2147 EXCALIBUR cotton hook,
14,5 cm




2155 EXCALIBUR gouge,
14 cm

Instrument Sets

Premium

A two-part base set with all important tools for daily tasks at the foot. Made of stainless steel in Solingen premium quality.

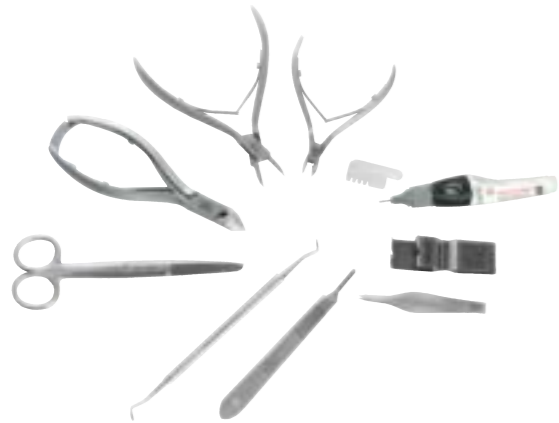

PREMIUM INSTRUMENT SET

7099.5 Set

SCOPE OF DELIVERY

2111	PREMIUM Head Cutters standard, 14 cm
2116.1	PREMIUM Corner nippers flame-pointed, 13 cm
2114	PREMIUM Skin nippers 10 cm, Cutting 7 mm
7030.3	PREMIUM Double probe extra narrow, 17 cm
2118.8	PREMIUM Doppelsonde extra schmal, 17 cm
2105.3	PREMIUM Splinter tweezers "Feilchenfeld" 8 cm
2169	Scalpel holder no. 3
2172.1	Blade pullers Scalpel-Ex
2196.5	SÜDA instrument care oil, 12 ml

Excalibur

A two-part base set with all important tools for daily tasks at the foot. Made of stainless steel in proven Excalibur Quality.



EXCALIBUR INSTRUMENT SET

7099.4 Set

SCOPE OF DELIVERY

2128	EXCALIBUR Head Cutters standard, 14 cm
2131.13	EXCALIBUR Corner nippers standard, 13 cm
2120	EXCALIBUR Skin nippers 10 cm, Cutting 5 mm
2138	EXCALIBUR Bandage scissors sharp / blunt, 13 cm
2189	EXCALIBUR Double probe with spoon, 17 cm
2175.1	EXCALIBUR Splinter tweezers Feilchenfeld, 7.5 cm
2169	Scalpel holder no. 3
2172.1	Blade pullers Scalpel-Ex
2196.5	SÜDA instrument care oil, 12 ml

Starter-Set callused skin

A two-part base set with all the main instruments for nail and callus treatment. Made of stainless steel in proven Excalibur Quality.



STARTER-SET CALLUSED SKIN


2179.7 Set

SCOPE OF DELIVERY

2128	EXCALIBUR Head Cutters standard, 14 cm
2135.11	EXCALIBUR Corner nippers flame-pointed, 11,5 cm
2120	EXCALIBUR Skin nippers 10 cm, Cutting 5 mm
2144.1	EXCALIBUR Double probe 17 cm
2175	EXCALIBUR splinter tweezers "Feilchenfeld" 11 cm
2141.5	SÜDA Corn planer Set
2178	OR scalpel holder
2178.00	OR Blades size 0, 10 pieces
2178.10	OR Blade size 00, 10 pieces

Starter-Set Nails

Eight-part base set with all major instruments specifically for nail treatment. Made of stainless steel in proven Excalibur Quality.

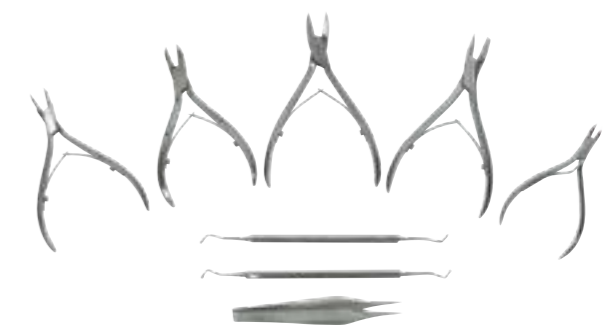

STARTER-SET NAILS

2179.8 Set

SCOPE OF DELIVERY

2130	EXCALIBUR corner nippers standard, 10 cm
2131.11	EXCALIBUR corner nippers standard, 11,5 cm
2135.11	EXCALIBUR corner nippers flame-pointed, 11,5 cm
2132.1	EXCALIBUR Nail nippers curved, 12 cm
2118	EXCALIBUR nail splitter 13 cm
2144.1	EXCALIBUR double probe 17 cm
2189	EXCALIBUR double probe with spoon, 17 cm
2175	EXCALIBUR splinter tweezers "Feilchenfeld" 11 cm

Corneal processing



CREDO-CORN PLANER

2141

Longitudinal slicer in metal version, sterilized, incl. blade and rasp clip.



CREDO CORN PLANER POPART

2141.6

Longitudinal planer with plastic handle incl. blade and rasp clip. Assorted colors, in blister. Colour not selectable.



SÜDA CORN PLANER SET

2141.5

Longitudinal planer in metal version, sterilized, incl. 10 blades.



CREDO BLADE PROFİ

2141.4 10 pieces

2141.40 100 pieces



CREDO DUOSOFT FUSSFEILE POPART

2180.2 1 piece

Stylish and practical foot file to remove excess corneal and delicate fine grinding of the feet. By interchangeable file blades, the file can be used repeatedly. Simply peel off old file blade and apply the new. Ideal for sales to end users. Assorted colors, in blister. Colour not selectable.



CREDO DUOSOFT SPARE RASP SURFACES

2180.3 2 coarse and 2 fine spare rasp surfaces

2180.30 fine, 25 pieces

2180.31 coarse, 25 pieces



CREDO CERAMIC SPARE RASP SURFACES POPART

2180.4 20 cm

For gentle callus removal on foot (Suitable for diabetics). With a fine and a coarse side. Long life and easy to clean. Ideal for sales to end users. Assorted colors in blister, color not selectable.



BIOSANA CALLUS FILE

2190 natur

2190.5 lacquered

From impregnated beech wood, FSC certified. File papers coarse / fine. Total length 25 cm. In the brightly lacquered files no color selection is possible.



CALLUS FILE FOOT FORM

2193 25 cm

Made of white plastic. File papers coarse / fine.



PUMICE STONE SPONGE

9844 1 piece

Proven scratch sponge for hands and feet. Removes calluses gently and naturally. Handy size (about the size of a 100 g soap bar). Colour not selectable.

UMICE STONE SPONGE ANTIBACTERIAL

9844.1 1 piece

For hygienic removal of the cornea. Size: L 9.9 x W 4.8 x H 2.5 cm. Colour not selectable.



MICROPLANE® CALLUS RASPER

2199 white, 1 piece

Dry skin and cornea are removed quickly and effortlessly. Application to the dry skin. By ergonomic shape easy access to all foot areas. Protective cap for safe storage and easy transport. Collecting cover for convenient and clean operation.



MICROPLANE® CALLUS SMOOTHER

2195 white, 1 piece

Smooths rough, dry skin. Recommended after Microplane® Callus file to achieve soft feet. Designed for a specific use on dry skin. Ideal to take along on trips. Protective cap for safe storage and easy transport.



Hygienic application.



Removing and disposing of the contaminated abrasive paper after the treatment.



Insert the steel file in a disinfectant bath, then rinse under running water, rinse, dry.



Load into the sterilizer



The fresh, self-adhesive surface is peeled from the carrier film.



New sandpaper easy stick back.

Callus file SK1

Unlike traditional corneal files, in which a commercial use already is difficult just because of the hygiene problem, the cornea file SK-1 allows unlimited use in constantly changing clientele. The solution to the hygiene problem is surprisingly simple, because the SK-1 consists of two parts.

Self-adhesive surface, which is cut to size and is re-stucked for each customer easily on the body.



Stainless steel body, which survives each required disinfection and sterilization routine without prejudice.
Length 25 cm

2191.3 Callus file SK-1

Steel file without abrasive surfaces.

RASP SURFACES COARSE FOR CALLUS FILE SK-1

2191.1 self-adhesive, 10 pieces

2191.10 self-adhesive, 100 pieces

RASP SURFACES FINE FOR CALLUS FILE SK-1

2191.2 self-adhesive, 10 pieces

2191.20 self-adhesive, 100 pieces

CALLUS FILE SK-1 SET

2191

Steel file incl. 10 rasp surfaces coarse and 10 rasp surfaces fine.

Hygiene

Clean with security

The complete hygiene process - from the thorough disinfection and cleaning to sterilization through to appropriate storage of instruments - is essential for each cosmetic institute. Compliance with the country-specific hygiene regulations play a major role. Only in this way you protect yourself and your customers from possible infections.

204 - 205

Ultrasonic device

206 - 215

Sterilization
technology

216 - 217

Hygiene sets

Ultrasonic devices

Fast disinfection and strong cleaning

Ultra 1431

For all medical instruments and cutlery.

1431 Ultra 1431

1431.001 Stainless steel basket



FEATURES

- Digital timer with 5 time settings (90/180/280/380/480s)
- Including deposit screen and lid
- New improved operator panel
- Innovative, patented design and ultra quiet operation

TECHNICAL SPECIFICATIONS

2.5 liter stainless steel tank, max. 2.2 liters
 Tank size: approx 246 x 146 x 80 mm
 Weight: 2.5 kg
 Cleaning power: 70 watts
 Ultrasonic Frequency: 42.000 Hz

DEGAS-FUNCTION

In the cleaning liquid, especially if it has been freshly filled, gas bubbles can reside. These prevent a clean formation of cavitation while conducting ultrasonic waves through the liquid and thus reduce the cleaning performance. By means of degas function, these gas bubbles are selectively discharged, the cleaning performance optimizes.

SWEEP-FUNCTION

Permanent change of the ultrasonic frequency to the better distribution of the cleaning performance in the liquid.

Ultra 1432

With switchable Degas function for medical instruments and cutlery.

1432 Ultra 1432

1432.001 Inner basket



FEATURES

- 1-30 min. timing
- Including inner basket and lid
- 2-colored LED display
- Switchable degas function
- Built 80 Watt heating (autom. Shutdown at approx 60° C)

TECHNICAL SPECIFICATIONS

2.0-liter stainless steel tank, max. 1.7L
 Tank size: approx 183 x 153 x 77 mm
 Weight: 2.2 kg
 Cleaning power: 70 watts
 Ultrasonic Frequency: 35.000 Hz

ELMASONIC S podo basic

Microprocessor controlled ultrasonic cleaning and sweep technology mark the latest technology. The sweep function is always active in this device. Thereby operation becomes easier.

8325.5 Elmasonic S podo basic



FEATURES

- Low-priced ultrasound device for the continuous or short operation
- Variable time level from 1 to 30 minutes.
- Housing and transducer tank entirely of stainless steel
- Automatically activated SWEEP function
- LED display for ultrasonic function

TECHNICAL SPECIFICATIONS

Power consumption: 230V / 50 - 60Hz
 Power: 30 Watt
 Sink Interior Dimensions: H 8,5 x W 19 x D 6 cm
 Device External dimensions: H 19 x W 23 x D 10.5 cm
 Weight: 1.8 kg
 Frequency: 37 kHz
 Volume: 800 ml

SCOPE OF DELIVERY

- Plastic lid, can be used as drip tray for the basket
- Plastic basket, 2 small parts mugs, tweezers
- elma clean 10 Cleaning liquid, 1.000 ml

SONIC 3D

Powerful ultrasound unit for professional use in podology practices as well as in larger foot care and cosmetic institutes. Designed for continuous operation. With its 3-liter volume, it also receives larger amounts of cleaning material, for example whole instrument sets.

1437 Sonic 3D



FEATURES

- DEGAS function
- LCD display
- Digital time and temperature setting 5-60 min.
- Instrument basket
- Pan and cover in stainless steel

TECHNICAL SPECIFICATIONS

Electrical connection: 200-240 V / 50 - 60 Hz
 Power: 120 watts
 Ultrasonic frequency: 40.000 Hz
 Volume: 3 l
 Tub interior dimensions: B 23.8 x H 10 x D 13.8 cm
 Unit dimensions: W 26.5 x H 22.5 x D 16.5 cm
 Weight: 3.4 kg



ATTENTION!
 Elma clean 10 can only be used as a cleaning additive - no disinfectant!

Cleaning and disinfecting equipment

MELAtherm®10

The cleaning and disinfection of instruments and cutters prior to sterilization is an integral part of the treatment process of the instruments in podology and foot care practices. The innovative concept of the thermal disinfecter MELAtherm®10 offers all the advantages for larger practices for a fast, safe and above all cost-saving processing.

FEATURES

- Active drying to protect the instruments from recontamination, corrosion and damage
- Dosing module and process media easily accessible and user-friendly in a media drawer
- Documentation possibilities via the integrated CF card slot or via two Ethernet interfaces to the network, PC or printer
- Practice-optimized processing chamber for short operating times, ergonomic working
- low consumption of process media and optimal cleaning results
- Standard width of 60 cm for uncomplicated installation in a hygiene line
- Alternatively, the device can be set up freely and fitted with a stainless steel cover plate
- Optionally available base cabinet for ergonomic operating height and additional storage space
- Basic basket (not included in the price - see accessories) optionally with and without injector rail for the assembly with turbine and handpiece attachments
- Comprehensive basket assortment for individual requirements available upon request
- Medical device Class IIb acc. Directive 93/42/EEC
- Cleaning and disinfecting unit according to Euro-norm DIN EN ISO 15883-1 / -2



1461.1	MELAtherm®10 DTA Built-in unit for three-phase connection
1461.2	MELAtherm®10 DTB Built-in unit for AC connection

SCOPE OF DELIVERY

- Thermodisinfector substructure without baskets
- Drawer with dosing unit for process media
- MEtherm starter set process media

TECHNICAL SPECIFICATIONS

Electrical connection DTA: 400 V, 50 Hz, 3 x 16 A
 Electrical connection DTB: 220-240 V, 50 Hz, 1 x 16 A
 Power consumption: 9.3 kW (DTA) / 3.3 kW (DTB)
 Dimensions Housing: W 59.8 x H 81.8 x D 67.8 cm
 Dimensions of washing chamber:
 B 46.5 x H 40.5 x D 44.8 cm



Accessories for MELAtherm®10

BASE BASKETS

1461.30	Basic basket without injector rail B 42 x H 9 x D 43 cm
1461.31	Basic basket with injector rail With 11 injector nozzles and clamping springs B 42 x H 12 x D 43 cm



INSTRUMENT BASKETS

1461.32	Instrument basket standard, 2/6 segment cross
1461.33	Shopping cart G, 2/6 segment along
1461.34	Instrument basket compact, 1/6 segment



FLEX BASKETS STACKABLE

1461.35	Flex basket 1, 1/6 segment
1461.36	Flex basket 2, 2/6 segment along the length
1461.37	Flex basket 3, 3/6 segment



PROCESS MEDIA

1461.38	MEtherm 50, canister 5 liters Mild alkaline, enzymatic cleaner Minimizes the formation of documents
1461.39	MEtherm 55, canister 5 liters neutralizer C (with citric acid) particularly gentle and material-compatible
1461.40	MEtherm 56, canister 5 liters neutralizer P (with phosphoric acid) for a particularly effective cleaning
1461.41	MEtherm 60, bottle 940 ml Liquid rinse aid to support spot-free drying



COVER PLATE AND SOCKET

1461.10	Stainless steel cover plate Dimensions: W 59.8 x H 2 x D 67.8 cm, weight: 6 kg
1461.11	Base cabinet for MELAtherm®10 Dimensions: W 60 x H 40 x D 60 cm, weight: 21 kg



Hot air sterilization

The heat sterilization is carried out at a temperature of 180 ° C by hot air, generated with concealed heating coils.

Hot air sterilizer MAXI

For all sterilizable instruments to 180° C.

1412 Hot air sterilizer MAXI



FEATURES

- Temperature range: from 0 to 200° C, continuously adjustable over temperature selector.
- Fully automatic operation of the entire sterilization process up to the power off.
- Practical use basket for the instruments.
- Indication of indoor temperature (only for MAXI-version).

SPECIFICATIONS FOR 1412 MAXI

Power requirements: 230 V / 50 - 60 Hz, 150 Watt
 Device dimensions (cm): H 37.5 x W 49 x D 19
 Basket dimensions (cm): H 22 x W 16.5 x D 6 cm
 Duration: 1,5 hours
 Weight: 8.5 kg
 Capacity of use basket: 2-4 sets of instruments
 approx. 2 - 4 sets of instruments

Hot air sterilizer MINI

For all sterilizable instruments to 180° C.

1411 Hot air sterilizer MINI



FEATURES

- Temperature range: from 0 to 200 ° C, continuously adjustable over temperature selector.
- Fully automatic operation of the entire sterilization process up to the power off.
- Practical use basket for the instruments.

SPECIFICATIONS FOR 1411 MINI

Power requirements: 230 V / 50 - 60 Hz, 150 Watt
 Device dimensions: H 31 x W 20 x D 19 cm
 Basket dimensions: H 24 x W 11 x D 2.7 cm
 Duration: 2 hours
 Weight: 3 kg
 Capacity of use basket: 1-2 sets of instruments

Recirculated air sterilization

In addition to the thermal distribution of the hot air is further circulated by a fan. Thereby, the chamber of the device is heated absolutely evenly and the sterilization process is complete after 60 minutes.

Recirculated air sterilizer STERIL UL

Für alle sterilisierbaren Instrumente bis 180 °C.

8235.10 STERIL UL



SCOPE OF DELIVERY

- Recirculated air sterilizer
- 1 jack for trays (item 8235.104)
- 1 set of holders for trays (item 8235.105)
- 2 aluminum trays (item 8235.101)

ACCESSORIES

- 8235.108 1 set of holders for 2 tablets incl. Lifter
- 8235.101 Aluminum tray
- 8235.104 Cartridge for recirculating air sterilizer
- 8235.105 Lift for cartridge / tablet
- 8235.102 1 set of brackets for 1 cassette / tray incl. Lifter
- 8235.103 Cassette lifter
- 1415 Cutter bowl stainless steel: W 9 x H 1,7 x D 7,2 cm

FEATURES

- Temperature range: from 60 to 180 ° C, can be selected individually.
- Automatic temperature control.
- Recirculation mode, guaranteed all-around sterilization.
- Timer with automatic shutdown.
- Capacity with 2 trays: ca. 2 - 4 Instrument sets per sterilization process.

TECHNICAL SPECIFICATIONS

Power requirements: 230 V / 50 Hz, 650 Watt
 Device dimensions: H 26 x W 31 x D 39.5 cm
 Internal Dimensions: H 7.4 x W 19.5 x D 29.5 cm
 Weight: 12 kg
 Volume: 4 liters

Autoclave Class S

The MELAG autoclaves Class S work with a pre- and post-vacuum in combination with a fractionated gravity system. This ensures that the air is completely removed from the sterilization chamber, from possible cavities of the instruments, and from the packaging. So these autoclaves fulfill the strict European standard EN 13060 and represent the current state of steam sterilization technique.

MELAtronic 23 EN

FEATURES

- Easy care, enamelled metal housing.
- Simple, convenient operation.
- Fully automatic, microprocessor-controlled program flow for monitoring and logging the sterilization parameters.
- Serial interface for connection to the report printer MELAprint, the CF Card Printer MELAflash or alternatively directly to your office PC.
- 4 programs are available, including a quick program for unwrapped and solid instruments.
- Water quality control.
- Optional water circuit or instrument gentle way system.
- Safety through the Info Display and electric door lock.
- Drying function and preheating function for better

1447 MELAtronic 23 EN

ACCESSORIES

1447.1 Tray for MELAtronic 23 EN



TECHNICAL SPECIFICATIONS

Power requirements: 230 V / 50 - 60 Hz / 1.600 watts
 Housing dimensions (cm): H 38 x W 52 x D 59
 Chamber volume: 19 liters
 Dimensions tray (cm): H 2 x W 19 x D 42
 Weight: 32 kg

SCOPE OF DELIVERY

The device is supplied with 3 trays (Art. 1447.1).
 Overall, there is a possibility of receiving 5 trays.

MELAtronic 15 EN+

A "plus" for the automatic drying!



SCOPE OF DELIVERY

The device is supplied with 3 trays (Art. 8235.401).

TECHNICAL SPECIFICATIONS

Power requirements: 230 V / 50 - 60 Hz / 1,500 W
 Housing dimensions (cm): W 44 x H 33 x D 50
 Chamber volume: 7 liters
 Dimensions tray (cm): H 2 x W 12 x D 35
 Weight: 19 kg

8235.81 MELAtronic 15 EN+ without MELAflash CF-Card

8235.8 MELAtronic 15 EN+ with MELAflash CF-Card

8235.82 Retrofitting drying MELAtronic 15 EN+

ACCESSORIES

8235.401 Tray for MELAtronic 15 EN+

Drying:

An essential prerequisite for the long-term storage of wrapped instruments is the dry removal of packaged sterilized goods after sterilization. According to the European standard for small autoclaves EN 13060, the residual moisture content of 0.2% for instruments must not be exceeded. These requirements are fully met the new MELAtronic 15 EN + through effective drying pump.

Water System:

For protecting and maintaining the value of the instruments quality of distilled or demineralized water used is crucial. By the MELAtronic 15 EN + conductance integrated water quality is checked. If the quality falls below a certain level, a start of the autoclave is only possible when new clean water is filled. In addition, this autoclave offers the

choice between water saving circulatory system or the one-way system where fresh water is used for each sterilization process.

Programs and times:

MELAtronic 15EN + offers four sterilization programs. A quick program (134 ° C) may unwrapped instruments after 18min be removed again from the autoclave. Also, the universal program (134 ° C) and the pioneer program (134 ° C) allow quick reuse of the sterilized wrapped instruments. The delicate program operates at a lower temperature (121 ° C) and is particularly suitable for heat-sensitive instruments.

Euroklav 29 VS+

Powerful class S autoclave with many advantages for the sterilization of simply packaged instruments. Equipped with versatile programs and documentation options. Cost-effective alternative for podiatry and pedicure practices, which do not require Class B autoclaves because of the instruments and packaging used.

FEATURES

- Functional and beautiful design
- Convenient and safe operation
- 4 sterilization programs and 2 test programs for function testing and documentation
- Includes quick program (20 min.) For unpacked instruments
- Control and monitoring of the program sequence by modern microprocessors
- Documentation of cycle data via printer, MELAflash CF card or practice network
- Convenient filling of the demineralized or distilled water thanks to the large opening of the storage tank and integrated filling aid
- Effortless removal of used water via quick coupling hoses
- Integrated measurement of water quality
- Easily adjustable options and special functions
- Possibility of switching from internal to external supply water supply with optional connection of a water treatment plant
- Large door lock for safe and easy opening and closing of the door
- Medical device Class IIb acc. Directive 93/42/EEC



1459.3 Euroklav 29 VS+ with MELAflash CF-Card

ACCESSORIES

1459.100 Tray for Euroklav 29 VS+

TECHNICAL SPECIFICATIONS

Electrical connection: 230 V, 50 Hz
 Power consumption: 2.100 watts
 Dimensions Housing: B 42.5 x H 48.5 x D 61 cm
 Chamber volume: 18 liters, Ø 25 x L 35 cm
 Measurements tray: W 19 x H 2 x D 29 cm
 Loading capacity: 3 kg instruments / 0.8 kg textiles
 Weight: 42 kg

SCOPE OF DELIVERY

- Euroklav 29 VS +
- Support for 5 tablets or 3 standard tray cassettes
- 3 trays (1459,100)
- MELAflash CF card



Autoclave Class B

Vacuklav 31 B+

According to the European standard EN 13060, class B is the highest class to which the strictest requirements for sterilization are made. Devices of this class are therefore fully usable for all loads (for example with hollow body and transmission instruments, turbines, etc.) occurring in practice, as well as for all types of sterilizing packages. This is achieved by the fractionated pre-vacuum process. Class B autoclaves are therefore absolutely future-proof and a worthwhile investment.



1459.4 Vacuklav 31 B+ with MELAflash CF-Card

FEATURES

- Functional and beautiful design
- Comfortable and safe operation
- Economical and efficient operation
- Patented air cooling system
- Powerful diaphragm pumps
- 5 sterilization programs and 2 test programs
- Functional testing and documentation
- Includes 2 fast programs (15 min for unpacked instruments and charge dependent 28-50 min for packaged instruments)
- Fully automatic sterilization process including drying
- Control and monitoring of the program sequence by modern microprocessors
- Logging of the sterilization parameters
- Documentation of cycle data via printer, MELAflash CF card or practice network
- Convenient filling of the demineralized or distilled water thanks to the large opening of the storage tank and integrated filling aid
- Effortless removal of used water via quick coupling hoses
- Integrated measurement of water quality
- Easily adjustable options and special functions
- Possibility of switching from internal to external supply water supply with optional connection of a water treatment plant
- Automatic pre-heating of the boiler
- Large door lock for safe and easy opening and closing of the door
- Medical device Class IIb acc. Directive 93/42/EEC

The patented air cooling of the vacuum pump as well as a simple operation and reliability characterize this all-rounder. At the same time, the innovative concept of Vacuklav 31 B + stands for fast operating times, versatile documentation possibilities and reliable sterilization results. The integrated monitoring of the feedwater and the electronic control of the Vacuklav 31 B + avoid operator errors and ensure the high quality of the sterilization and the maximum protection of your instruments.

TECHNICAL SPECIFICATIONS

Electrical connection: 230 V, 50-60 Hz
 Power consumption: 2,100 watts
 Dimensions Housing: B 42.5 x H 48.5 x D 63 cm
 Chamber volume: 17 liters, Ø 25 x L 35 cm
 Measurements Tray: W 19 x H 2 x D 29 cm
 Loading capacity: 5 kg instruments / 1.8 kg textiles
 Weight: 45 kg

SCOPE OF DELIVERY

- Vacuklav 31 B +
- Support for 5 trays or 3 standard tray trays
- 5 trays (1459,100)
- MELAflash CF card

ACCESSORIES

1459.100 Tray for Vacuklav 31 B+

Accessories for Autoclave



PROTOCOL PRINTER MELAPRINT

8235.5 Protocol printer

8235.501 Paper rolls, 5 pieces

Handy needle printer, logged and documented the parameters that are acquired on every sterilization passage automatically with MELAG autoclaves.



MELAFASH CF-CARD

8235.6 MELAflash CF-Card

3 components for easy and comfortable documentation of MELAG autoclave. MELAflash CF Card Printer; MELAflash CF Card. Power cable, cable for autoclave; MELAflash card reader



STERILIZATION REPORT

1443.5 1 block of 50 sheets

For documentation of the sterilization process.



MELAFOL TRANSPARENT HOSES

1459.6 1 roll, w 100 mm / L 200 m

1459.7 1 roll, w 75 mm / L 200 m

Transparent pipes, foil / paper for packing and storing sterilized goods. With color change indicator. Only for autoclave!



ADHESIVE FILM BAGS

8235.206 60 mm x 100 mm

8235.203 90 mm x 250 mm

8235.207 140 mm x 260 mm

Time-saving alternative to weld: Insert tool, remove protective strip, close bag. Done! 1 Pack of 200 pieces. Only for autoclave!



TAPE INDICATOR

1449.3 Klebebandindikator 1 Rolle, B 19 mm/L 50 m

A color change on the tape indicates that was sterilized. Can also be used as an indicator strip on the cassette and closed containers. Only for autoclave!

1449.4 Indicator strips for steam, 250 pieces 28,90 €

Accessories for hot-air sterilizer



FILM TUBES

1412.1 1 roll, W 75 mm / L 100 m

1412.2 1 roll, W 150 mm / L 100 m

Film tubes of special polyamide film for packing and storage of sterilization. Only for hot air sterilizers!



TAPE INDICATOR

1412.3 1 roll, W 19 mm / L 50 m

A color change on the tape indicates that was sterilized. Can also be used as an indicator strip on the cassette and closed containers. Only for hot air sterilizers!

1412.4 Indicator strips Hot air, 250 pieces




CUTTER BOWL

1415 Cutter bowl made of stainless steel

Dimensions: W 9 x H 1,7 x D 7,2 cm



Sealing Machine

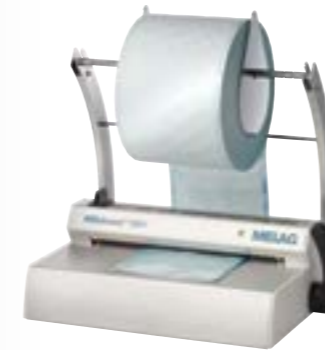


MELAseal 100+

Sealer for packaging sterilization. Efficient, reliable and easy to use. With visual indicator for monitoring temperature and sealing time.

1459.5 Sealing Machine MELAseal 100+

1459.003 Roll holder Comfort

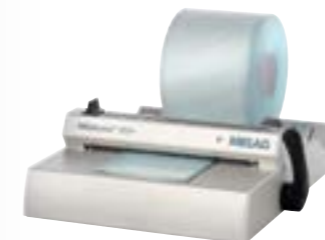

1459.5 Sealing Machine MELAseal 100+

1459.002 Roll holder Standard

Hygiene-Basic equipment

CLEAN-SET Basic

Basic equipment that is necessary to comply with the hygiene rules.



CLEAN-SET BASIS

8409.55 CLEAN-SET basic

SCOPE OF DELIVERY

- 8235.10 STERIL UL
- 8325.5 Elmasonic S podo basic
- 9220.100 Sterillium, 500 ml bottle
- 9220.200 Cutasept F, 250 ml-Spray bottle
- 9220.310 Bacillol plus, 500 ml-bottle
- 9220.403 Korsolex med AF, 2.000 ml-bottle

SCOPE OF DELIVERY

- 9220.601 Baktolin pure, 500 ml-bottle
- 9220.700 Mikrobac forte, 10 bags à 20 ml
- 9220.907 Instruments-RS-Spray, 200 ml
- 9220.914 Universal measuring cup 1.000 ml

Hygiene comfort equipment

CLEAN-SET Comfort

Comfortable full range for gentle sterilization of the instruments in a few minutes.



CLEAN-SET COMFORT

8409.56 CLEAN-SET Comfort



SCOPE OF DELIVERY

- 8235.8 MELAtronic 15 EN+ with MELAflash CF-Card
- 8325.5 Elmasonic S podo basic
- 9220.100 Sterillium, 500 ml-bottle
- 9220.200 Cutasept F, 250 ml-spray bottle
- 9220.310 Bacillol plus, 500 ml-bottle
- 9220.403 Korsolex med AF, 2000 ml-bottle
- 9220.502 Bacillol Tissues, Dispenser 100 wipes
- 9220.601 Baktolin pure, 500 ml-bottle
- 9220.613 Towel dispenser
- 9220.740 Quartamon Med, 2000 ml-bottle

SCOPE OF DELIVERY

- 9220.744 Gigasept Instru AF, 2.000 ml-bottle
- 9220.612 Battery-operated dispenser, Filling volume 1.400 m
- 9220.900 Disinfection tray
- 9220.902 Sterile cup, 200 ml
- 9220.905 Eurospender 3000, 500 ml-bottle
- 9220.907 Instruments-RS-Spray, 200 ml
- 9220.914 Universal-Measuring cup 1.000 ml

Masthead

EDITOR:

IONTO Health & Beauty GmbH
Greschbachstraße 3
76229 Karlsruhe
GERMANY

Section 1 General - Scope

1. Definitions

Entrepreneurs (commercial customers) within the meaning of the terms and conditions of business are natural or legal persons or partnerships with legal capacity with whom business relationships are entered into and who in the case of entering into a legal transaction act by way of pursuing a commercial of self-employed professional activity. Ultimate consumers within the meaning of the terms and conditions of business are natural persons with whom business relationships are entered into whereby such business relationships cannot be allocated to a commercial or independent professional activity. Customers within the meaning of the terms and conditions of business are both entrepreneurs and ultimate consumers. User within the meaning of the terms and conditions of business is IONTO Health & Beauty GmbH, Greschbachstraße 3, D-76229 Karlsruhe.

2. The terms and conditions of business apply to all current and future business relationships between users and customers as stated in their respective, current version. They apply to all details stated in this respect in brochures, price lists, advertisements etc. This applies irrespective of whether or not the details are made verbally, in writing, by facsimile or via the internet.

3. General terms and conditions of business that vary from, are contrary to or supplement our terms and conditions of business shall not become a part of the contract, including if known, unless the validity of such conditions is expressly approved in writing by the user.

4. The user is entitled to amend these general terms and conditions of business at any time with effect for the future.

Section 2 Entering into a contract

1. The user's offers are subject to change without notice and subject to another, written agreement. Technical modifications and alterations to the form, colour and / or weight are permitted as part of acceptable measures.

2. By way of placing an order, the customer states with binding force the intention to acquire the ordered goods. The user is entitled to accept the contractual offer that applies to the order within two weeks following receipt by the user. Acceptance may be stated wither in writing or by delivering the goods to the customer.

3. If the ultimate consumer orders the goods electronically, the user shall provide confirmation of the order without undue delay. Confirmation of receipt does not, at this stage, constitute binding acceptance of the order. The confirmation of receipt may be associated with the acceptance statement.

4. A contract shall be entered into subject to the correct and timely self-delivery by the user's supplier. This applies only in the event that non-delivery is not the user's responsibility. The customer shall be informed of the non-availability of the performance without undue delay. A potential counter-performance shall be reimbursed without undue delay.

5. Insofar as the ultimate consumer orders the goods electronically, the user's contract text shall be saved and on request forwarded to the customer including these general terms and conditions of business by e-mail.

6. Amendments to a contract that has been entered into shall be subject to the written form. This written form requirement may only be rescinded in writing.

7. With regard to the cancellation options for ultimate consumers, reference is expressly made to the separate cancellation policy (Section 4).

8. Section 454 and Section 455 BGB (German Civil Code) do not apply.

Section 3 Reservation of title

1. In the case of contracts entered into with end users, we reserve ownership of the goods up until payment in full of the purchase price. In the case of contracts entered into with entrepreneurs, the user reserves ownership of the goods up until settlement in full of all claims resulting from an ongoing business relationship.

2. The customer undertakes to treat the goods with due care. Insofar as servicing and inspection work is required, the customer is to regularly perform such work at its own cost. Similarly, the customer undertakes to provide notification of any changes in the place of residence or principal place of business.

3. If the value of the securities in place for the user exceed its claims in total by more than 20 %, at the customer's request the user undertakes to release securities in the sum of the excessive value up to the limit of 120 % of the outstanding receivables at the customer's discretion.

4. The customer undertakes to inform the user without undue delay of third party access to the goods, e.g. in the event of seizure and potential damage or destruction of the goods without undue delay.

5. In the event of conduct in breach of contract by the customer, in particular default in payment or in the event of violation of an obligation in accordance with sub-section 4 of this provision, the user shall be entitled to withdraw from the contract and demand the return of the goods.

6. The customer is entitled to sell the goods in the ordinary course of business. The customer assigns to the user at this point in time all claims in the sum of the invoice amount to which the customer is entitled against third parties from the sale. The user accepts the assignment. Following the assignment, the entrepreneur shall be authorised to collect the claim. The user reserves the right to collect the claim as soon as the entrepreneur fails to honour its payment obligations or does not properly honour these and defaults in payment. If the entrepreneur is wholly or partially in default with one or more payments, if the entrepreneur discontinues payment or if an application is filed for the institution of insolvency proceedings regarding its assets, it shall no longer be entitled to dispose of the reserved goods. The user shall then be entitled to withdraw

from the contract, take back the reserved goods or withdraw the authority granted to the customer to collect receivables from the selling on of goods and demand information about the recipients of the reserved goods and inform these of the assignment of the corresponding claim and collect the claim itself. The user is entitled to utilise the reserved goods, which may be reclaimed, on the open market. In the event of further processing, the user shall remain the manufacturer within the meaning of user Section 950 BGB.

7. Neither the assertion of the reservation of title nor the seizure of the delivery item by the user is deemed to constitute a withdrawal from contract.

Section 4 Right of cancellation and return

The ultimate user has a cancellation right in accordance with Section 355 BGB. Insofar reference is expressly made to the separate cancellation policy. Cancellation policy (for ultimate consumers only) Ultimate consumers may cancel their declaration to enter into a contract within two weeks, without stating a reason for such action, in text form (e.g. letter, facsimile or e-mail) or - if the item is surrendered to them prior to expiry of the period - by returning the item. The period shall commence following receipt of these directions in text form. Sending the cancellation or the item in good time shall be deemed sufficient to safeguard the cancellation period. The cancellation is to be directed to: IONTO Health & Beauty GmbH, Greschbachstraße 3, D-76229 Karlsruhe, facsimile 0721 9770-290, e-mail: info@ionto.de

In the event of a valid cancellation, mutually received performance is to be returned and, where applicable, benefits received are to be surrendered. If the received performance cannot be returned in full or in part by the ultimate consumer, or merely in a poorer condition, the ultimate consumer shall, where applicable, be required to provide compensation for lost value. In the case of surrendering items, this shall not apply if the deterioration in the item is solely attributable to reviewing it - as would be possible for example in a retail store. In other respects, the ultimate consumer may avoid the obligation to provide compensation for lost value involving the deterioration caused by the intended utilisation of the item by not taking possession of the item and avoiding any action that causes the value of the item to deteriorate. Items suitable for shipment as parcels from a sales price of EUR 40.00 are to be returned at our cost and risk. Items that are not suitable for shipment as parcels are to be collected from the ultimate consumer. Obligations to provide compensation for payments must be honoured within 30 days. The period shall commence for the ultimate consumer on sending the cancellation notice or the item, while for the user it shall commence upon receipt.

Section 5 Remuneration

1. The offered purchase price is binding. In dealings with ultimate consumers the turnover tax is included in the purchase price. In other respects, the stated (net) product price does not include shipping or insurance costs or expenses and taxes; these costs shall be charged separately.

2. The user reserves the right to adjust prices accordingly with effect for future business transactions in the event of changes to exchange rates, customs duties, taxes, freight and insurance costs or prime costs.

3. In the case of contracts entered into with companies, the user shall be entitled to amend the prices after entering into the contract too if the cost factors that are authoritative for the pricing are proven to have changed or the user's supplier has demonstrably increased the prices.

4. In the absence of agreements to the contrary, in each case the prices in accordance with the user's valid price list on the day of the delivery are deemed applicable (list price).

5. By way of the contract of purchase the customer undertakes to pay to the user the agreed purchase price and to accept the purchased item. If the customer refuses to accept the delivery for a reason that is not the user's responsibility, or if the customer withdraws from the contract without justification, the customer is to pay 25 % of the agreed purchase price, if the user does not insist on performance of the contract, as compensation for the incurred costs and loss of expected profit provided the customer does not furnish proof that the user did not sustain any damage at all or sustained significantly less damage than that reflected in the flat-rate. The user reserves the right to request compensation of the specific loss instead of contractual performance or the flat-rate.

Section 6 Terms and conditions of payment

1. In the absence of agreements to the contrary, the purchase price as well as other remuneration, ancillary performance prices and disbursed costs shall fall due for payment upon hand over of the object of purchase.

2. The customer may pay the purchase price via cash on delivery (plus EUR 2 delivery fee), in cash (2% trade discount), invoice (14 days without deduction) or via the SEPA direct debit procedure (2% trade discount). The direct debit shall be collected 3 days following the invoice date. The pre-notification is reduced to 3 days. The buyer assures that its bank account will have adequate cover. The user does not undertake to accept cheques or bills of exchange. If these are accepted, this shall only apply on account of performance. The customer shall only be permitted to deduct a trade discount in the case of receipt of payment in full and in good time by the user.

3. Following consultation, new customers are to place initial orders subject to advance payment. This also applies to followup deliveries up to a total net order volume of EUR 300.00.

4. The user reserves the right to hold back deliveries and / or services in the event of default in payment by the customer up until payment of potential outstanding claims, demand default interest as well as

compensation for additional damage sustained by the user as a result of the default.

5. The customer may only exercise a right of retention if its counter-claim is based on the same contractual relationship.

6. In dealings with the user the customer shall not be entitled to set off counter-claims provided the counter-claims are not undisputed or have not become res judicata.

7. If the customer has deceived the user upon placing the order regarding its credit rating, or if it is lacking, and if the user was unaware of such a situation, the user may withdraw from the contract without the setting of an additional period. If such circumstances occur after the placing of an order, the user only undertakes to render additional services subject to appropriate advance payment.

8. During the default, the ultimate consumer is to impose an interest rate of 5 percentage points above the current base lending rate on the obligation to pay money. During the default, the entrepreneur is to charge interest at 8 percentage points above the base lending rate on the obligation to pay money.

Section 7 Delivery

1. Delivery times are only agreed on an approximate basis. Agreed delivery periods shall only commence if the customer has honoured a duty to collaborate that may apply.

2. The user is entitled to provide partial deliveries. The location to which the goods are to be delivered is to be stated in the order confirmation. The goods shall be sent in suitable packaging that is customary in the trade for normal shipping. In the event of default in acceptance, all costs in that respect shall be borne by the customer, in particular warehousing and working hours costs.

3. In the absence of agreements to the contrary, the user reserves the right to select the transport route and mode of transport. As a matter of principle, all deliveries shall be made ex warehouse and at the customer's risk and cost. The customer shall also carry the transport risk if the shipping costs are, in an exceptional case, borne by the user. The goods may be insured against transport damage subject to reimbursement of expenses in that respect. If the customer is an ultimate consumer, the risk of accidental loss and accidental deterioration of the sold goods, including in the case of sales shipment, shall pass to the customer on hand over of the item.

4. The hand over shall be deemed to have taken place if the customer is in default.

5. Used equipment and equipment that is acquired in the online auction may only be delivered in Germany.

Section 8 Warranty

1. If the customer is an entrepreneur, the user shall provide subsequent improvement, at its own discretion, regarding faulty goods by way of subsequent performance or replacement delivery within a reasonable period. All articles and products that are stated separately in an invoice/order confirmation with an individual price are deemed independent items to which potential warranty rights apply separately.

2. If the customer is an ultimate consumer, the following procedure is agreed upon with consideration given to the user's economic interests to rectify a defect regarding the goods: In the case of products with a value of less than EUR 500 the ultimate consumer may initially only request replacement. If the value of the object of sale exceeds EUR 500.00, the user shall initially be entitled to an attempt to provide subsequent improvement within a reasonable period. A subsequent improvement period of 20 days is deemed reasonable. If the subsequent improvement is unacceptable in economic terms, subsequent performance shall be provided by way of replacement.

3. If the subsequent performance fails, as a matter of principle the customer shall be entitled, at its own discretion, to request a reduction of the remuneration or rescission of the contract. In the event of merely a minor breach of contract, in particular in the case of minor defects only, the customer shall not, however, be entitled to withdraw.

4. Entrepreneurs must provide the user with notification of obvious defects without undue delay, at the latest within a period of two weeks from receipt of the goods; otherwise, asserting a warranty claim is excluded. Providing notification in good time is deemed sufficient with regard to honouring the period. The burden of proof lies completely with the entrepreneur regarding any claim preconditions, in particular regarding the defect itself, the time of identifying the defect and providing notification of the defect in good time. Ultimate consumers must notify the user in writing of obvious defects within a period of two months from the time at which the condition of goods in breach of contract was identified. Receipt of notification by the user is deemed authoritative for safeguarding the period. If the ultimate consumer fails to provide such notification, the warranty rights shall expire two months after the ultimate consumer identified the defect. This does not apply in the event of fraudulent intent on the part of the user. The burden of proof lies with the ultimate consumer regarding the time at which the defect is identified. If the ultimate consumer was persuaded to purchase the item on the basis of inaccurate manufacturers details, the burden of proof regarding its purchasing decision shall lie with the ultimate consumer. In the event of used goods, the burden of proof lies with the ultimate consumer regarding the faulty nature of the item. With regard to complaints which, following a review, prove not to contain a defect that is the user's responsibility, or defects that are attributable to operating errors or inappropriate handling by the customer, the customer undertakes to reimburse all the user's costs incurred as a result of the review of the alleged defects.

5. If the customer intends to withdraw from the contract because of a defect in title or material defect

after subsequent performance has failed, in that respect the customer shall not be entitled to additionally claim for damages regarding the defect. If the customer intends to claim for compensatory damages after subsequent performance has failed, the goods shall remain in the customer's possession if this is acceptable. The compensatory damages shall be limited to the difference between the purchase price and the value of the faulty item. This shall not apply if the user caused the breach of contract as a result of fraudulent intent.

6. The guarantee period for entrepreneurs is one year from delivery. If the delivery items are used goods, and if these are sold expressly as used items, the guarantee shall be excluded. The user shall not be liable for wear-and-tear that is typically to be expected or loss in value or other malfunctions or limitations in use based on an item's age. A claim for replacement is excluded. The period of limitations for ultimate consumers is two years from delivery. In the case of used items, the period of limitations is one year from delivery of the goods. This shall not apply if the customer fails to provide notification of the defect in good time (sub-section 4 of this provision).

7. If the customer is an entrepreneur, as a matter of principle only the manufacturer's product description is deemed applicable as the quality of the goods. In addition, public statements, praise, goods and representative samples or advertising by the manufacturer do not constitute any details with regard to the quality of the goods as per agreement.

8. The user shall not provide the customer with any guarantees in a legal sense. This does not affect manufacturers' guarantees.

Section 9 Limitation on liability

1. In the case of minor negligent violations of obligations, the user's liability shall be limited to the average damage according to type of goods that is foreseeable, typical for such contracts and directly affects the goods. This also applies in the case of minor negligent violations of obligations on the part of our legal representatives or vicarious agents. The following applies in dealings with entrepreneurs:

2.1 In the case of minor negligence, the user shall only be liable in the case of violation of key contractual obligations (cardinal obligations) and in the case of personal damage and in accordance with the German Product Liability Act. In other respects, the user's pre-contractual, contractual and non-contractual liability is limited to intent and gross negligence, whereby the limitation on liability also applies in the event of culpability on the part of a vicarious agent or legal representative of the user.

2.2 Insofar as key contractual obligations are not violated as a result of gross negligence or intent, the user's liability for the occurred damage (including indirect damage) shall be limited to damage that is average, typical and foreseeable at the time of entering into the contract and in terms of amount shall be limited to the order value for each case of damage and damage associated with series. All further-reaching rights and claims, irrespective of the legal reason, are excluded, in particular compensation for indirect damage such as the loss of expected profits.

2.3 Liability shall be excluded insofar as the sustained damage could have been avoided by acceptable measures by the customer aimed at limiting damage.

3. The above-mentioned limitation on liability shall not apply to the customer's claims resulting from product liability. Furthermore, the limitation on liability shall not apply in the case of physical injury, detrimental effects on health or death of the customer attributable to the user.

4. The customer's claims for compensatory damages regarding a defect shall fall under the statute of limitations one year from delivery of the goods. This does not apply if the user is proven to have acted with fraudulent intent.

Section 10 Final provisions

1. The laws of the Federal Republic of Germany shall apply to this contract. The provisions of the UN Contract Law as well as other international laws do not apply.

2. If the customer is a merchant, legal person under public law or special federal funds, the user's registered office is deemed the place of jurisdiction for all disputes resulting from this contract. The same applies if the customer does not have a general place of jurisdiction in Germany, or the customer's place of residence or usual place of abode are not known at the time legal action is brought. The user reserves the right to bring legal action at the customer's principal place of business. The user's registered office in Karlsruhe is deemed the place of performance.

3. Customer data are subject to electronic data processing as part of the order processing. The user shall comply with the provisions of the German Data Protection Act when using personal data.

4. In the event that individual provisions of the contract entered into with the customer, including these general terms and conditions of business, are wholly or partially invalid, this shall not affect the other provisions. The wholly or partially invalid regulation is to be replaced by a regulation the economic success of which comes closest to that of the invalid provision.

IONTO Health & Beauty GmbH

Greschbachstrasse 3 . D-76229 Karlsruhe . Germany

IONTO Health & Beauty GmbH
Greschbachstraße 3
76229 Karlsruhe
Germany

FOOT CARE / PODIATRY

- > Phone 0721 9770 828
- > export@ionto.de
- > www.sueda.de/en/